

Titration report



Method

Method ID	SERCALKSAMPLES	Release state	Released
Name	TA of Seawater	Modified by	SystemInternal
Modified on	2/6/2024 11:26:27 AM		

Summary

Samples

Sample ID	Name	Result	State
1	GW5 8/7 0850 rep3, 11(1)	Sample size 60.069 g	
	Content R1	2.1109 mmol/L ⁽⁹⁾	OK
2	GW6 8/7 0815 rep3, 12(2)	Sample size 60.677 g	
	Content R1	1.9435 mmol/L ⁽⁹⁾	OK
3	CRM 209-0435, 13(3)	Sample size 56.581 g	
	Content R1	2.2111 mmol/L ⁽⁹⁾	OK
4	FTS 2023 8/16 0530-rep1, 14(4)	Sample size 58.572 g	
	Content R1	1.4967 mmol/L ⁽⁹⁾	OK
5	FTS 2023 8/16 0530-rep2, 15(5)	Sample size 59.903 g	
	Content R1	1.4949 mmol/L ⁽⁹⁾	OK
6	FTS 2023 8/16 0630 rep1, 16(6)	Sample size 58.182 g	
	Content R1	1.5103 mmol/L ⁽¹⁾	OK
7	FTS 2023 8/16 0630 rep2, 17(7)	Sample size 63.279 g	
	Content R1	1.5085 mmol/L ⁽⁹⁾	OK
8	FTS 2023 8/16 0730 rep1, 1(8)	Sample size 56.609 g	
	Content R1	1.5673 mmol/L ⁽¹⁾	OK
9	FTS 2023 8/16 0730 rep2, 2(9)	Sample size 57.505 g	
	Content R1	1.5598 mmol/L ⁽¹⁾	OK
10	FTS 2023 8/16 0830 0830 rep1, 3(10)	Sample size 58.970 g	
	Content R1	1.6619 mmol/L ⁽⁹⁾	OK
11	FTS 2023 8/16 0830 0830 rep2, 4(11)	Sample size 58.791 g	
	Content R1	1.6634 mmol/L ⁽⁹⁾	OK
12	FTS 2023 8/16 0930 rep1, 5(12)	Sample size 56.574 g	
	Content R1	1.8306 mmol/L ⁽⁹⁾	OK
13	FTS 2023 8/16 0930 rep2, 6(13)	Sample size 57.219 g	
	Content R1	1.8295 mmol/L ⁽⁹⁾	OK
14	FTS 2023 8/16 1030 (rep1), 7(14)	Sample size 58.001 g	
	Content R1	2.0248 mmol/L ⁽⁹⁾	OK
15	FTS 2023 8/16 1030 (rep2), 8(15)	Sample size 57.048 g	
	Content R1	2.0282 mmol/L ⁽⁹⁾	OK
16	FTS 2023 8/16 1130 rep1, 9(16)	Sample size 56.341 g	
	Content R1	2.1835 mmol/L ⁽⁹⁾	OK
17	FTS 2023 8/16 1130 rep2, 10(17)	Sample size 59.213 g	
	Content R1	2.1806 mmol/L ⁽¹⁾	OK
18	FTS 2023 8/16 1230 rep1, 11(18)	Sample size 59.157 g	
	Content R1	2.2690 mmol/L ⁽¹⁾	OK
19	FTS 2023 8/16 1230 rep2, 12(19)	Sample size 56.693 g	
	Content R1	2.2701 mmol/L ⁽⁹⁾	OK

Titration report



20	FTS 2023 8/16 13:30 rep1, 13(20)	Sample size	57.225 g	
		Content R1	2.2898 mmol/L ⁽⁹⁾	OK
21	FTS 2023 8/16 13:30 rep2, 14(21)	Sample size	58.141 g	
		Content R1	2.2862 mmol/L ⁽⁹⁾	OK
22	CRM 209-0435, 15(22)	Sample size	56.828 g	
		Content R1	2.2069 mmol/L ⁽⁹⁾	OK
23	FTS 2023 8/16 14:30 rep1, 16(23)	Sample size	56.073 g	
		Content R1	2.1798 mmol/L ⁽¹⁾	OK
24	FTS 2023 8/16 14:30 rep2, 17(24)	Sample size	61.350 g	
		Content R1	2.1424 mmol/L ⁽⁹⁾	OK
25	FTS 2023 8/16 15:30 rep1, 1(25)	Sample size	60.436 g	
		Content R1	2.0962 mmol/L ⁽⁹⁾	OK
26	FTS 2023 8/16 15:30 rep2, 2(26)	Sample size	60.202 g	
		Content R1	2.0971 mmol/L ⁽⁹⁾	OK
27	FTS 2023 8/16 16:30 rep1, 3(27)	Sample size	56.136 g	
		Content R1	2.0528 mmol/L ⁽⁹⁾	OK
28	FTS 2023 8/16 16:30 rep2, 4(28)	Sample size	58.294 g	
		Content R1	2.0554 mmol/L ⁽⁹⁾	OK
29	FTS 2023 8/16 17:30 rep1, 5(29)	Sample size	62.119 g	
		Content R1	1.9278 mmol/L ⁽⁹⁾	OK
30	FTS 2023 8/16 17:30 rep2, 6(30)	Sample size	57.120 g	
		Content R1	1.9278 mmol/L ⁽¹⁾	OK

Statistics

R1

Minimum	1.9536 mmol/L ^(1·9)	Minimum	1.4949 mmol/L ⁽⁹⁾
Standard deviation	0.2730 ^(1·9)	Rel. standard deviation	13.973 % ^(1·9)
Total samples	30	Samples excluded	0

Main Result Details

Result title	Content	Sample ID1	GW5 8/7 0850 rep3
Result name	R1	Sample ID2	11(1)
Result	2.1109 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content	Sample ID1	GW6 8/7 0815 rep3
Result name	R1	Sample ID2	12(2)
Result	1.9435 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	

Titration report

LabX

Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 2.2111 mmol/L ⁽⁹⁾

Sample ID1 CRM 209-0435
Sample ID2 13(3)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 1.4967 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 0530-rep1
Sample ID2 14(4)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 1.4949 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 0530-rep2
Sample ID2 15(5)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 1.5103 mmol/L ⁽¹⁾

Sample ID1 FTS 2023 8/16 0630 rep1
Sample ID2 16(6)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Titration report

LabX

Result title Content
Result name R1
Result 1.5085 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 0630 rep2
Sample ID2 17(7)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 1.5673 mmol/L ⁽¹⁾

Sample ID1 FTS 2023 8/16 0730 rep1
Sample ID2 1(8)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 1.5598 mmol/L ⁽¹⁾

Sample ID1 FTS 2023 8/16 0730 rep2
Sample ID2 2(9)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 1.6619 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 0830 0830 rep1
Sample ID2 3(10)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title Content
Result name R1
Result 1.6634 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 0830 0830 rep2
Sample ID2 4(11)
Result comment
Result state OK

Titration report

LabX

Result state comment

Nominal value

Compliance state

Released

Is result deleted

No

Calculation formula R1=Q*C/m

Result title Content

Result name R1

Result 1.8306 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 0930 rep1

Sample ID2 5(12)

Result comment

Result state OK

Result state comment

Nominal value

Compliance state

Released

Is result deleted

No

Calculation formula R1=Q*C/m

Result title Content

Result name R1

Result 1.8295 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 0930 rep2

Sample ID2 6(13)

Result comment

Result state OK

Result state comment

Nominal value

Compliance state

Released

Is result deleted

No

Calculation formula R1=Q*C/m

Result title Content

Result name R1

Result 2.0248 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 1030 (rep1)

Sample ID2 7(14)

Result comment

Result state OK

Result state comment

Nominal value

Compliance state

Released

Is result deleted

No

Calculation formula R1=Q*C/m

Result title Content

Result name R1

Result 2.0282 mmol/L ⁽⁹⁾

Sample ID1 FTS 2023 8/16 1030 (rep2)

Sample ID2 8(15)

Result comment

Result state OK

Result state comment

Nominal value

Compliance state

Released

Is result deleted

No

Titration report

LabX

Calculation formula R1=Q*C/m

Result title	Content	Sample ID1	FTS 2023 8/16 1130 rep1
Result name	R1	Sample ID2	9(16)
Result	2.1835 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content	Sample ID1	FTS 2023 8/16 1130 rep2
Result name	R1	Sample ID2	10(17)
Result	2.1806 mmol/L ⁽¹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content	Sample ID1	FTS 2023 8/16 1230 rep1
Result name	R1	Sample ID2	11(18)
Result	2.2690 mmol/L ⁽¹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content	Sample ID1	FTS 2023 8/16 1230 rep2
Result name	R1	Sample ID2	12(19)
Result	2.2701 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content	Sample ID1	FTS 2023 8/16 13:30 rep1
Result name	R1	Sample ID2	13(20)
Result	2.2898 mmol/L ⁽⁹⁾	Result comment	

Titration report



Result state	OK
Result state comment	
Nominal value	
Compliance state	Released
Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content
Result name	R1
Result	2.2862 mmol/L ⁽⁹⁾

Sample ID1	FTS 2023 8/16 13:30 rep2
Sample ID2	14(21)
Result comment	
Result state	OK
Result state comment	
Nominal value	
Compliance state	Released
Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content
Result name	R1
Result	2.2069 mmol/L ⁽⁹⁾

Sample ID1	CRM 209-0435
Sample ID2	15(22)
Result comment	
Result state	OK
Result state comment	
Nominal value	
Compliance state	Released
Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content
Result name	R1
Result	2.1798 mmol/L ⁽¹⁾

Sample ID1	FTS 2023 8/16 14:30 rep1
Sample ID2	16(23)
Result comment	
Result state	OK
Result state comment	
Nominal value	
Compliance state	Released
Is result deleted	No

Calculation formula R1=Q*C/m

Result title	Content
Result name	R1
Result	2.1424 mmol/L ⁽⁹⁾

Sample ID1	FTS 2023 8/16 14:30 rep2
Sample ID2	17(24)
Result comment	
Result state	OK
Result state comment	
Nominal value	
Compliance state	Released

Titration report



Is result deleted No

Calculation formula $R1=Q \cdot C/m$

Result title	Content	Sample ID1	FTS 2023 8/16 15:30 rep1
Result name	R1	Sample ID2	1(25)
Result	2.0962 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula $R1=Q \cdot C/m$

Result title	Content	Sample ID1	FTS 2023 8/16 15:30 rep2
Result name	R1	Sample ID2	2(26)
Result	2.0971 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula $R1=Q \cdot C/m$

Result title	Content	Sample ID1	FTS 2023 8/16 16:30 rep1
Result name	R1	Sample ID2	3(27)
Result	2.0528 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula $R1=Q \cdot C/m$

Result title	Content	Sample ID1	FTS 2023 8/16 16:30 rep2
Result name	R1	Sample ID2	4(28)
Result	2.0554 mmol/L ⁽⁹⁾	Result comment	
		Result state	OK
		Result state comment	
		Nominal value	
		Compliance state	Released
		Is result deleted	No

Calculation formula $R1=Q \cdot C/m$

Result title	Content	Sample ID1	FTS 2023 8/16 17:30 rep1
--------------	---------	------------	--------------------------

Titration report

LabX

Result name R1
Result 1.9278 mmol/L ⁽⁹⁾

Sample ID2 5(29)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula R1=Q*C/m

Result title Content
Result name R1
Result 1.9278 mmol/L ⁽¹⁾

Sample ID1 FTS 2023 8/16 17:30 rep2
Sample ID2 6(30)
Result comment
Result state OK
Result state comment
Nominal value
Compliance state Released
Is result deleted No

Calculation formula R1=Q*C/m

Raw Results

Scope 1/1, Sample 1/30 (GW5 8/7 0850 rep3; 11(1))

Sample (S1)

Creation time 2/21/2024 11:22:55 AM
Sample start time 2/21/2024 11:22:55 AM
Sample size 60.069 g
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 11:22:55 AM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 11:22:56 AM
T 24.9 °C
E 8.111 pH
tMe 01:13 min
ESTMe 3.076 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 11:24:10 AM
Terminated at Potential
EEQ1 4.529 pH ⁽⁹⁾
EHNv1 6.229 pH
Q1 0.126797 mmol ⁽⁹⁾

Titration report

LabX

VEQ1	1.266701 mL ⁽⁹⁾
EST	8.112 pH
VEX	0.794799 mL
VEND	2.06150 mL
t	23:52 min
QEND	0.206356 mmol
QEX	0.079559 mmol
EPD	6.752 pH
nEQ	2 ⁽⁹⁾

Content (R1)

Creation time	2/21/2024 11:48:10 AM ⁽⁹⁾
Result	2.1109 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 11:48:10 AM
---------------	-----------------------

Scope 1/1, Sample 2/30 (GW6 8/7 0815 rep3; 12(2))

Sample (S1)

Creation time	2/21/2024 11:48:35 AM
Sample start time	2/21/2024 11:48:35 AM
Sample size	60.677 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 11:48:36 AM
---------------	-----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 11:48:36 AM
T	25.1 °C
E	8.092 pH
tMe	01:12 min
ESTMe	4.812 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 11:49:48 AM
Terminated at	Potential
EEQ1	4.571 pH ⁽⁹⁾
EHN1	6.230 pH
Q1	0.117926 mmol ⁽⁹⁾
VEQ1	1.178086 mL ⁽⁹⁾
EST	8.094 pH
VEX	0.815664 mL
VEND	1.99375 mL

Titration report

LabX

t	23:24 min
QEND	0.199574 mmol
QEX	0.081648 mmol
EPD	6.710 pH
nEQ	2 ⁽⁹⁾

Content (R1)

Creation time	2/21/2024 12:13:18 PM ⁽⁹⁾
Result	1.9435 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 12:13:18 PM
---------------	-----------------------

Scope 1/1, Sample 3/30 (CRM 209-0435; 13(3))

Sample (S1)

Creation time	2/21/2024 12:13:29 PM
Sample start time	2/21/2024 12:13:29 PM
Sample size	56.581 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.022 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 12:13:29 PM
---------------	-----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 12:13:29 PM
T	25.1 °C
E	8.034 pH
tMe	04:16 min
ESTMe	6.742 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 12:17:46 PM
Terminated at	Potential
EEQ1	4.476 pH ⁽⁹⁾
EHN1	6.119 pH
Q1	0.125108 mmol ⁽⁹⁾
VEQ1	1.249835 mL ⁽⁹⁾
EST	8.027 pH
VEX	0.737665 mL
VEND	1.98750 mL
t	22:52 min
QEND	0.198949 mmol
QEX	0.073840 mmol
EPD	6.662 pH

Titration report

LabX

nEQ 1

Content (R1)

Creation time 2/21/2024 12:40:43 PM ⁽⁹⁾

Result 2.2111 mmol/L ⁽⁹⁾

Decimal places 4 ⁽⁹⁾

Substance None ⁽⁹⁾

M 1 g/mol ⁽⁹⁾

z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 12:40:44 PM

Scope 1/1, Sample 4/30 (FTS 2023 8/16 0530-rep1; 14(4))

Sample (S1)

Creation time 2/21/2024 12:40:55 PM

Sample start time 2/21/2024 12:40:55 PM

Sample size 58.572 g

Correction factor 1.0

Temperature 25.0 °C

Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 12:40:55 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 12:40:55 PM

T 25.3 °C

E 7.482 pH

tMe 01:12 min

ESTMe 5.144 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 12:42:08 PM

Terminated at Potential

EEQ1 4.629 pH ⁽⁹⁾

EHN1 6.217 pH

Q1 0.087662 mmol ⁽⁹⁾

VEQ1 0.875744 mL ⁽⁹⁾

EST 7.495 pH

VEX 0.795756 mL

VEND 1.67150 mL

t 20:47 min

QEND 0.167317 mmol

QEX 0.079655 mmol

EPD 6.463 pH

nEQ 1

Content (R1)

Creation time 2/21/2024 1:03:01 PM ⁽⁹⁾

Titration report

LabX

Result	1.4967 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 1:03:01 PM
---------------	----------------------

Scope 1/1, Sample 5/30 (FTS 2023 8/16 0530-rep2; 15(5))

Sample (S1)

Creation time	2/21/2024 1:03:12 PM
Sample start time	2/21/2024 1:03:12 PM
Sample size	59.903 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 1:03:13 PM
---------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 1:03:13 PM
T	25.3 °C
E	7.489 pH
tMe	01:12 min
ESTMe	5.489 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 1:04:26 PM
Terminated at	Potential
EEQ1	4.628 pH ⁽⁹⁾
EHN1	6.200 pH
Q1	0.089546 mmol ⁽⁹⁾
VEQ1	0.894569 mL ⁽⁹⁾
EST	7.501 pH
VEX	0.791181 mL
VEND	1.68575 mL
t	20:48 min
QEND	0.168744 mmol
QEX	0.079197 mmol
EPD	6.461 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 1:25:19 PM ⁽⁹⁾
Result	1.4949 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾

Titration report

LabX

z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 1:25:19 PM

Scope 1/1, Sample 6/30 (FTS 2023 8/16 0630 rep1; 16(6))

Sample (S1)

Creation time 2/21/2024 1:25:30 PM

Sample start time 2/21/2024 1:25:30 PM

Sample size 58.182 g

Correction factor 1.0

Temperature 25.0 °C

Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 1:25:30 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 1:25:30 PM

T 25.2 °C

E 7.332 pH

tMe 01:12 min

ESTMe 4.295 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 1:26:43 PM

Terminated at Potential

EEQ1 4.592 pH

EHN1 6.183 pH

Q1 0.087870 mmol

VEQ1 0.877820 mL

EST 7.347 pH

VEX 0.750074 mL

VEND 1.64575 mL

t 20:17 min

QEND 0.164740 mmol

QEX 0.075082 mmol

EPD 6.424 pH

nEQ 2 ⁽²²⁾

Content (R1)

Creation time 2/21/2024 1:47:06 PM ⁽¹⁾

Result 1.5103 mmol/L ⁽¹⁾

Decimal places 4 ⁽¹⁾

Substance None ⁽¹⁾

M 1 g/mol ⁽¹⁾

z 1 ⁽¹⁾

Rinse (Rinse1)

Creation time 2/21/2024 1:47:06 PM

Titration report

LabX

Scope 1/1, Sample 7/30 (FTS 2023 8/16 0630 rep2; 17(7))

Sample (S1)

Creation time	2/21/2024 1:47:17 PM
Sample start time	2/21/2024 1:47:17 PM
Sample size	63.279 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 1:47:17 PM
----------------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 1:47:18 PM
T	25.2 °C
E	7.373 pH
tMe	01:11 min
ESTMe	4.398 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 1:48:30 PM
Terminated at	Potential
EEQ1	4.602 pH ⁽⁹⁾
EHN1	6.178 pH
Q1	0.095454 mmol ⁽⁹⁾
VEQ1	0.953589 mL ⁽⁹⁾
EST	7.387 pH
VEX	0.872411 mL
VEND	1.82600 mL
t	21:49 min
QEND	0.182783 mmol
QEX	0.087328 mmol
EPD	6.476 pH
nEQ	2 ⁽⁹⁾

Content (R1)

Creation time	2/21/2024 2:10:24 PM ⁽⁹⁾
Result	1.5085 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 2:10:25 PM
----------------------	----------------------

Scope 1/1, Sample 8/30 (FTS 2023 8/16 0730 rep1; 1(8))

Sample (S1)

Creation time	2/21/2024 2:34:14 PM
----------------------	----------------------

Titration report

LabX

Sample start time 2/21/2024 2:34:14 PM
Sample size 56.609 g
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 2:34:14 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 2:34:14 PM
T 25.2 °C
E 7.313 pH
tMe 01:12 min
ESTMe 4.350 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 2:35:27 PM
Terminated at Potential
EEQ1 4.600 pH
EHN1 6.192 pH
Q1 0.088724 mmol
VEQ1 0.886350 mL
EST 7.326 pH
VEX 0.771005 mL
VEND 1.64250 mL
t 20:17 min
QEND 0.164414 mmol
QEX 0.077178 mmol
EPD 6.440 pH
nEQ 2

Content (R1)

Creation time 2/21/2024 2:55:49 PM ⁽¹⁾
Result 1.5673 mmol/L ⁽¹⁾
Decimal places 4 ⁽¹⁾
Substance None ⁽¹⁾
M 1 g/mol ⁽¹⁾
z 1 ⁽¹⁾

Rinse (Rinse1)

Creation time 2/21/2024 2:55:49 PM

Scope 1/1, Sample 9/30 (FTS 2023 8/16 0730 rep2; 2(9))

Sample (S1)

Creation time 2/21/2024 2:56:00 PM
Sample start time 2/21/2024 2:56:00 PM
Sample size 57.505 g
Correction factor 1.0
Temperature 25.0 °C

Titration report

LabX

Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 2:56:01 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 2:56:01 PM

T 25.3 °C

E 7.359 pH

tMe 01:12 min

ESTMe 4.578 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 2:57:14 PM

Terminated at Potential

EEQ1 4.612 pH

EHN1 6.194 pH

Q1 0.089697 mmol

VEQ1 0.896071 mL

EST 7.374 pH

VEX 0.771683 mL

VEND 1.69325 mL

t 20:47 min

QEND 0.169494 mmol

QEX 0.077245 mmol

EPD 6.445 pH

nEQ 2 ⁽²²⁾

Content (R1)

Creation time 2/21/2024 3:18:07 PM ⁽¹⁾

Result 1.5598 mmol/L ⁽¹⁾

Decimal places 4 ⁽¹⁾

Substance None ⁽¹⁾

M 1 g/mol ⁽¹⁾

z 1 ⁽¹⁾

Rinse (Rinse1)

Creation time 2/21/2024 3:18:07 PM

Scope 1/1, Sample 10/30 (FTS 2023 8/16 0830 0830 rep1; 3(10))

Sample (S1)

Creation time 2/21/2024 3:18:18 PM

Sample start time 2/21/2024 3:18:18 PM

Sample size 58.970 g

Correction factor 1.0

Temperature 25.0 °C

Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 3:18:18 PM

Titration report

LabX

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 3:18:18 PM
T	25.2 °C
E	7.460 pH
tMe	01:12 min
ESTMe	4.525 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 3:19:31 PM
Terminated at	Potential
EEQ1	4.579 pH ⁽⁹⁾
EHN1	6.207 pH
Q1	0.098003 mmol ⁽⁹⁾
VEQ1	0.979050 mL ⁽⁹⁾
EST	7.473 pH
VEX	0.782950 mL
VEND	1.76200 mL
t	21:20 min
QEND	0.176376 mmol
QEX	0.078373 mmol
EPD	6.517 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 3:40:56 PM ⁽⁹⁾
Result	1.6619 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 3:40:57 PM
----------------------	----------------------

Scope 1/1, Sample 11/30 (FTS 2023 8/16 0830 0830 rep2; 4(11))

Sample (S1)

Creation time	2/21/2024 3:41:08 PM
Sample start time	2/21/2024 3:41:08 PM
Sample size	58.791 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 3:41:08 PM
----------------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 3:41:08 PM
T	25.2 °C

Titration report

LabX

E 7.478 pH
tMe 01:12 min
ESTMe 4.496 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 3:42:21 PM
Terminated at Potential
EEQ1 4.575 pH ⁽⁹⁾
EHNv1 6.205 pH
Q1 0.097791 mmol ⁽⁹⁾
VEQ1 0.976929 mL ⁽⁹⁾
EST 7.490 pH
VEX 0.785821 mL
VEND 1.76275 mL
t 21:19 min
QEND 0.176451 mmol
QEX 0.078661 mmol
EPD 6.515 pH
nEQ 1

Content (R1)

Creation time 2/21/2024 4:03:45 PM ⁽⁹⁾
Result 1.6634 mmol/L ⁽⁹⁾
Decimal places 4 ⁽⁹⁾
Substance None ⁽⁹⁾
M 1 g/mol ⁽⁹⁾
z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 4:03:45 PM

Scope 1/1, Sample 12/30 (FTS 2023 8/16 0930 rep1; 5(12))

Sample (S1)

Creation time 2/21/2024 4:03:56 PM
Sample start time 2/21/2024 4:03:56 PM
Sample size 56.574 g
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 4:03:56 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 4:03:56 PM
T 25.2 °C
E 7.633 pH
tMe 01:12 min
ESTMe 4.615 pH

Titration report

LabX

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 4:05:10 PM
Terminated at	Potential
EEQ1	4.570 pH ⁽⁹⁾
EHN1	6.227 pH
Q1	0.103566 mmol ⁽⁹⁾
VEQ1	1.034627 mL ⁽⁹⁾
EST	7.643 pH
VEX	0.773123 mL
VEND	1.80775 mL
t	21:49 min
QEND	0.180956 mmol
QEX	0.077390 mmol
EPD	6.573 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 4:27:03 PM ⁽⁹⁾
Result	1.8306 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 4:27:04 PM
----------------------	----------------------

Scope 1/1, Sample 13/30 (FTS 2023 8/16 0930 rep2; 6(13))

Sample (S1)

Creation time	2/21/2024 4:27:15 PM
Sample start time	2/21/2024 4:27:15 PM
Sample size	57.219 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 4:27:15 PM
----------------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 4:27:15 PM
T	25.2 °C
E	7.649 pH
tMe	01:12 min
ESTMe	4.320 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 4:28:28 PM
Terminated at	Potential

Titration report

LabX

EEQ1	4.579 pH ⁽⁹⁾
EHN1	6.230 pH
Q1	0.104685 mmol ⁽⁹⁾
VEQ1	1.045807 mL ⁽⁹⁾
EST	7.660 pH
VEX	0.766193 mL
VEND	1.81200 mL
t	21:49 min
QEND	0.181381 mmol
QEX	0.076696 mmol
EPD	6.584 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 4:50:22 PM ⁽⁹⁾
Result	1.8295 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 4:50:23 PM
---------------	----------------------

Scope 1/1, Sample 14/30 (FTS 2023 8/16 1030 (rep1); 7(14))

Sample (S1)

Creation time	2/21/2024 4:50:34 PM
Sample start time	2/21/2024 4:50:34 PM
Sample size	58.001 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 4:50:34 PM
---------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 4:50:34 PM
T	25.2 °C
E	7.519 pH
tMe	01:12 min
ESTMe	4.710 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 4:51:47 PM
Terminated at	Potential
EEQ1	4.563 pH ⁽⁹⁾
EHN1	6.215 pH
Q1	0.117440 mmol ⁽⁹⁾
VEQ1	1.173228 mL ⁽⁹⁾

Titration report

LabX

EST	7.529 pH
VEX	0.773272 mL
VEND	1.94650 mL
t	22:49 min
QEND	0.194845 mmol
QEX	0.077405 mmol
EPD	6.621 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 5:14:42 PM ⁽⁹⁾
Result	2.0248 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 5:14:42 PM
---------------	----------------------

Scope 1/1, Sample 15/30 (FTS 2023 8/16 1030 (rep2); 8(15))

Sample (S1)

Creation time	2/21/2024 5:14:53 PM
Sample start time	2/21/2024 5:14:53 PM
Sample size	57.048 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 5:14:53 PM
---------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 5:14:54 PM
T	25.2 °C
E	7.569 pH
tMe	01:12 min
ESTMe	4.743 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 5:16:06 PM
Terminated at	Potential
EEQ1	4.555 pH ⁽⁹⁾
EHN1	6.225 pH
Q1	0.115702 mmol ⁽⁹⁾
VEQ1	1.155861 mL ⁽⁹⁾
EST	7.579 pH
VEX	0.782139 mL
VEND	1.93800 mL
t	22:51 min

Titration report

LabX

QEND	0.193994 mmol
QEX	0.078292 mmol
EPD	6.625 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 5:39:03 PM ⁽⁹⁾
Result	2.0282 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 5:39:03 PM
---------------	----------------------

Scope 1/1, Sample 16/30 (FTS 2023 8/16 1130 rep1; 9(16))

Sample (S1)

Creation time	2/21/2024 5:39:14 PM
Sample start time	2/21/2024 5:39:14 PM
Sample size	56.341 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 5:39:14 PM
---------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 5:39:14 PM
T	25.2 °C
E	7.624 pH
tMe	01:12 min
ESTMe	4.695 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 5:40:27 PM
Terminated at	Potential
EEQ1	4.560 pH ⁽⁹⁾
EHN1	6.228 pH
Q1	0.123018 mmol ⁽⁹⁾
VEQ1	1.228954 mL ⁽⁹⁾
EST	7.633 pH
VEX	0.766546 mL
VEND	1.99550 mL
t	23:20 min
QEND	0.199750 mmol
QEX	0.076731 mmol
EPD	6.664 pH
nEQ	1

Titration report

LabX

Content (R1)

Creation time	2/21/2024 6:03:53 PM ⁽⁹⁾
Result	2.1835 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 6:03:53 PM
----------------------	----------------------

Scope 1/1, Sample 17/30 (FTS 2023 8/16 1130 rep2; 10(17))

Sample (S1)

Creation time	2/21/2024 6:04:04 PM
Sample start time	2/21/2024 6:04:04 PM
Sample size	59.213 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 6:04:05 PM
----------------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 6:04:05 PM
T	25.2 °C
E	7.662 pH
tMe	01:12 min
ESTMe	4.737 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 6:05:18 PM
Terminated at	Potential
EEQ1	4.561 pH
EHN1	6.226 pH
Q1	0.129119 mmol
VEQ1	1.289905 mL
EST	7.670 pH
VEX	0.804619 mL
VEND	2.11600 mL
t	24:21 min
QEND	0.211812 mmol
QEX	0.080542 mmol
EPD	6.693 pH
nEQ	2 ⁽²²⁾

Content (R1)

Creation time	2/21/2024 6:29:44 PM ⁽¹⁾
Result	2.1806 mmol/L ⁽¹⁾

Titration report

LabX

Decimal places	4 ⁽¹⁾
Substance	None ⁽¹⁾
M	1 g/mol ⁽¹⁾
z	1 ⁽¹⁾

Rinse (Rinse1)

Creation time	2/21/2024 6:29:44 PM
---------------	----------------------

Scope 1/1, Sample 18/30 (FTS 2023 8/16 1230 rep1; 11(18))

Sample (S1)

Creation time	2/21/2024 6:29:55 PM
Sample start time	2/21/2024 6:29:55 PM
Sample size	59.157 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 6:29:56 PM
---------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 6:29:56 PM
T	25.3 °C
E	7.742 pH
tMe	01:12 min
ESTMe	4.718 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 6:31:09 PM
Terminated at	Potential
EEQ1	4.554 pH
EHN1	6.236 pH
Q1	0.134230 mmol
VEQ1	1.340956 mL
EST	7.747 pH
VEX	0.799870 mL
VEND	2.16900 mL
t	24:52 min
QEND	0.217117 mmol
QEX	0.080067 mmol
EPD	6.726 pH
nEQ	2 ⁽²²⁾

Content (R1)

Creation time	2/21/2024 6:56:06 PM ⁽¹⁾
Result	2.2690 mmol/L ⁽¹⁾
Decimal places	4 ⁽¹⁾
Substance	None ⁽¹⁾
M	1 g/mol ⁽¹⁾
z	1 ⁽¹⁾

Titration report

LabX

Rinse (Rinse1)

Creation time 2/21/2024 6:56:07 PM

Scope 1/1, Sample 19/30 (FTS 2023 8/16 1230 rep2; 12(19))

Sample (S1)

Creation time 2/21/2024 6:56:45 PM

Sample start time 2/21/2024 6:56:45 PM

Sample size 56.693 g

Correction factor 1.0

Temperature 25.0 °C

Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 6:56:46 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 6:56:46 PM

T 25.2 °C

E 7.792 pH

tMe 01:12 min

ESTMe 4.731 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 6:57:59 PM

Terminated at Potential

EEQ1 4.562 pH ⁽⁹⁾

EHN1 6.244 pH

Q1 0.128699 mmol ⁽⁹⁾

VEQ1 1.285706 mL ⁽⁹⁾

EST 7.796 pH

VEX 0.766294 mL

VEND 2.05200 mL

t 23:50 min

QEND 0.205405 mmol

QEX 0.076706 mmol

EPD 6.716 pH

nEQ 1

Content (R1)

Creation time 2/21/2024 7:21:55 PM ⁽⁹⁾

Result 2.2701 mmol/L ⁽⁹⁾

Decimal places 4 ⁽⁹⁾

Substance None ⁽⁹⁾

M 1 g/mol ⁽⁹⁾

z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 7:21:55 PM

Scope 1/1, Sample 20/30 (FTS 2023 8/16 13:30 rep1; 13(20))

Titration report

LabX

Sample (S1)

Creation time	2/21/2024 7:22:06 PM
Sample start time	2/21/2024 7:22:06 PM
Sample size	57.225 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 7:22:06 PM
----------------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 7:22:06 PM
T	25.3 °C
E	7.923 pH
tMe	01:12 min
ESTMe	4.727 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 7:23:19 PM
Terminated at	Potential
EEQ1	4.536 pH ⁽⁹⁾
EHN1	6.251 pH
Q1	0.131033 mmol ⁽⁹⁾
VEQ1	1.309022 mL ⁽⁹⁾
EST	7.926 pH
VEX	0.794728 mL
VEND	2.10375 mL
t	24:24 min
QEND	0.210585 mmol
QEX	0.079552 mmol
EPD	6.749 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 7:47:48 PM ⁽⁹⁾
Result	2.2898 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 7:47:49 PM
----------------------	----------------------

Scope 1/1, Sample 21/30 (FTS 2023 8/16 13:30 rep2; 14(21))

Sample (S1)

Creation time	2/21/2024 7:48:00 PM
Sample start time	2/21/2024 7:48:00 PM

Titration report

LabX

Sample size 58.141 g
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 7:48:00 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 7:48:00 PM
T 25.4 °C
E 7.929 pH
tMe 01:12 min
ESTMe 4.736 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 7:49:13 PM
Terminated at Potential
EEQ1 4.547 pH ⁽⁹⁾
EHN1 6.250 pH
Q1 0.132920 mmol ⁽⁹⁾
VEQ1 1.327869 mL ⁽⁹⁾
EST 7.933 pH
VEX 0.782881 mL
VEND 2.11075 mL
t 24:23 min
QEND 0.211286 mmol
QEX 0.078366 mmol
EPD 6.760 pH
nEQ 1

Content (R1)

Creation time 2/21/2024 8:13:42 PM ⁽⁹⁾
Result 2.2862 mmol/L ⁽⁹⁾
Decimal places 4 ⁽⁹⁾
Substance None ⁽⁹⁾
M 1 g/mol ⁽⁹⁾
z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 8:13:42 PM

Scope 1/1, Sample 22/30 (CRM 209-0435; 15(22))

Sample (S1)

Creation time 2/21/2024 8:13:53 PM
Sample start time 2/21/2024 8:13:53 PM
Sample size 56.828 g
Correction factor 1.0
Temperature 25.0 °C
Density 1.022 g/mL

Titration report

LabX

Titration stand (Titrationstand1)

Creation time 2/21/2024 8:13:53 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 8:13:54 PM

T 25.4 °C

E 7.899 pH

tMe 01:12 min

ESTMe 4.416 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 8:15:07 PM

Terminated at Potential

EEQ1 4.477 pH ⁽⁹⁾

EHNV1 6.118 pH

Q1 0.125414 mmol ⁽⁹⁾

VEQ1 1.252884 mL ⁽⁹⁾

EST 7.900 pH

VEX 0.740866 mL

VEND 1.99375 mL

t 22:50 min

QEND 0.199574 mmol

QEX 0.074161 mmol

EPD 6.655 pH

nEQ 1

Content (R1)

Creation time 2/21/2024 8:38:02 PM ⁽⁹⁾

Result 2.2069 mmol/L ⁽⁹⁾

Decimal places 4 ⁽⁹⁾

Substance None ⁽⁹⁾

M 1 g/mol ⁽⁹⁾

z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 8:38:02 PM

Scope 1/1, Sample 23/30 (FTS 2023 8/16 14:30 rep1; 16(23))

Sample (S1)

Creation time 2/21/2024 8:38:13 PM

Sample start time 2/21/2024 8:38:13 PM

Sample size 56.073 g

Correction factor 1.0

Temperature 25.0 °C

Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 8:38:13 PM

Measure (normal) [1] (Measurenormal1)

Titration report

LabX

Creation time	2/21/2024 8:38:14 PM
T	25.3 °C
E	8.010 pH
tMe	01:12 min
ESTMe	4.945 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 8:39:26 PM
Terminated at	Potential
EEQ1	4.581 pH
EHN1	6.288 pH
Q1	0.122227 mmol
VEQ1	1.221052 mL
EST	8.012 pH
VEX	0.733688 mL
VEND	1.97700 mL
t	23:20 min
QEND	0.197898 mmol
QEX	0.073442 mmol
EPD	6.785 pH
nEQ	2 ⁽²²⁾

Content (R1)

Creation time	2/21/2024 9:02:52 PM ⁽¹⁾
Result	2.1798 mmol/L ⁽¹⁾
Decimal places	4 ⁽¹⁾
Substance	None ⁽¹⁾
M	1 g/mol ⁽¹⁾
z	1 ⁽¹⁾

Rinse (Rinse1)

Creation time	2/21/2024 9:02:52 PM
----------------------	----------------------

Scope 1/1, Sample 24/30 (FTS 2023 8/16 14:30 rep2; 17(24))

Sample (S1)

Creation time	2/21/2024 9:03:03 PM
Sample start time	2/21/2024 9:03:03 PM
Sample size	61.350 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 9:03:03 PM
----------------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 9:03:04 PM
T	25.2 °C
E	7.971 pH
tMe	01:12 min

Titration report

LabX

ESTMe 4.776 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 9:04:17 PM
Terminated at Potential
EEQ1 4.542 pH ⁽⁹⁾
EHNv1 6.247 pH
Q1 0.131434 mmol ⁽⁹⁾
VEQ1 1.313026 mL ⁽⁹⁾
EST 7.974 pH
VEX 0.804224 mL
VEND 2.11725 mL
t 24:25 min
QEND 0.211937 mmol
QEX 0.080503 mmol
EPD 6.762 pH
nEQ 2 ⁽⁹⁾

Content (R1)

Creation time 2/21/2024 9:28:47 PM ⁽⁹⁾
Result 2.1424 mmol/L ⁽⁹⁾
Decimal places 4 ⁽⁹⁾
Substance None ⁽⁹⁾
M 1 g/mol ⁽⁹⁾
z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 9:28:47 PM

Scope 1/1, Sample 25/30 (FTS 2023 8/16 15:30 rep1; 1(25))

Sample (S1)

Creation time 2/21/2024 9:28:59 PM
Sample start time 2/21/2024 9:28:59 PM
Sample size 60.436 g
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 9:28:59 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 9:28:59 PM
T 25.3 °C
E 7.977 pH
tMe 01:11 min
ESTMe 4.735 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 9:30:12 PM

Titration report

LabX

Terminated at	Potential
EEQ1	4.552 pH ⁽⁹⁾
EHN1	6.247 pH
Q1	0.126683 mmol ⁽⁹⁾
VEQ1	1.265567 mL ⁽⁹⁾
EST	7.980 pH
VEX	0.841183 mL
VEND	2.10675 mL
t	24:25 min
QEND	0.210886 mmol
QEX	0.084202 mmol
EPD	6.740 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 9:54:42 PM ⁽⁹⁾
Result	2.0962 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 9:54:42 PM
----------------------	----------------------

Scope 1/1, Sample 26/30 (FTS 2023 8/16 15:30 rep2; 2(26))

Sample (S1)

Creation time	2/21/2024 9:54:53 PM
Sample start time	2/21/2024 9:54:53 PM
Sample size	60.202 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 9:54:54 PM
----------------------	----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 9:54:54 PM
T	25.1 °C
E	7.995 pH
tMe	01:12 min
ESTMe	4.667 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 9:56:07 PM
Terminated at	Potential
EEQ1	4.554 pH ⁽⁹⁾
EHN1	6.250 pH
Q1	0.126248 mmol ⁽⁹⁾

Titration report

LabX

VEQ1	1.261216 mL ⁽⁹⁾
EST	7.997 pH
VEX	0.791034 mL
VEND	2.05225 mL
t	23:54 min
QEND	0.205430 mmol
QEX	0.079183 mmol
EPD	6.744 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 10:20:06 PM ⁽⁹⁾
Result	2.0971 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 10:20:06 PM
---------------	-----------------------

Scope 1/1, Sample 27/30 (FTS 2023 8/16 16:30 rep1; 3(27))

Sample (S1)

Creation time	2/21/2024 10:25:32 PM
Sample start time	2/21/2024 10:25:32 PM
Sample size	56.136 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 10:25:32 PM
---------------	-----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 10:25:32 PM
T	25.0 °C
E	7.897 pH
tMe	01:11 min
ESTMe	4.713 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 10:26:45 PM
Terminated at	Potential
EEQ1	4.560 pH ⁽⁹⁾
EHN1	6.248 pH
Q1	0.115236 mmol ⁽⁹⁾
VEQ1	1.151204 mL ⁽⁹⁾
EST	7.901 pH
VEX	0.772796 mL
VEND	1.92400 mL

Titration report

LabX

t	22:54 min
QEND	0.192592 mmol
QEX	0.077357 mmol
EPD	6.677 pH
nEQ	1

Content (R1)

Creation time	2/21/2024 10:49:45 PM ⁽⁹⁾
Result	2.0528 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 10:49:45 PM
---------------	-----------------------

Scope 1/1, Sample 28/30 (FTS 2023 8/16 16:30 rep2; 4(28))

Sample (S1)

Creation time	2/21/2024 10:55:17 PM
Sample start time	2/21/2024 10:55:17 PM
Sample size	58.294 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 10:55:18 PM
---------------	-----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 10:55:18 PM
T	25.0 °C
E	7.914 pH
tMe	01:12 min
ESTMe	4.712 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 10:56:31 PM
Terminated at	Potential
EEQ1	4.551 pH ⁽⁹⁾
EHN1	6.246 pH
Q1	0.119815 mmol ⁽⁹⁾
VEQ1	1.196956 mL ⁽⁹⁾
EST	7.919 pH
VEX	0.792544 mL
VEND	1.98950 mL
t	23:23 min
QEND	0.199149 mmol
QEX	0.079334 mmol
EPD	6.702 pH

Titration report

LabX

nEQ 1

Content (R1)

Creation time 2/21/2024 11:19:59 PM ⁽⁹⁾
Result 2.0554 mmol/L ⁽⁹⁾
Decimal places 4 ⁽⁹⁾
Substance None ⁽⁹⁾
M 1 g/mol ⁽⁹⁾
z 1 ⁽⁹⁾

Rinse (Rinse1)

Creation time 2/21/2024 11:20:00 PM

Scope 1/1, Sample 29/30 (FTS 2023 8/16 17:30 rep1; 5(29))

Sample (S1)

Creation time 2/21/2024 11:20:11 PM
Sample start time 2/21/2024 11:20:11 PM
Sample size 62.119 g
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 2/21/2024 11:20:11 PM

Measure (normal) [1] (Measurenormal1)

Creation time 2/21/2024 11:20:11 PM
T 25.2 °C
E 7.867 pH
tMe 01:11 min
ESTMe 4.701 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/21/2024 11:21:24 PM
Terminated at Potential
EEQ1 4.555 pH ⁽⁹⁾
EHN1 6.233 pH
Q1 0.119756 mmol ⁽⁹⁾
VEQ1 1.196363 mL ⁽⁹⁾
EST 7.873 pH
VEX 0.848637 mL
VEND 2.04500 mL
t 23:54 min
QEND 0.204704 mmol
QEX 0.084949 mmol
EPD 6.685 pH
nEQ 2 ⁽⁹⁾

Content (R1)

Creation time 2/21/2024 11:45:24 PM ⁽⁹⁾

Titration report

LabX

Result	1.9278 mmol/L ⁽⁹⁾
Decimal places	4 ⁽⁹⁾
Substance	None ⁽⁹⁾
M	1 g/mol ⁽⁹⁾
z	1 ⁽⁹⁾

Rinse (Rinse1)

Creation time	2/21/2024 11:45:24 PM
----------------------	-----------------------

Scope 1/1, Sample 30/30 (FTS 2023 8/16 17:30 rep2; 6(30))

Sample (S1)

Creation time	2/21/2024 11:45:49 PM
Sample start time	2/21/2024 11:45:49 PM
Sample size	57.120 g
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/21/2024 11:45:50 PM
----------------------	-----------------------

Measure (normal) [1] (Measurenormal1)

Creation time	2/21/2024 11:45:50 PM
T	25.2 °C
E	7.901 pH
tMe	01:12 min
ESTMe	4.678 pH

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/21/2024 11:47:03 PM
Terminated at	Potential
EEQ1	4.583 pH
EHN1	6.243 pH
Q1	0.110118 mmol
VEQ1	1.100077 mL
EST	7.906 pH
VEX	0.765581 mL
VEND	1.85975 mL
t	22:25 min
QEND	0.186161 mmol
QEX	0.076635 mmol
EPD	6.645 pH
nEQ	2

Content (R1)

Creation time	2/22/2024 12:09:33 AM ⁽¹⁾
Result	1.9278 mmol/L ⁽¹⁾
Decimal places	4 ⁽¹⁾
Substance	None ⁽¹⁾
M	1 g/mol ⁽¹⁾

Titration report



z

1 ⁽¹⁾

Rinse (Rinse1)

Creation time 2/22/2024 12:09:34 AM

After Sample Scopes

Park (Park1)

Creation time 2/22/2024 10:01:57 AM

Measure (normal) [1] Scope 1/1, Sample 1/30

Measured Values Scope 1/1, Sample 1/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	3.076	19.3	
2	4.841	19.9	1.765
3	5.333	20.7	0.492
4	5.678	21.5	0.345
5	5.975	22.2	0.297
6	6.245	22.8	0.270
7	6.482	23.3	0.237
8	6.673	23.6	0.191
9	6.860	23.9	0.187
10	7.018	24.1	0.158
11	7.165	24.3	0.147
12	7.291	24.4	0.126
13	7.402	24.5	0.111
14	7.488	24.6	0.086
15	7.571	24.7	0.083
16	7.641	24.7	0.070
17	7.703	24.8	0.062
18	7.756	24.8	0.053
19	7.803	24.8	0.047
20	7.839	24.8	0.036
21	7.874	24.9	0.035
22	7.904	24.9	0.030
23	7.932	24.9	0.028
24	7.955	24.9	0.023
25	7.974	24.9	0.019
26	7.992	24.9	0.018
27	8.006	24.9	0.014
28	8.019	24.9	0.013
29	8.030	24.9	0.011
30	8.041	24.9	0.011
31	8.050	24.9	0.009
32	8.058	24.9	0.008
33	8.064	24.9	0.006
34	8.070	24.9	0.006
35	8.075	24.9	0.005


Titration report

LabX

36	8.080	24.9	0.005
37	8.083	24.9	0.003
38	8.087	24.9	0.004
39	8.090	24.9	0.003
40	8.092	24.9	0.002
41	8.095	24.9	0.003
42	8.096	24.9	0.001
43	8.098	24.9	0.002
44	8.099	24.9	0.001
45	8.101	24.9	0.002
46	8.102	24.9	0.001
47	8.103	24.9	0.001
48	8.104	24.9	0.001
49	8.105	24.9	0.001
50	8.106	24.9	0.001
51	8.106	24.9	0.000
52	8.107	24.9	0.001
53	8.108	24.9	0.001
54	8.108	24.9	0.000
55	8.108	24.9	0.000
56	8.109	24.9	0.001
57	8.110	24.9	0.001
58	8.110	24.9	0.000
59	8.110	24.9	0.000
60	8.111	24.9	0.001

Titration (EQP) [1] Scope 1/1, Sample 1/30


E - V Curve Scope 1/1, Sample 1/30

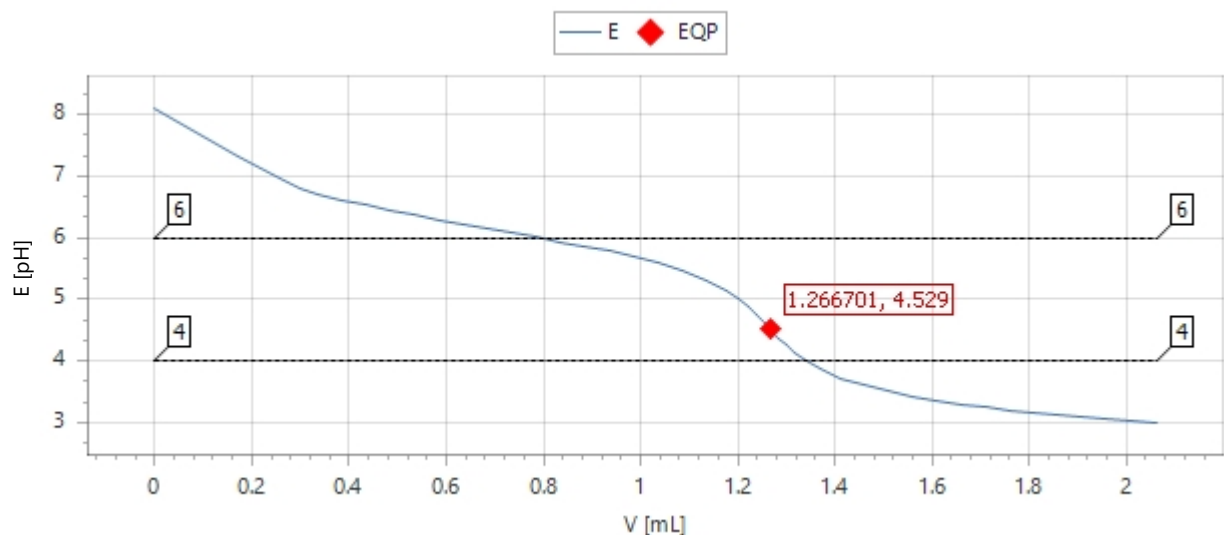

dE/dV - V Curve Scope 1/1, Sample 1/30

Titration report

LabX

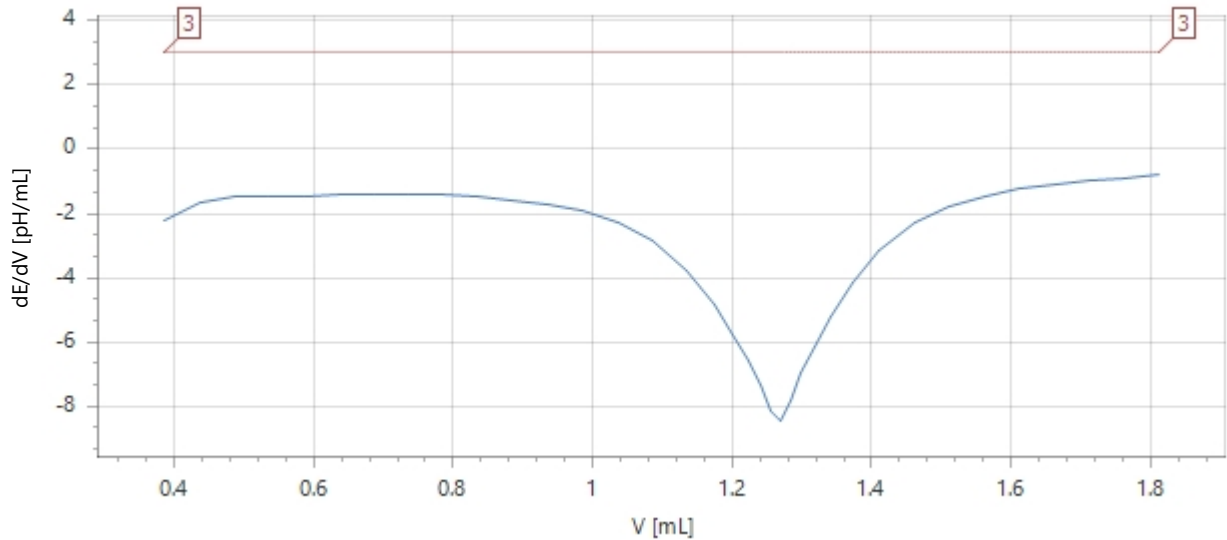

Measured Values Scope 1/1, Sample 1/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	8.112		24.9			
29	0.17150	7.313	0.17150	24.9	-0.799		
59	0.25725	6.969	0.08575	24.9	-0.344		
90	0.30000	6.802	0.04275	24.9	-0.167		
180	0.33600	6.697	0.03600	24.9	-0.105		
211	0.38600	6.614	0.05000	24.9	-0.083	-2.20	
241	0.43600	6.528	0.05000	24.9	-0.086	-1.65	
271	0.48600	6.447	0.05000	24.9	-0.081	-1.46	
302	0.53600	6.370	0.05000	24.9	-0.077	-1.47	
332	0.58600	6.296	0.05000	24.9	-0.074	-1.45	
362	0.63600	6.225	0.05000	24.9	-0.071	-1.41	
392	0.68600	6.157	0.05000	24.9	-0.068	-1.39	
422	0.73600	6.084	0.05000	24.9	-0.073	-1.39	
453	0.78600	6.016	0.05000	24.9	-0.068	-1.43	
483	0.83600	5.941	0.05000	24.9	-0.075	-1.48	
513	0.88600	5.866	0.05000	24.9	-0.075	-1.57	
543	0.93600	5.783	0.05000	24.9	-0.083	-1.70	
574	0.98600	5.690	0.05000	24.9	-0.093	-1.93	
604	1.03600	5.587	0.05000	24.9	-0.103	-2.27	
634	1.08600	5.463	0.05000	24.9	-0.124	-2.84	
664	1.13600	5.310	0.05000	24.9	-0.153	-3.76	
695	1.17500	5.153	0.03900	24.9	-0.157	-4.85	
725	1.20075	5.024	0.02575	24.9	-0.129	-5.76	
755	1.22250	4.889	0.02175	24.9	-0.135	-6.61	
785	1.24000	4.760	0.01750	24.9	-0.129	-7.35	
816	1.25550	4.628	0.01550	24.9	-0.132	-8.14	
	1.266701 ⁽⁹⁾	4.529 ⁽⁹⁾					EQP
846	1.26925	4.506	0.01375	24.9	-0.122	-8.43	
876	1.28400	4.375	0.01475	24.9	-0.131	-7.82	
906	1.29925	4.253	0.01525	24.9	-0.122	-6.98	
936	1.31825	4.123	0.01900	24.9	-0.130	-6.15	

Titration report

LabX

967	1.34125	3.995	0.02300	24.9	-0.128	-5.21
997	1.37175	3.863	0.03050	24.9	-0.132	-4.15
1027	1.41150	3.730	0.03975	24.9	-0.133	-3.13
1057	1.46150	3.604	0.05000	24.9	-0.126	-2.29
1088	1.51150	3.506	0.05000	24.9	-0.098	-1.78
1118	1.56150	3.424	0.05000	24.9	-0.082	-1.46
1148	1.61150	3.357	0.05000	24.9	-0.067	-1.24
1178	1.66150	3.297	0.05000	24.9	-0.060	-1.09
1209	1.71150	3.245	0.05000	24.9	-0.052	-0.98
1239	1.76150	3.199	0.05000	24.9	-0.046	-0.88
1269	1.81150	3.156	0.05000	24.9	-0.043	-0.80
1299	1.86150	3.118	0.05000	24.9	-0.038	
1330	1.91150	3.083	0.05000	24.9	-0.035	
1360	1.96150	3.050	0.05000	24.9	-0.033	
1390	2.01150	3.020	0.05000	24.9	-0.030	
1420	2.06150	2.992	0.05000	24.9	-0.028	

Measure (normal) [1] Scope 1/1, Sample 2/30

Measured Values Scope 1/1, Sample 2/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.812	23.1	
2	5.091	23.3	0.279
3	5.388	23.6	0.297
4	5.679	23.9	0.291
5	5.946	24.1	0.267
6	6.193	24.3	0.247
7	6.422	24.5	0.229
8	6.605	24.6	0.183
9	6.788	24.8	0.183
10	6.951	24.8	0.163
11	7.090	24.9	0.139
12	7.219	25.0	0.129
13	7.326	25.0	0.107
14	7.419	25.0	0.093
15	7.496	25.1	0.077
16	7.572	25.1	0.076
17	7.638	25.1	0.066
18	7.694	25.1	0.056
19	7.746	25.1	0.052
20	7.790	25.1	0.044
21	7.827	25.1	0.037
22	7.857	25.1	0.030
23	7.887	25.1	0.030
24	7.913	25.1	0.026
25	7.934	25.1	0.021
26	7.954	25.1	0.020
27	7.970	25.1	0.016
28	7.985	25.1	0.015

Titration report

LabX


29	7.997	25.1	0.012
30	8.008	25.1	0.011
31	8.017	25.1	0.009
32	8.026	25.1	0.009
33	8.034	25.1	0.008
34	8.041	25.1	0.007
35	8.047	25.1	0.006
36	8.052	25.1	0.005
37	8.056	25.1	0.004
38	8.061	25.1	0.005
39	8.064	25.1	0.003
40	8.067	25.1	0.003
41	8.070	25.1	0.003
42	8.073	25.1	0.003
43	8.075	25.1	0.002
44	8.077	25.1	0.002
45	8.079	25.1	0.002
46	8.080	25.1	0.001
47	8.082	25.1	0.002
48	8.083	25.1	0.001
49	8.085	25.1	0.002
50	8.085	25.1	0.000
51	8.087	25.1	0.002
52	8.088	25.1	0.001
53	8.088	25.1	0.000
54	8.089	25.1	0.001
55	8.090	25.1	0.001
56	8.090	25.1	0.000
57	8.091	25.1	0.001
58	8.091	25.1	0.000
59	8.092	25.1	0.001
60	8.092	25.1	0.000

Titration (EQP) [1] Scope 1/1, Sample 2/30


E - V Curve Scope 1/1, Sample 2/30

Titration report

LabX

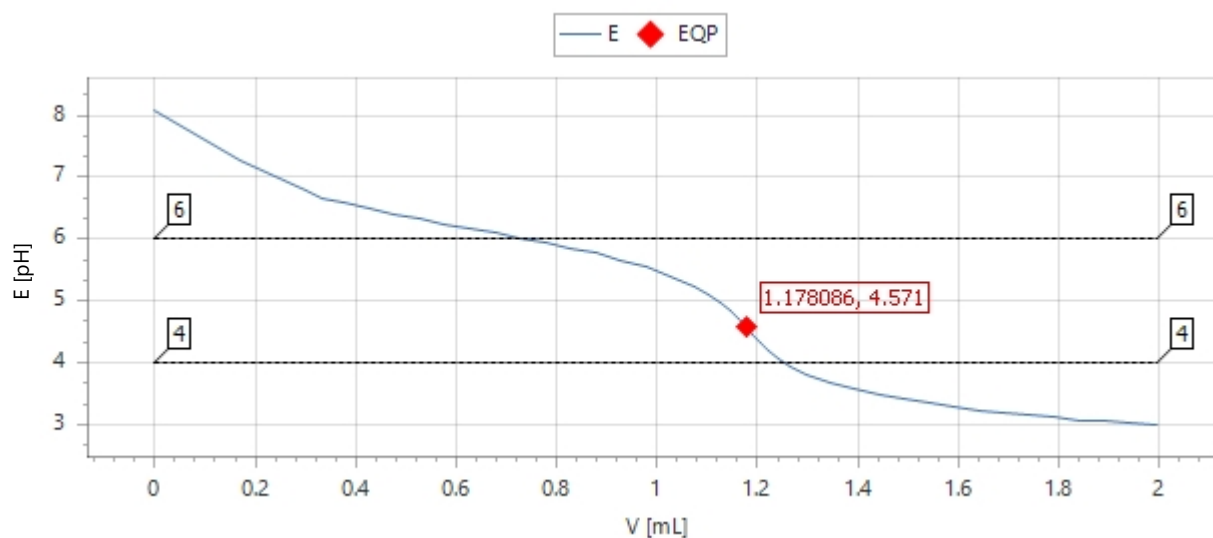

dE/dV - V Curve Scope 1/1, Sample 2/30

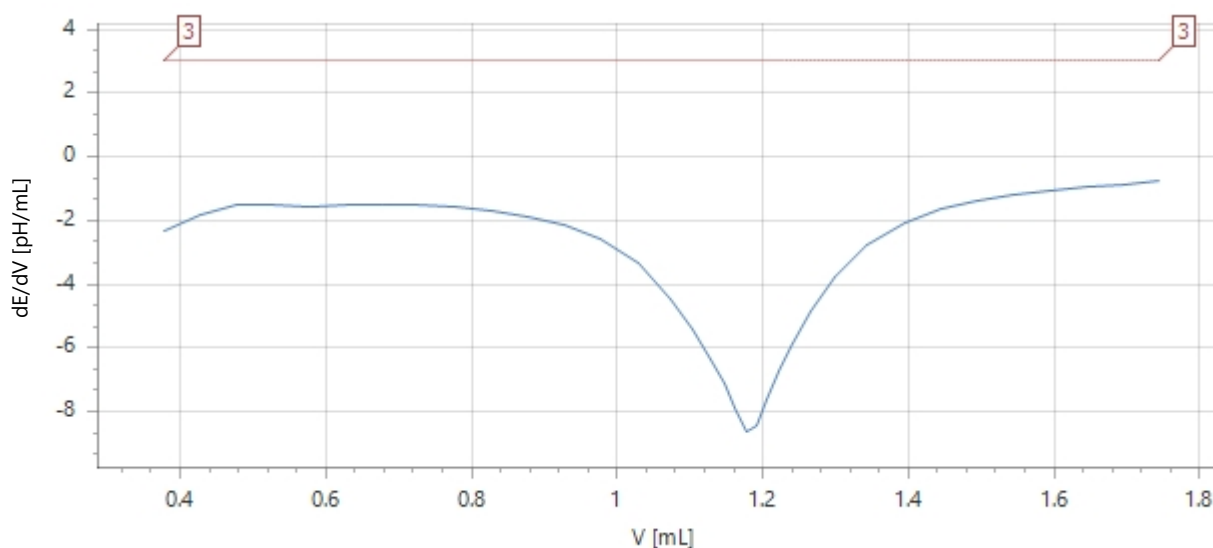

Measured Values Scope 1/1, Sample 2/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	8.094		25.1			
32	0.17150	7.275	0.17150	25.1	-0.819		
62	0.25725	6.944	0.08575	25.1	-0.331		
92	0.30000	6.774	0.04275	25.1	-0.170		
183	0.33325	6.659	0.03325	25.1	-0.115		
213	0.37825	6.583	0.04500	25.1	-0.076	-2.34	
243	0.42825	6.492	0.05000	25.1	-0.091	-1.80	
273	0.47825	6.410	0.05000	25.1	-0.082	-1.54	
304	0.52825	6.327	0.05000	25.1	-0.083	-1.53	
334	0.57825	6.246	0.05000	25.1	-0.081	-1.55	
364	0.62825	6.174	0.05000	25.1	-0.072	-1.52	
394	0.67825	6.099	0.05000	25.1	-0.075	-1.52	
425	0.72825	6.018	0.05000	25.1	-0.081	-1.54	
455	0.77825	5.940	0.05000	25.1	-0.078	-1.60	

Titration report

LabX

485	0.82825	5.857	0.05000	25.1	-0.083	-1.70	EQP
515	0.87825	5.769	0.05000	25.1	-0.088	-1.87	
546	0.92825	5.668	0.05000	25.1	-0.101	-2.14	
576	0.97825	5.552	0.05000	25.1	-0.116	-2.60	
606	1.02825	5.410	0.05000	25.1	-0.142	-3.34	
636	1.07400	5.246	0.04575	25.1	-0.164	-4.44	
667	1.10375	5.111	0.02975	25.1	-0.135	-5.42	
697	1.12725	4.976	0.02350	25.1	-0.135	-6.33	
727	1.14625	4.844	0.01900	25.1	-0.132	-7.11	
757	1.16225	4.714	0.01600	25.1	-0.130	-7.90	
788	1.17650	4.586	0.01425	25.1	-0.128	-8.60	
	1.178086 ⁽⁹⁾	4.571 ⁽⁹⁾					
818	1.19025	4.460	0.01375	25.1	-0.126	-8.43	
848	1.20475	4.336	0.01450	25.1	-0.124	-7.59	
878	1.22150	4.206	0.01675	25.1	-0.130	-6.75	
908	1.24100	4.081	0.01950	25.1	-0.125	-5.88	
939	1.26650	3.950	0.02550	25.1	-0.131	-4.85	
969	1.29925	3.821	0.03275	25.1	-0.129	-3.80	
999	1.34375	3.686	0.04450	25.1	-0.135	-2.79	
1029	1.39375	3.572	0.05000	25.1	-0.114	-2.09	
1060	1.44375	3.480	0.05000	25.1	-0.092	-1.65	
1090	1.49375	3.404	0.05000	25.1	-0.076	-1.38	
1120	1.54375	3.340	0.05000	25.1	-0.064	-1.18	
1150	1.59375	3.283	0.05000	25.1	-0.057	-1.05	
1181	1.64375	3.233	0.05000	25.1	-0.050	-0.93	
1211	1.69375	3.189	0.05000	25.1	-0.044	-0.85	
1241	1.74375	3.148	0.05000	25.1	-0.041	-0.77	
1272	1.79375	3.111	0.05000	25.1	-0.037		
1302	1.84375	3.077	0.05000	25.1	-0.034		
1332	1.89375	3.045	0.05000	25.1	-0.032		
1362	1.94375	3.016	0.05000	25.1	-0.029		
1392	1.99375	2.988	0.05000	25.1	-0.028		

Measure (normal) [1] Scope 1/1, Sample 3/30

Measured Values Scope 1/1, Sample 3/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	6.742	18.6	
2	6.806	19.3	0.064
3	6.931	20.2	0.125
4	7.131	21.1	0.200
5	7.344	21.9	0.213
6	7.532	22.6	0.188
7	7.692	23.2	0.160
8	7.813	23.6	0.121
9	7.922	23.9	0.109
10	8.008	24.2	0.086
11	8.078	24.4	0.070
12	8.136	24.5	0.058

Titration report

LabX

13	8.177	24.6	0.041
14	8.214	24.7	0.037
15	8.245	24.8	0.031
16	8.259	24.9	0.014
17	8.270	24.9	0.011
18	8.273	25.0	0.003
19	8.278	25.0	0.005
20	8.278	25.0	0.000
21	8.272	25.0	-0.006
22	8.263	25.0	-0.009
23	8.264	25.1	0.001
24	8.255	25.1	-0.009
25	8.246	25.1	-0.009
26	8.237	25.1	-0.009
27	8.229	25.1	-0.008
28	8.220	25.1	-0.009
29	8.211	25.1	-0.009
30	8.202	25.1	-0.009
31	8.193	25.1	-0.009
32	8.184	25.1	-0.009
33	8.176	25.1	-0.008
34	8.167	25.1	-0.009
35	8.163	25.1	-0.004
36	8.154	25.1	-0.009
37	8.146	25.1	-0.008
38	8.138	25.1	-0.008
39	8.131	25.1	-0.007
40	8.125	25.1	-0.006
41	8.119	25.1	-0.006
42	8.113	25.1	-0.006
43	8.106	25.1	-0.007
44	8.100	25.1	-0.006
45	8.094	25.1	-0.006
46	8.090	25.1	-0.004
47	8.087	25.1	-0.003
48	8.079	25.1	-0.008
49	8.075	25.1	-0.004
50	8.071	25.1	-0.004
51	8.066	25.1	-0.005
52	8.064	25.1	-0.002
53	8.059	25.1	-0.005
54	8.055	25.1	-0.004
55	8.052	25.1	-0.003
56	8.048	25.1	-0.004
57	8.044	25.1	-0.004
58	8.040	25.1	-0.004
59	8.038	25.1	-0.002
60	8.034	25.1	-0.004


Titration (EQP) [1]

Scope 1/1, Sample 3/30


Titration report

LabX

E - V Curve Scope 1/1, Sample 3/30

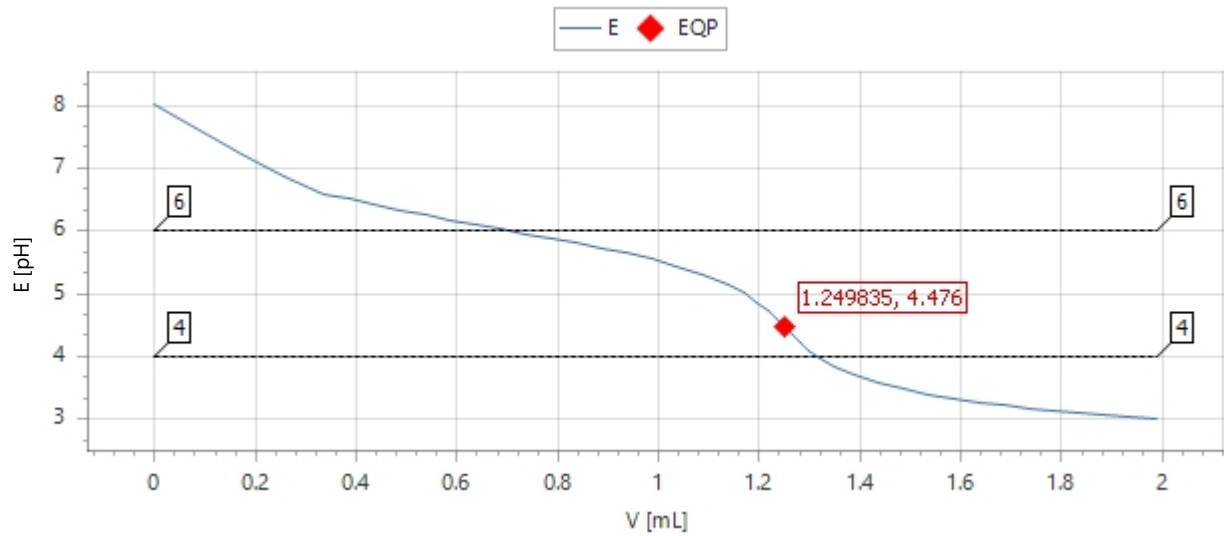

dE/dV - V Curve Scope 1/1, Sample 3/30

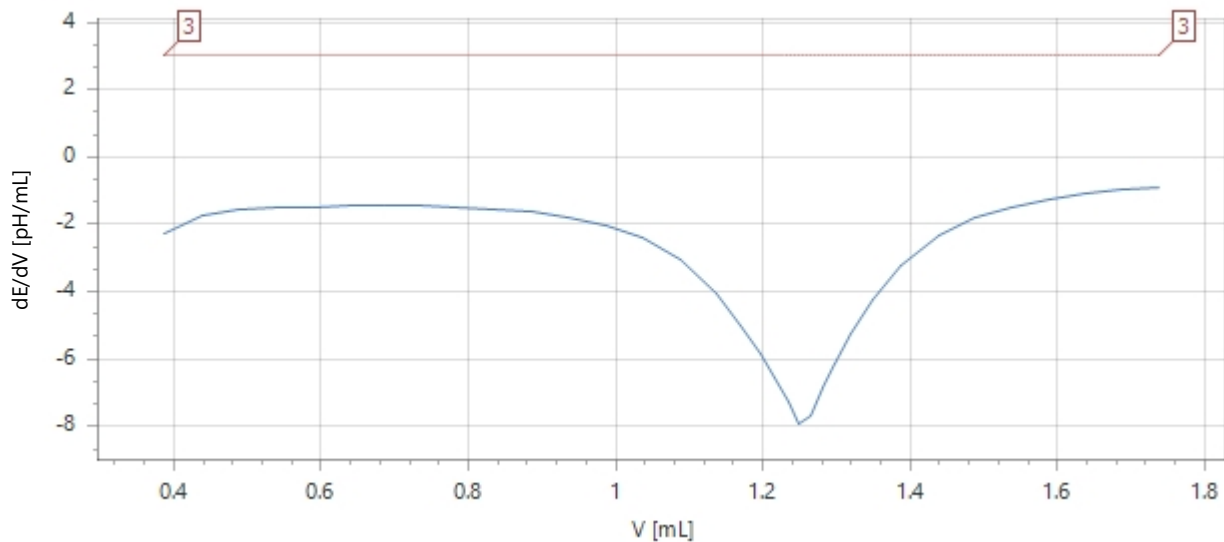

Measured Values Scope 1/1, Sample 3/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	8.027		25.1			
29	0.17150	7.217	0.17150	25.1	-0.810		
59	0.25725	6.869	0.08575	25.1	-0.348		
90	0.30000	6.704	0.04275	25.1	-0.165		
180	0.33725	6.597	0.03725	25.1	-0.107		
211	0.38725	6.507	0.05000	25.1	-0.090	-2.27	
241	0.43725	6.415	0.05000	25.1	-0.092	-1.74	
271	0.48725	6.326	0.05000	25.1	-0.089	-1.55	
302	0.53725	6.248	0.05000	25.1	-0.078	-1.52	
332	0.58725	6.174	0.05000	25.1	-0.074	-1.48	
362	0.63725	6.101	0.05000	25.1	-0.073	-1.43	
392	0.68725	6.031	0.05000	25.1	-0.070	-1.42	

Titration report

LabX

422	0.73725	5.958	0.05000	25.1	-0.073	-1.44
453	0.78725	5.885	0.05000	25.1	-0.073	-1.48
483	0.83725	5.810	0.05000	25.1	-0.075	-1.54
513	0.88725	5.729	0.05000	25.1	-0.081	-1.63
543	0.93725	5.644	0.05000	25.1	-0.085	-1.78
574	0.98725	5.548	0.05000	25.1	-0.096	-2.02
604	1.03725	5.440	0.05000	25.0	-0.108	-2.40
634	1.08725	5.310	0.05000	25.0	-0.130	-3.02
664	1.13725	5.144	0.05000	25.0	-0.166	-4.04
695	1.17075	5.004	0.03350	25.0	-0.140	-4.99
725	1.19625	4.868	0.02550	25.0	-0.136	-5.86
755	1.21625	4.740	0.02000	25.0	-0.128	-6.56
785	1.23400	4.607	0.01775	25.0	-0.133	-7.28
816	1.24925	4.481	0.01525	25.0	-0.126	-7.90
	1.249835 ⁽⁹⁾	4.476 ⁽⁹⁾				
846	1.26425	4.355	0.01500	25.0	-0.126	-7.66
876	1.28000	4.228	0.01575	25.0	-0.127	-6.85
906	1.29775	4.106	0.01775	25.0	-0.122	-6.11
937	1.32050	3.976	0.02275	25.0	-0.130	-5.25
967	1.34950	3.845	0.02900	25.0	-0.131	-4.26
997	1.38750	3.715	0.03800	25.0	-0.130	-3.25
1027	1.43750	3.585	0.05000	25.0	-0.130	-2.36
1058	1.48750	3.484	0.05000	25.0	-0.101	-1.82
1088	1.53750	3.402	0.05000	25.0	-0.082	-1.48
1118	1.58750	3.333	0.05000	25.0	-0.069	-1.26
1148	1.63750	3.274	0.05000	25.0	-0.059	-1.09
1178	1.68750	3.221	0.05000	25.0	-0.053	-0.98
1209	1.73750	3.175	0.05000	25.0	-0.046	-0.88
1239	1.78750	3.132	0.05000	25.0	-0.043	
1269	1.83750	3.094	0.05000	25.0	-0.038	
1300	1.88750	3.058	0.05000	25.0	-0.036	
1330	1.93750	3.026	0.05000	25.0	-0.032	
1360	1.98750	2.996	0.05000	25.0	-0.030	

EQP

Measure (normal) [1] Scope 1/1, Sample 4/30

Measured Values Scope 1/1, Sample 4/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	5.144	23.7	
2	5.200	23.9	0.056
3	5.376	24.1	0.176
4	5.540	24.3	0.164
5	5.682	24.5	0.142
6	5.801	24.7	0.119
7	5.908	24.8	0.107
8	5.993	24.9	0.085
9	6.079	25.0	0.086
10	6.158	25.0	0.079
11	6.231	25.1	0.073

Titration report

LabX

12	6.300	25.1	0.069
13	6.366	25.2	0.066
14	6.427	25.2	0.061
15	6.478	25.2	0.051
16	6.532	25.2	0.054
17	6.584	25.2	0.052
18	6.634	25.2	0.050
19	6.680	25.2	0.046
20	6.724	25.3	0.044
21	6.767	25.3	0.043
22	6.804	25.3	0.037
23	6.844	25.3	0.040
24	6.883	25.3	0.039
25	6.915	25.3	0.032
26	6.948	25.3	0.033
27	6.980	25.3	0.032
28	7.009	25.3	0.029
29	7.028	25.3	0.019
30	7.060	25.3	0.032
31	7.087	25.3	0.027
32	7.113	25.3	0.026
33	7.137	25.3	0.024
34	7.161	25.3	0.024
35	7.181	25.3	0.020
36	7.199	25.3	0.018
37	7.218	25.3	0.019
38	7.237	25.3	0.019
39	7.254	25.3	0.017
40	7.272	25.3	0.018
41	7.288	25.3	0.016
42	7.304	25.3	0.016
43	7.318	25.3	0.014
44	7.332	25.3	0.014
45	7.345	25.3	0.013
46	7.357	25.3	0.012
47	7.369	25.3	0.012
48	7.381	25.3	0.012
49	7.391	25.3	0.010
50	7.401	25.3	0.010
51	7.410	25.3	0.009
52	7.420	25.3	0.010
53	7.429	25.3	0.009
54	7.438	25.3	0.009
55	7.446	25.3	0.008
56	7.453	25.3	0.007
57	7.461	25.3	0.008
58	7.468	25.3	0.007
59	7.475	25.3	0.007
60	7.481	25.3	0.006


Titration (EQP) [1]

Scope 1/1, Sample 4/30


Titration report

LabX

E - V Curve Scope 1/1, Sample 4/30

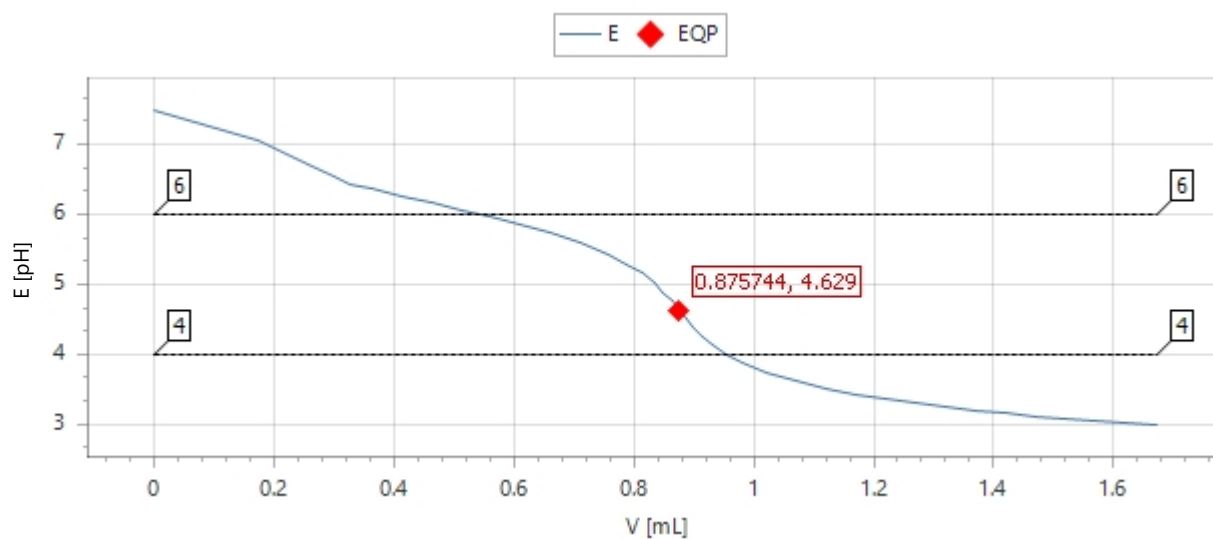

dE/dV - V Curve Scope 1/1, Sample 4/30

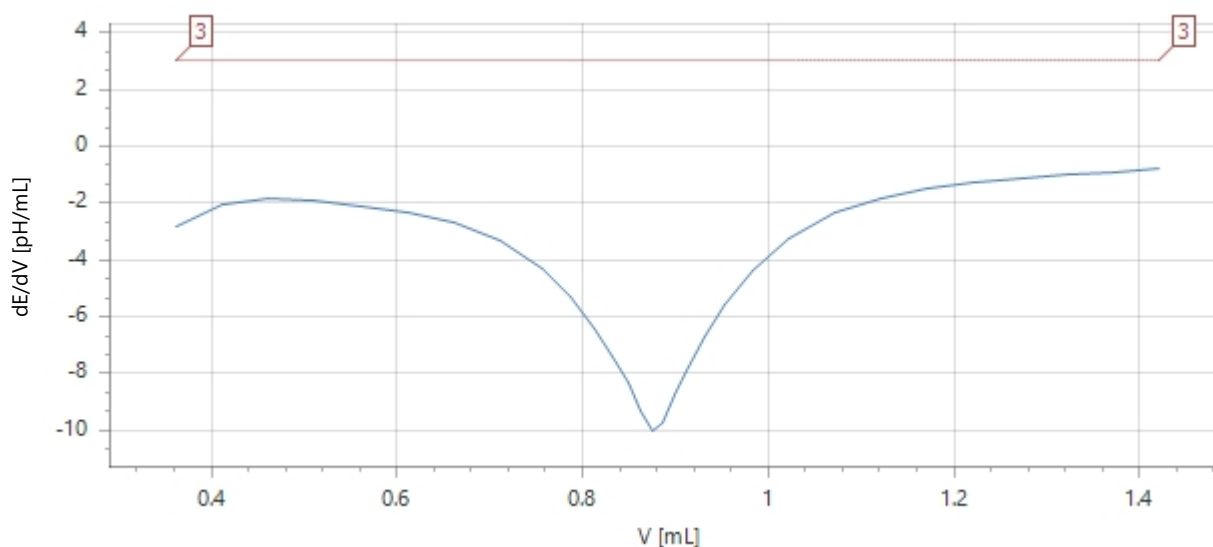

Measured Values Scope 1/1, Sample 4/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.495		25.3			
27	0.17150	7.057	0.17150	25.2	-0.438		
58	0.25725	6.724	0.08575	25.2	-0.333		
88	0.30000	6.533	0.04275	25.2	-0.191		
178	0.32625	6.423	0.02625	25.2	-0.110		
209	0.36175	6.360	0.03550	25.2	-0.063	-2.85	
239	0.41175	6.267	0.05000	25.2	-0.093	-2.07	
269	0.46175	6.171	0.05000	25.2	-0.096	-1.81	
300	0.51175	6.070	0.05000	25.2	-0.101	-1.92	
330	0.56175	5.966	0.05000	25.2	-0.104	-2.11	
360	0.61175	5.855	0.05000	25.2	-0.111	-2.32	
390	0.66175	5.735	0.05000	25.2	-0.120	-2.67	

Titration report

LabX

420	0.71175	5.591	0.05000	25.2	-0.144	-3.30	EQP
451	0.75650	5.434	0.04475	25.2	-0.157	-4.28	
481	0.78775	5.297	0.03125	25.2	-0.137	-5.30	
511	0.81275	5.161	0.02500	25.2	-0.136	-6.39	
542	0.83275	5.027	0.02000	25.2	-0.134	-7.42	
572	0.84900	4.895	0.01625	25.2	-0.132	-8.34	
602	0.86275	4.762	0.01375	25.2	-0.133	-9.33	
632	0.87475	4.640	0.01200	25.2	-0.122	-9.98	
	0.875744 ⁽⁹⁾	4.629 ⁽⁹⁾					
662	0.88725	4.504	0.01250	25.2	-0.136	-9.70	
693	0.89900	4.388	0.01175	25.2	-0.116	-8.74	
723	0.91400	4.258	0.01500	25.2	-0.130	-7.74	
753	0.93175	4.127	0.01775	25.2	-0.131	-6.71	
783	0.95350	4.000	0.02175	25.2	-0.127	-5.58	
814	0.98250	3.869	0.02900	25.2	-0.131	-4.38	
844	1.02150	3.736	0.03900	25.1	-0.133	-3.26	
874	1.07150	3.610	0.05000	25.1	-0.126	-2.36	
904	1.12150	3.506	0.05000	25.1	-0.104	-1.82	
935	1.17150	3.424	0.05000	25.1	-0.082	-1.50	
965	1.22150	3.354	0.05000	25.1	-0.070	-1.27	
995	1.27150	3.294	0.05000	25.1	-0.060	-1.11	
1025	1.32150	3.242	0.05000	25.1	-0.052	-0.98	
1056	1.37150	3.194	0.05000	25.1	-0.048	-0.89	
1086	1.42150	3.153	0.05000	25.1	-0.041	-0.80	
1116	1.47150	3.113	0.05000	25.1	-0.040		
1146	1.52150	3.078	0.05000	25.1	-0.035		
1176	1.57150	3.045	0.05000	25.1	-0.033		
1207	1.62150	3.015	0.05000	25.1	-0.030		
1237	1.67150	2.986	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 5/30

Measured Values Scope 1/1, Sample 5/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	5.489	23.7	
2	5.085	23.9	-0.404
3	5.248	24.1	0.163
4	5.422	24.3	0.174
5	5.571	24.5	0.149
6	5.698	24.7	0.127
7	5.807	24.8	0.109
8	5.897	24.9	0.090
9	5.989	25.0	0.092
10	6.074	25.1	0.085
11	6.152	25.1	0.078
12	6.227	25.2	0.075
13	6.296	25.2	0.069
14	6.361	25.2	0.065
15	6.415	25.2	0.054

Titration report

LabX


16	6.474	25.3	0.059
17	6.531	25.3	0.057
18	6.586	25.3	0.055
19	6.638	25.3	0.052
20	6.687	25.3	0.049
21	6.733	25.3	0.046
22	6.772	25.3	0.039
23	6.815	25.3	0.043
24	6.857	25.3	0.042
25	6.895	25.3	0.038
26	6.929	25.3	0.034
27	6.963	25.3	0.034
28	6.995	25.3	0.032
29	7.019	25.3	0.024
30	7.048	25.3	0.029
31	7.080	25.3	0.032
32	7.107	25.3	0.027
33	7.133	25.3	0.026
34	7.157	25.3	0.024
35	7.180	25.3	0.023
36	7.199	25.3	0.019
37	7.219	25.3	0.020
38	7.239	25.3	0.020
39	7.258	25.3	0.019
40	7.275	25.3	0.017
41	7.292	25.3	0.017
42	7.308	25.3	0.016
43	7.322	25.3	0.014
44	7.337	25.3	0.015
45	7.352	25.3	0.015
46	7.365	25.3	0.013
47	7.378	25.3	0.013
48	7.389	25.3	0.011
49	7.399	25.3	0.010
50	7.410	25.3	0.011
51	7.419	25.3	0.009
52	7.429	25.3	0.010
53	7.438	25.3	0.009
54	7.447	25.3	0.009
55	7.455	25.3	0.008
56	7.462	25.3	0.007
57	7.469	25.3	0.007
58	7.476	25.3	0.007
59	7.482	25.3	0.006
60	7.488	25.3	0.006

Titration (EQP) [1] Scope 1/1, Sample 5/30


E - V Curve Scope 1/1, Sample 5/30

Titration report

LabX

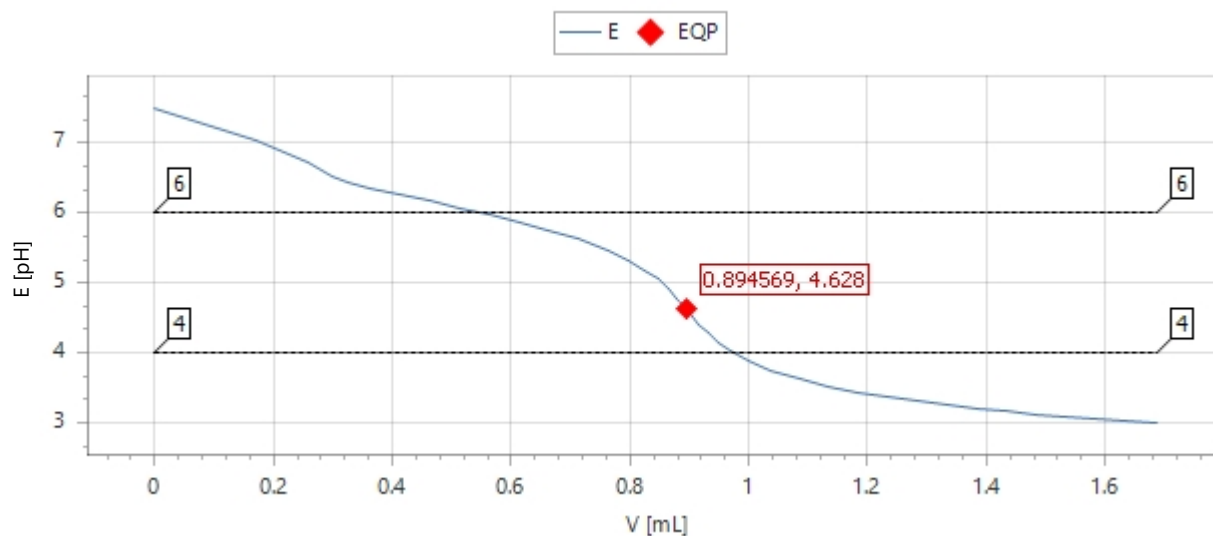

dE/dV - V Curve Scope 1/1, Sample 5/30

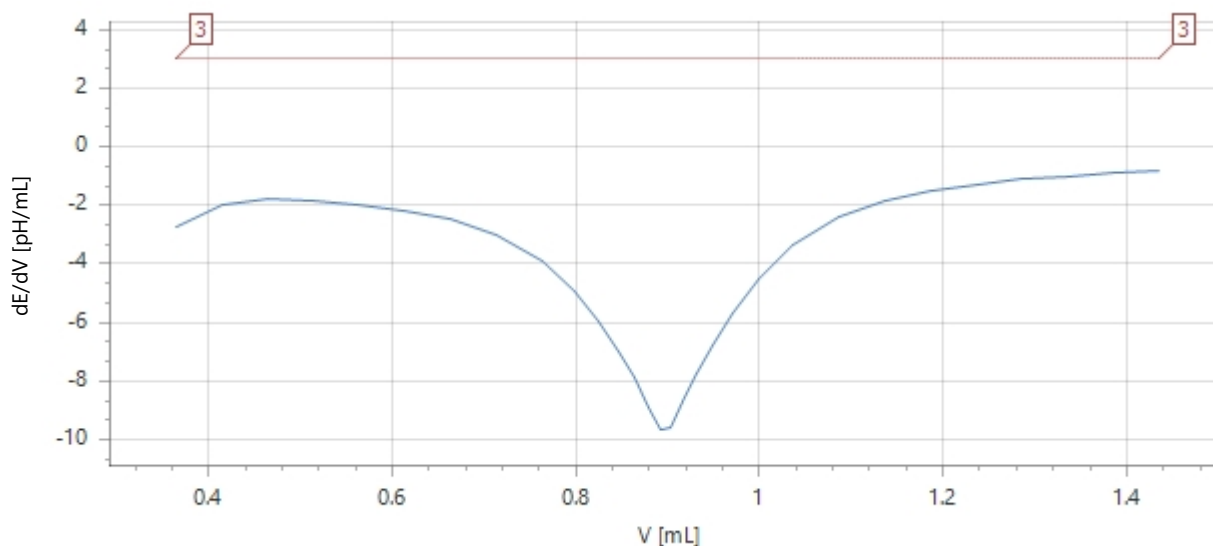

Measured Values Scope 1/1, Sample 5/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.501		25.3			
28	0.17150	7.040	0.17150	25.3	-0.461		
58	0.25725	6.717	0.08575	25.3	-0.323		
88	0.30000	6.528	0.04275	25.3	-0.189		
179	0.32650	6.423	0.02650	25.2	-0.105		
209	0.36375	6.356	0.03725	25.2	-0.067	-2.75	
239	0.41375	6.268	0.05000	25.2	-0.088	-2.01	
269	0.46375	6.167	0.05000	25.2	-0.101	-1.75	
300	0.51375	6.074	0.05000	25.2	-0.093	-1.85	
330	0.56375	5.977	0.05000	25.2	-0.097	-2.01	
360	0.61375	5.869	0.05000	25.2	-0.108	-2.18	
390	0.66375	5.755	0.05000	25.2	-0.114	-2.47	
421	0.71375	5.624	0.05000	25.1	-0.131	-3.01	
451	0.76375	5.460	0.05000	25.1	-0.164	-3.93	

Titration report

LabX

481	0.79800	5.322	0.03425	25.1	-0.138	-4.90	EQP
511	0.82525	5.186	0.02725	25.1	-0.136	-5.98	
542	0.84700	5.049	0.02175	25.1	-0.137	-7.01	
572	0.86425	4.918	0.01725	25.1	-0.131	-7.90	
602	0.87900	4.786	0.01475	25.1	-0.132	-8.87	
632	0.89200	4.655	0.01300	25.1	-0.131	-9.63	
	0.894569 ⁽⁹⁾	4.628 ⁽⁹⁾					
662	0.90400	4.531	0.01200	25.1	-0.124	-9.58	
693	0.91700	4.402	0.01300	25.1	-0.129	-8.71	
723	0.93125	4.278	0.01425	25.1	-0.124	-7.73	
753	0.94875	4.147	0.01750	25.1	-0.131	-6.77	
784	0.97025	4.018	0.02150	25.1	-0.129	-5.68	
814	0.99800	3.889	0.02775	25.1	-0.129	-4.52	
844	1.03575	3.756	0.03775	25.1	-0.133	-3.38	
874	1.08575	3.622	0.05000	25.1	-0.134	-2.42	
904	1.13575	3.519	0.05000	25.1	-0.103	-1.86	
935	1.18575	3.436	0.05000	25.1	-0.083	-1.51	
965	1.23575	3.366	0.05000	25.1	-0.070	-1.27	
995	1.28575	3.305	0.05000	25.1	-0.061	-1.11	
1025	1.33575	3.252	0.05000	25.1	-0.053	-0.99	
1056	1.38575	3.205	0.05000	25.1	-0.047	-0.89	
1086	1.43575	3.163	0.05000	25.1	-0.042	-0.81	
1116	1.48575	3.123	0.05000	25.1	-0.040		
1146	1.53575	3.088	0.05000	25.1	-0.035		
1177	1.58575	3.055	0.05000	25.1	-0.033		
1207	1.63575	3.024	0.05000	25.1	-0.031		
1237	1.68575	2.995	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 6/30

Measured Values Scope 1/1, Sample 6/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.295	23.7	
2	4.886	23.8	0.591
3	5.162	24.1	0.276
4	5.359	24.3	0.197
5	5.513	24.5	0.154
6	5.641	24.6	0.128
7	5.750	24.8	0.109
8	5.838	24.9	0.088
9	5.925	24.9	0.087
10	6.005	25.0	0.080
11	6.077	25.0	0.072
12	6.146	25.1	0.069
13	6.211	25.1	0.065
14	6.271	25.1	0.060
15	6.323	25.1	0.052
16	6.376	25.2	0.053
17	6.426	25.2	0.050

Titration report

LabX


18	6.472	25.2	0.046
19	6.518	25.2	0.046
20	6.562	25.2	0.044
21	6.603	25.2	0.041
22	6.638	25.2	0.035
23	6.676	25.2	0.038
24	6.712	25.2	0.036
25	6.745	25.2	0.033
26	6.777	25.2	0.032
27	6.808	25.2	0.031
28	6.839	25.2	0.031
29	6.868	25.2	0.029
30	6.895	25.2	0.027
31	6.918	25.2	0.023
32	6.943	25.2	0.025
33	6.967	25.2	0.024
34	6.990	25.2	0.023
35	7.011	25.2	0.021
36	7.023	25.2	0.012
37	7.045	25.2	0.022
38	7.068	25.2	0.023
39	7.086	25.2	0.018
40	7.105	25.2	0.019
41	7.123	25.2	0.018
42	7.140	25.2	0.017
43	7.154	25.2	0.014
44	7.168	25.2	0.014
45	7.182	25.2	0.014
46	7.196	25.2	0.014
47	7.209	25.2	0.013
48	7.220	25.2	0.011
49	7.232	25.2	0.012
50	7.243	25.2	0.011
51	7.254	25.2	0.011
52	7.264	25.2	0.010
53	7.274	25.2	0.010
54	7.283	25.2	0.009
55	7.292	25.2	0.009
56	7.300	25.2	0.008
57	7.308	25.2	0.008
58	7.316	25.2	0.008
59	7.324	25.2	0.008
60	7.331	25.2	0.007

Titration (EQP) [1] Scope 1/1, Sample 6/30


E - V Curve Scope 1/1, Sample 6/30

Titration report

LabX

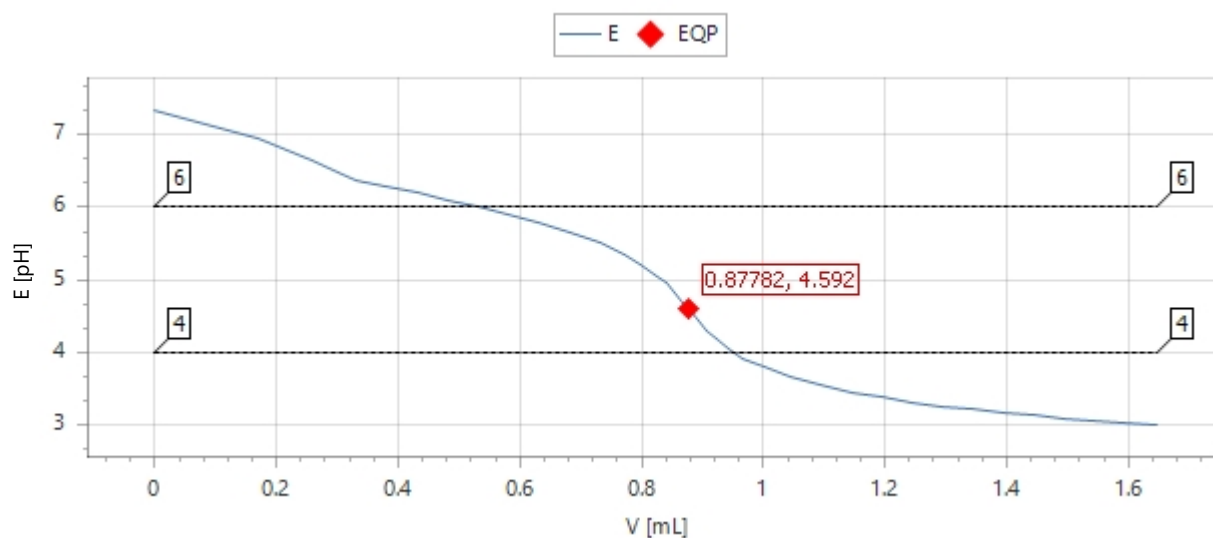

dE/dV - V Curve Scope 1/1, Sample 6/30

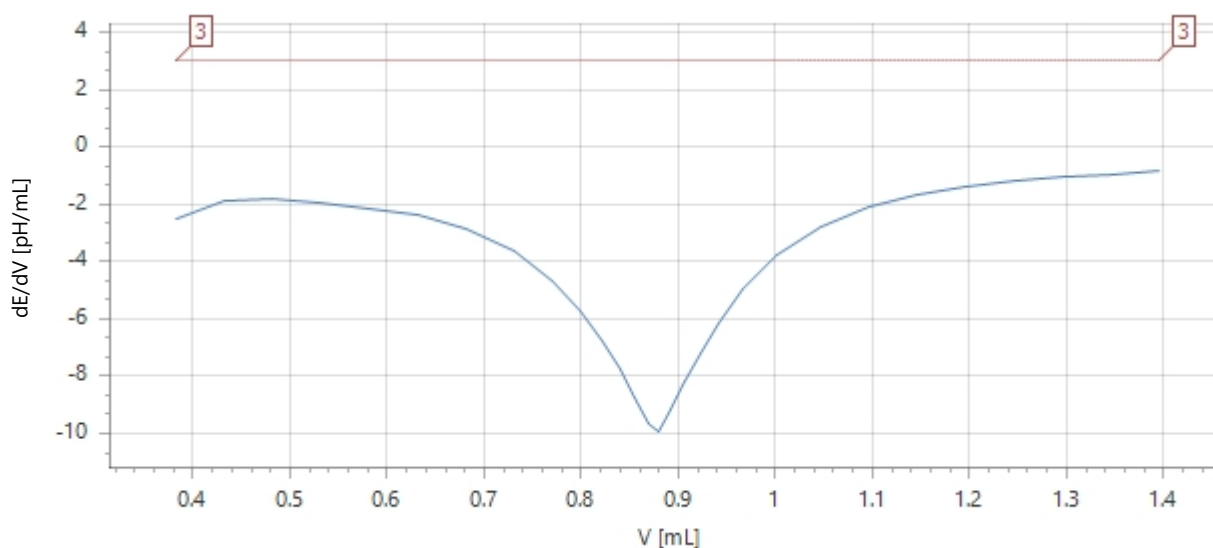

Measured Values Scope 1/1, Sample 6/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.347		25.2			
28	0.17150	6.955	0.17150	25.2	-0.392		
58	0.25725	6.639	0.08575	25.2	-0.316		
88	0.30000	6.470	0.04275	25.2	-0.169		
179	0.33200	6.379	0.03200	25.2	-0.091		
209	0.38200	6.292	0.05000	25.2	-0.087	-2.51	
239	0.43200	6.196	0.05000	25.2	-0.096	-1.89	
270	0.48200	6.101	0.05000	25.2	-0.095	-1.80	
300	0.53200	6.002	0.05000	25.2	-0.099	-1.95	
330	0.58200	5.895	0.05000	25.2	-0.107	-2.15	
360	0.63200	5.786	0.05000	25.2	-0.109	-2.39	
390	0.68200	5.659	0.05000	25.2	-0.127	-2.85	
421	0.73200	5.505	0.05000	25.2	-0.154	-3.65	
451	0.77050	5.358	0.03850	25.2	-0.147	-4.66	

Titration report

LabX

481	0.79950	5.219	0.02900	25.2	-0.139	-5.76	EQP
511	0.82200	5.086	0.02250	25.2	-0.133	-6.81	
542	0.84050	4.952	0.01850	25.2	-0.134	-7.79	
572	0.85575	4.816	0.01525	25.2	-0.136	-8.73	
602	0.86825	4.692	0.01250	25.2	-0.124	-9.63	
	0.877820	4.592					
632	0.88050	4.564	0.01225	25.2	-0.128	-9.90	
662	0.89275	4.436	0.01225	25.2	-0.128	-9.19	
693	0.90600	4.312	0.01325	25.2	-0.124	-8.19	
723	0.92225	4.183	0.01625	25.2	-0.129	-7.22	
753	0.94225	4.052	0.02000	25.2	-0.131	-6.13	
784	0.96725	3.923	0.02500	25.2	-0.129	-4.96	
814	1.00100	3.790	0.03375	25.2	-0.133	-3.78	
844	1.04575	3.659	0.04475	25.2	-0.131	-2.75	
874	1.09575	3.544	0.05000	25.2	-0.115	-2.05	
904	1.14575	3.455	0.05000	25.2	-0.089	-1.63	
935	1.19575	3.380	0.05000	25.2	-0.075	-1.36	
965	1.24575	3.315	0.05000	25.2	-0.065	-1.17	
995	1.29575	3.260	0.05000	25.2	-0.055	-1.04	
1026	1.34575	3.210	0.05000	25.2	-0.050	-0.93	
1056	1.39575	3.166	0.05000	25.2	-0.044	-0.84	
1086	1.44575	3.125	0.05000	25.2	-0.041		
1116	1.49575	3.089	0.05000	25.2	-0.036		
1146	1.54575	3.054	0.05000	25.1	-0.035		
1177	1.59575	3.023	0.05000	25.1	-0.031		
1207	1.64575	2.993	0.05000	25.1	-0.030		

Measure (normal) [1] Scope 1/1, Sample 7/30

Measured Values Scope 1/1, Sample 7/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.398	23.6	
2	4.919	23.7	0.521
3	5.179	24.0	0.260
4	5.372	24.2	0.193
5	5.528	24.4	0.156
6	5.657	24.6	0.129
7	5.768	24.7	0.111
8	5.856	24.8	0.088
9	5.947	24.9	0.091
10	6.029	25.0	0.082
11	6.105	25.0	0.076
12	6.175	25.1	0.070
13	6.241	25.1	0.066
14	6.303	25.1	0.062
15	6.355	25.1	0.052
16	6.410	25.1	0.055
17	6.462	25.1	0.052
18	6.512	25.2	0.050

Titration report

LabX


19	6.558	25.2	0.046
20	6.602	25.2	0.044
21	6.644	25.2	0.042
22	6.680	25.2	0.036
23	6.718	25.2	0.038
24	6.755	25.2	0.037
25	6.790	25.2	0.035
26	6.824	25.2	0.034
27	6.857	25.2	0.033
28	6.889	25.2	0.032
29	6.912	25.2	0.023
30	6.938	25.2	0.026
31	6.964	25.2	0.026
32	6.990	25.2	0.026
33	7.012	25.2	0.022
34	7.028	25.2	0.016
35	7.056	25.2	0.028
36	7.077	25.2	0.021
37	7.097	25.2	0.020
38	7.116	25.2	0.019
39	7.135	25.2	0.019
40	7.152	25.2	0.017
41	7.169	25.2	0.017
42	7.185	25.2	0.016
43	7.199	25.2	0.014
44	7.213	25.2	0.014
45	7.227	25.2	0.014
46	7.240	25.2	0.013
47	7.253	25.2	0.013
48	7.265	25.2	0.012
49	7.277	25.2	0.012
50	7.287	25.2	0.010
51	7.298	25.2	0.011
52	7.308	25.2	0.010
53	7.318	25.2	0.010
54	7.327	25.2	0.009
55	7.335	25.2	0.008
56	7.343	25.2	0.008
57	7.350	25.2	0.007
58	7.358	25.2	0.008
59	7.366	25.2	0.008
60	7.373	25.2	0.007

Titration (EQP) [1] Scope 1/1, Sample 7/30


E - V Curve Scope 1/1, Sample 7/30

Titration report

LabX

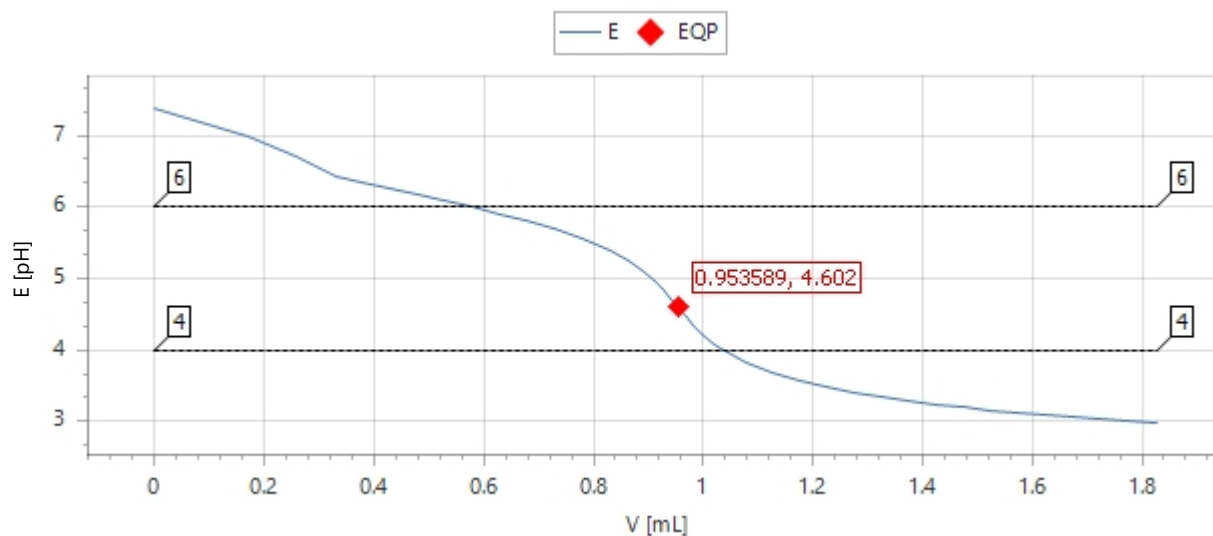

dE/dV - V Curve Scope 1/1, Sample 7/30

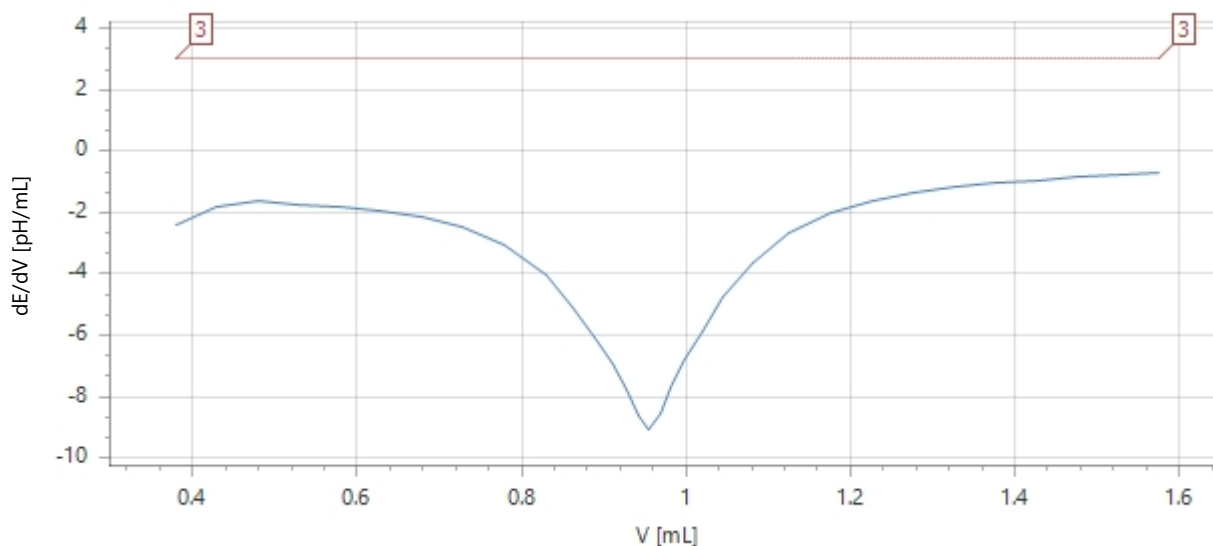

Measured Values Scope 1/1, Sample 7/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.387		25.2			
28	0.17150	6.995	0.17150	25.2	-0.392		
58	0.25725	6.698	0.08575	25.2	-0.297		
88	0.30000	6.527	0.04275	25.2	-0.171		
179	0.32925	6.433	0.02925	25.2	-0.094		
209	0.37925	6.354	0.05000	25.1	-0.079	-2.43	
239	0.42925	6.266	0.05000	25.1	-0.088	-1.79	
269	0.47925	6.173	0.05000	25.1	-0.093	-1.62	
300	0.52925	6.089	0.05000	25.1	-0.084	-1.73	
330	0.57925	5.997	0.05000	25.1	-0.092	-1.84	
360	0.62925	5.902	0.05000	25.1	-0.095	-1.93	
390	0.67925	5.801	0.05000	25.1	-0.101	-2.13	
421	0.72925	5.687	0.05000	25.1	-0.114	-2.48	
451	0.77925	5.553	0.05000	25.1	-0.134	-3.05	

Titration report

LabX

481	0.82925	5.388	0.05000	25.1	-0.165	-4.05	EQP
511	0.86300	5.247	0.03375	25.1	-0.141	-5.05	
542	0.88875	5.112	0.02575	25.1	-0.135	-6.03	
572	0.90950	4.978	0.02075	25.1	-0.134	-6.92	
602	0.92650	4.848	0.01700	25.1	-0.130	-7.73	
632	0.94150	4.716	0.01500	25.1	-0.132	-8.61	
	0.953589 ⁽⁹⁾	4.602 ⁽⁹⁾					
663	0.95475	4.591	0.01325	25.1	-0.125	-9.04	
693	0.96825	4.461	0.01350	25.1	-0.130	-8.56	
723	0.98200	4.341	0.01375	25.1	-0.120	-7.68	
753	0.99950	4.207	0.01750	25.1	-0.134	-6.77	
784	1.01950	4.081	0.02000	25.1	-0.126	-5.84	
814	1.04575	3.952	0.02625	25.1	-0.129	-4.75	
844	1.08075	3.817	0.03500	25.1	-0.135	-3.64	
874	1.12600	3.687	0.04525	25.1	-0.130	-2.67	
904	1.17600	3.576	0.05000	25.1	-0.111	-2.01	
935	1.22600	3.487	0.05000	25.1	-0.089	-1.61	
965	1.27600	3.413	0.05000	25.1	-0.074	-1.35	
995	1.32600	3.349	0.05000	25.1	-0.064	-1.17	
1025	1.37600	3.293	0.05000	25.1	-0.056	-1.03	
1056	1.42600	3.244	0.05000	25.1	-0.049	-0.93	
1086	1.47600	3.200	0.05000	25.1	-0.044	-0.84	
1116	1.52600	3.160	0.05000	25.1	-0.040	-0.76	
1146	1.57600	3.123	0.05000	25.1	-0.037	-0.70	
1177	1.62600	3.089	0.05000	25.1	-0.034		
1207	1.67600	3.058	0.05000	25.1	-0.031		
1237	1.72600	3.028	0.05000	25.1	-0.030		
1267	1.77600	3.002	0.05000	25.1	-0.026		
1298	1.82600	2.975	0.05000	25.1	-0.027		

Measure (normal) [1] Scope 1/1, Sample 8/30

Measured Values Scope 1/1, Sample 8/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.350	22.2	
2	5.079	22.5	0.729
3	5.322	22.9	0.243
4	5.493	23.3	0.171
5	5.629	23.7	0.136
6	5.742	24.0	0.113
7	5.831	24.2	0.089
8	5.920	24.4	0.089
9	5.999	24.6	0.079
10	6.073	24.7	0.074
11	6.139	24.8	0.066
12	6.200	24.9	0.061
13	6.258	24.9	0.058
14	6.307	25.0	0.049
15	6.357	25.0	0.050

Titration report

LabX


16	6.405	25.1	0.048
17	6.452	25.1	0.047
18	6.496	25.1	0.044
19	6.539	25.1	0.043
20	6.575	25.1	0.036
21	6.614	25.1	0.039
22	6.651	25.1	0.037
23	6.687	25.1	0.036
24	6.721	25.1	0.034
25	6.753	25.1	0.032
26	6.780	25.1	0.027
27	6.811	25.1	0.031
28	6.840	25.1	0.029
29	6.868	25.1	0.028
30	6.893	25.1	0.025
31	6.915	25.1	0.022
32	6.938	25.1	0.023
33	6.958	25.1	0.020
34	6.980	25.1	0.022
35	7.000	25.1	0.020
36	7.016	25.1	0.016
37	7.030	25.1	0.014
38	7.054	25.1	0.024
39	7.075	25.1	0.021
40	7.090	25.1	0.015
41	7.107	25.1	0.017
42	7.123	25.1	0.016
43	7.139	25.1	0.016
44	7.153	25.1	0.014
45	7.167	25.1	0.014
46	7.179	25.1	0.012
47	7.191	25.1	0.012
48	7.203	25.1	0.012
49	7.214	25.1	0.011
50	7.225	25.1	0.011
51	7.236	25.1	0.011
52	7.246	25.1	0.010
53	7.254	25.1	0.008
54	7.264	25.1	0.010
55	7.273	25.1	0.009
56	7.281	25.1	0.008
57	7.289	25.1	0.008
58	7.297	25.1	0.008
59	7.304	25.1	0.007
60	7.311	25.1	0.007

Titration (EQP) [1] Scope 1/1, Sample 8/30


E - V Curve Scope 1/1, Sample 8/30

Titration report

LabX

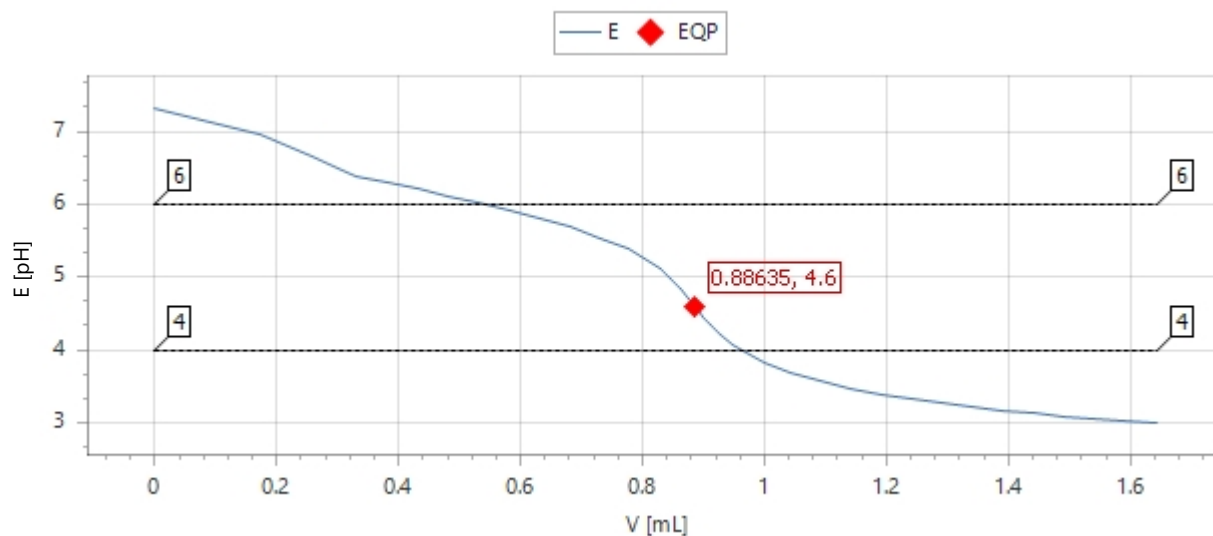

dE/dV - V Curve Scope 1/1, Sample 8/30

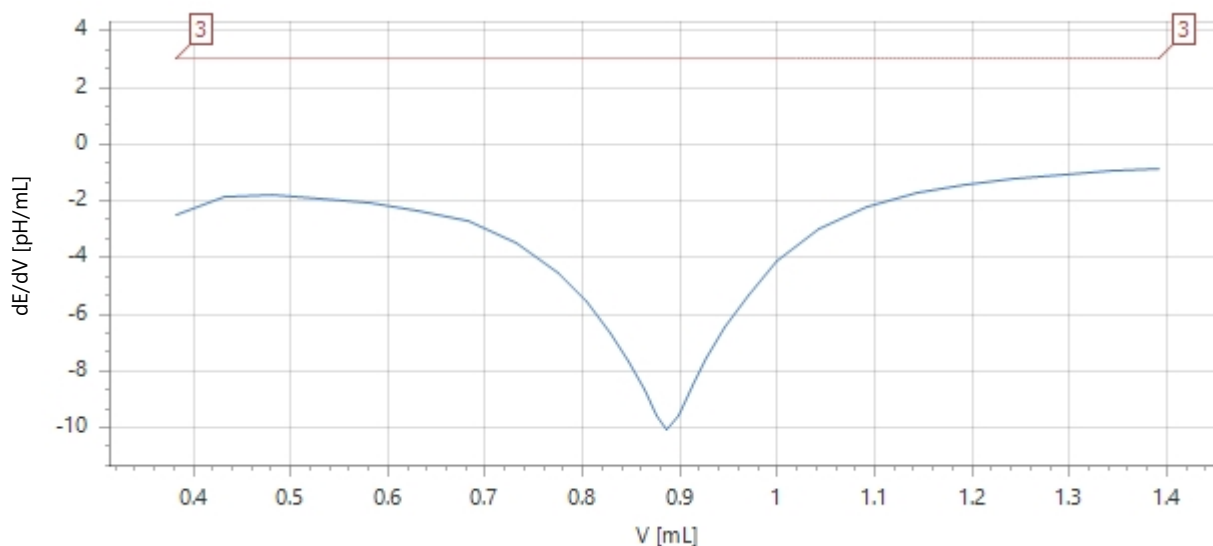

Measured Values Scope 1/1, Sample 8/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.326		25.1			
28	0.17150	6.948	0.17150	25.1	-0.378		
58	0.25725	6.650	0.08575	25.1	-0.298		
88	0.30000	6.486	0.04275	25.1	-0.164		
179	0.33200	6.395	0.03200	25.2	-0.091		
209	0.38200	6.307	0.05000	25.2	-0.088	-2.46	
239	0.43200	6.213	0.05000	25.2	-0.094	-1.88	
269	0.48200	6.117	0.05000	25.2	-0.096	-1.77	
300	0.53200	6.023	0.05000	25.2	-0.094	-1.91	
330	0.58200	5.920	0.05000	25.2	-0.103	-2.10	
360	0.63200	5.811	0.05000	25.2	-0.109	-2.32	
390	0.68200	5.686	0.05000	25.2	-0.125	-2.74	
421	0.73200	5.540	0.05000	25.2	-0.146	-3.47	
451	0.77550	5.379	0.04350	25.2	-0.161	-4.53	

Titration report

LabX

481	0.80425	5.249	0.02875	25.2	-0.130	-5.56	EQP
511	0.82850	5.108	0.02425	25.2	-0.141	-6.70	
542	0.84675	4.980	0.01825	25.2	-0.128	-7.66	
572	0.86300	4.841	0.01625	25.2	-0.139	-8.62	
602	0.87575	4.712	0.01275	25.2	-0.129	-9.60	
	0.886350	4.600					
632	0.88725	4.591	0.01150	25.2	-0.121	-10.04	
663	0.89950	4.463	0.01225	25.2	-0.128	-9.54	
693	0.91275	4.334	0.01325	25.2	-0.129	-8.50	
723	0.92750	4.210	0.01475	25.2	-0.124	-7.55	
753	0.94625	4.078	0.01875	25.2	-0.132	-6.48	
784	0.96950	3.949	0.02325	25.2	-0.129	-5.31	
814	1.00025	3.819	0.03075	25.2	-0.130	-4.13	
844	1.04250	3.684	0.04225	25.2	-0.135	-3.02	
874	1.09250	3.562	0.05000	25.2	-0.122	-2.21	
904	1.14250	3.466	0.05000	25.2	-0.096	-1.73	
935	1.19250	3.387	0.05000	25.2	-0.079	-1.43	
965	1.24250	3.320	0.05000	25.2	-0.067	-1.22	
995	1.29250	3.262	0.05000	25.2	-0.058	-1.07	
1026	1.34250	3.211	0.05000	25.2	-0.051	-0.96	
1056	1.39250	3.165	0.05000	25.2	-0.046	-0.86	
1086	1.44250	3.124	0.05000	25.2	-0.041		
1116	1.49250	3.086	0.05000	25.2	-0.038		
1146	1.54250	3.051	0.05000	25.2	-0.035		
1177	1.59250	3.019	0.05000	25.1	-0.032		
1207	1.64250	2.989	0.05000	25.1	-0.030		

Measure (normal) [1] Scope 1/1, Sample 9/30

Measured Values Scope 1/1, Sample 9/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.578	23.7	
2	4.965	23.9	0.387
3	5.209	24.1	0.244
4	5.397	24.3	0.188
5	5.548	24.5	0.151
6	5.675	24.7	0.127
7	5.784	24.8	0.109
8	5.870	24.9	0.086
9	5.958	25.0	0.088
10	6.038	25.1	0.080
11	6.111	25.1	0.073
12	6.178	25.1	0.067
13	6.242	25.2	0.064
14	6.300	25.2	0.058
15	6.351	25.2	0.051
16	6.404	25.2	0.053
17	6.454	25.2	0.050
18	6.502	25.2	0.048

Titration report

LabX


19	6.548	25.3	0.046
20	6.591	25.3	0.043
21	6.633	25.3	0.042
22	6.672	25.3	0.039
23	6.706	25.3	0.034
24	6.741	25.3	0.035
25	6.776	25.3	0.035
26	6.810	25.3	0.034
27	6.842	25.3	0.032
28	6.874	25.3	0.032
29	6.899	25.3	0.025
30	6.925	25.3	0.026
31	6.952	25.3	0.027
32	6.977	25.3	0.025
33	7.000	25.3	0.023
34	7.019	25.3	0.019
35	7.038	25.3	0.019
36	7.061	25.3	0.023
37	7.083	25.3	0.022
38	7.103	25.3	0.020
39	7.122	25.3	0.019
40	7.138	25.3	0.016
41	7.155	25.3	0.017
42	7.170	25.3	0.015
43	7.184	25.3	0.014
44	7.198	25.3	0.014
45	7.212	25.3	0.014
46	7.224	25.3	0.012
47	7.237	25.3	0.013
48	7.249	25.3	0.012
49	7.260	25.3	0.011
50	7.271	25.3	0.011
51	7.281	25.3	0.010
52	7.290	25.3	0.009
53	7.300	25.3	0.010
54	7.310	25.3	0.010
55	7.319	25.3	0.009
56	7.327	25.3	0.008
57	7.335	25.3	0.008
58	7.343	25.3	0.008
59	7.350	25.3	0.007
60	7.357	25.3	0.007

Titration (EQP) [1] Scope 1/1, Sample 9/30


E - V Curve Scope 1/1, Sample 9/30

Titration report

LabX

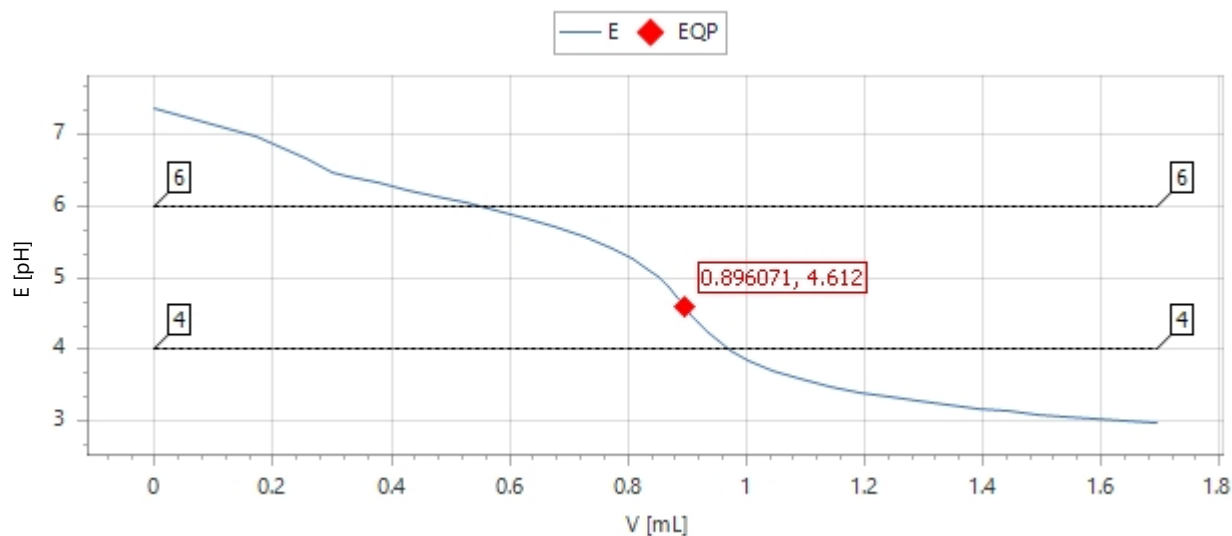

dE/dV - V Curve Scope 1/1, Sample 9/30

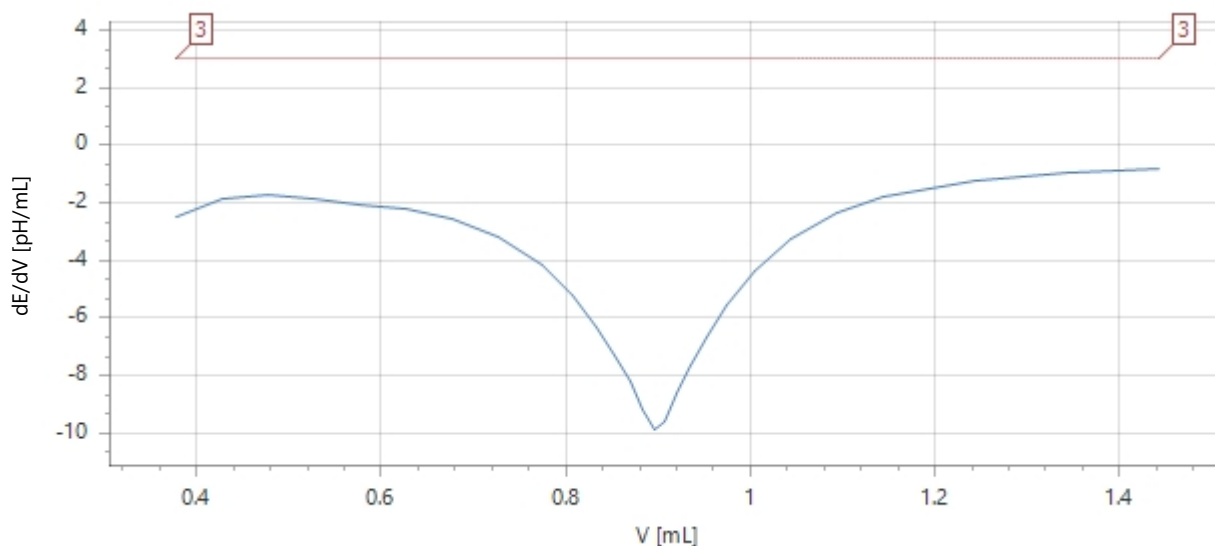

Measured Values Scope 1/1, Sample 9/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.374		25.3			
28	0.17150	6.969	0.17150	25.2	-0.405		
58	0.25725	6.663	0.08575	25.2	-0.306		
88	0.30000	6.485	0.04275	25.2	-0.178		
179	0.32725	6.407	0.02725	25.2	-0.078		
209	0.37725	6.323	0.05000	25.2	-0.084	-2.47	
239	0.42725	6.232	0.05000	25.2	-0.091	-1.83	
269	0.47725	6.141	0.05000	25.2	-0.091	-1.72	
300	0.52725	6.044	0.05000	25.2	-0.097	-1.88	
330	0.57725	5.945	0.05000	25.2	-0.099	-2.04	
360	0.62725	5.838	0.05000	25.2	-0.107	-2.23	
390	0.67725	5.721	0.05000	25.2	-0.117	-2.58	
420	0.72725	5.582	0.05000	25.2	-0.139	-3.18	
451	0.77475	5.419	0.04750	25.2	-0.163	-4.19	

Titration report

LabX

481	0.80700	5.282	0.03225	25.2	-0.137	-5.22	EQP
511	0.83250	5.145	0.02550	25.2	-0.137	-6.31	
542	0.85250	5.013	0.02000	25.2	-0.132	-7.30	
572	0.86925	4.877	0.01675	25.2	-0.136	-8.22	
602	0.88300	4.750	0.01375	25.2	-0.127	-9.19	
632	0.89575	4.615	0.01275	25.2	-0.135	-9.85	
	0.896071	4.612					
662	0.90700	4.500	0.01125	25.2	-0.115	-9.61	
693	0.92075	4.366	0.01375	25.2	-0.134	-8.62	
723	0.93525	4.238	0.01450	25.2	-0.128	-7.65	
753	0.95225	4.115	0.01700	25.2	-0.123	-6.70	
783	0.97475	3.984	0.02250	25.2	-0.131	-5.55	
814	1.00400	3.853	0.02925	25.2	-0.131	-4.36	
844	1.04325	3.718	0.03925	25.1	-0.135	-3.24	
874	1.09325	3.590	0.05000	25.1	-0.128	-2.34	
904	1.14325	3.490	0.05000	25.1	-0.100	-1.81	
935	1.19325	3.408	0.05000	25.1	-0.082	-1.48	
965	1.24325	3.339	0.05000	25.1	-0.069	-1.26	
995	1.29325	3.279	0.05000	25.1	-0.060	-1.10	
1025	1.34325	3.227	0.05000	25.1	-0.052	-0.98	
1056	1.39325	3.180	0.05000	25.1	-0.047	-0.88	
1086	1.44325	3.138	0.05000	25.1	-0.042	-0.80	
1116	1.49325	3.099	0.05000	25.1	-0.039		
1146	1.54325	3.064	0.05000	25.1	-0.035		
1177	1.59325	3.031	0.05000	25.1	-0.033		
1207	1.64325	3.001	0.05000	25.1	-0.030		
1237	1.69325	2.972	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 10/30

Measured Values Scope 1/1, Sample 10/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.525	23.6	
2	5.069	23.8	0.544
3	5.326	24.0	0.257
4	5.512	24.2	0.186
5	5.661	24.4	0.149
6	5.786	24.6	0.125
7	5.894	24.7	0.108
8	5.982	24.8	0.088
9	6.070	24.9	0.088
10	6.150	25.0	0.080
11	6.226	25.0	0.076
12	6.297	25.1	0.071
13	6.362	25.1	0.065
14	6.423	25.1	0.061
15	6.474	25.1	0.051
16	6.529	25.2	0.055
17	6.581	25.2	0.052

Titration report

LabX


18	6.631	25.2	0.050
19	6.677	25.2	0.046
20	6.720	25.2	0.043
21	6.762	25.2	0.042
22	6.803	25.2	0.041
23	6.838	25.2	0.035
24	6.875	25.2	0.037
25	6.907	25.2	0.032
26	6.938	25.2	0.031
27	6.970	25.2	0.032
28	7.000	25.2	0.030
29	7.019	25.2	0.019
30	7.047	25.2	0.028
31	7.077	25.2	0.030
32	7.102	25.2	0.025
33	7.126	25.2	0.024
34	7.147	25.2	0.021
35	7.168	25.2	0.021
36	7.186	25.2	0.018
37	7.205	25.2	0.019
38	7.224	25.2	0.019
39	7.241	25.2	0.017
40	7.257	25.2	0.016
41	7.272	25.2	0.015
42	7.287	25.2	0.015
43	7.300	25.2	0.013
44	7.313	25.2	0.013
45	7.327	25.2	0.014
46	7.339	25.2	0.012
47	7.351	25.2	0.012
48	7.362	25.2	0.011
49	7.372	25.2	0.010
50	7.381	25.2	0.009
51	7.391	25.2	0.010
52	7.400	25.2	0.009
53	7.408	25.2	0.008
54	7.417	25.2	0.009
55	7.425	25.2	0.008
56	7.433	25.2	0.008
57	7.439	25.2	0.006
58	7.446	25.2	0.007
59	7.452	25.2	0.006
60	7.459	25.2	0.007

Titration (EQP) [1] Scope 1/1, Sample 10/30


E - V Curve Scope 1/1, Sample 10/30

Titration report

LabX

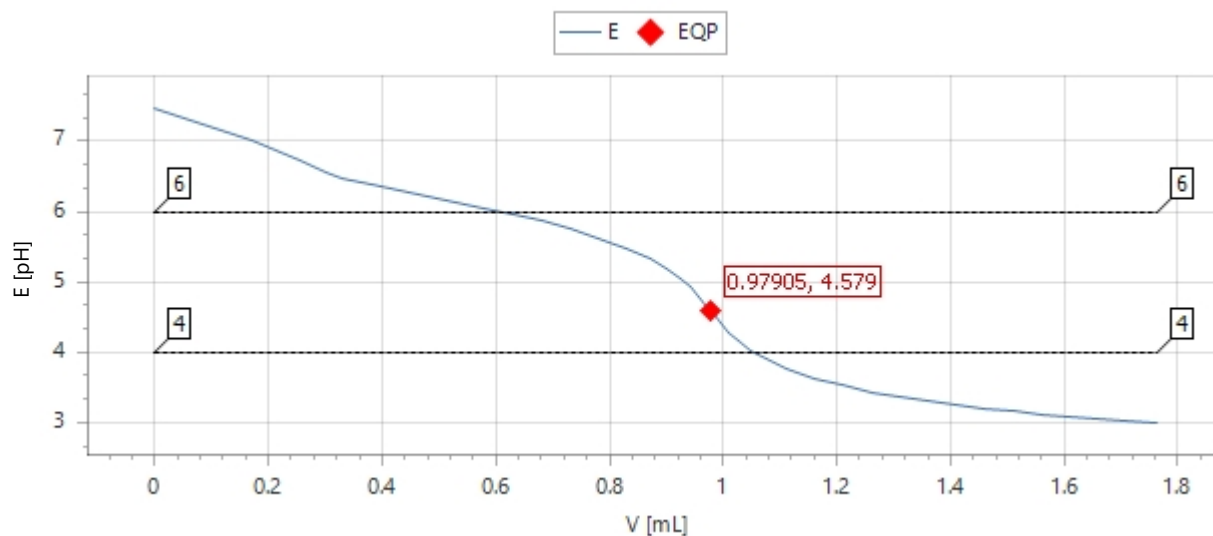

dE/dV - V Curve Scope 1/1, Sample 10/30

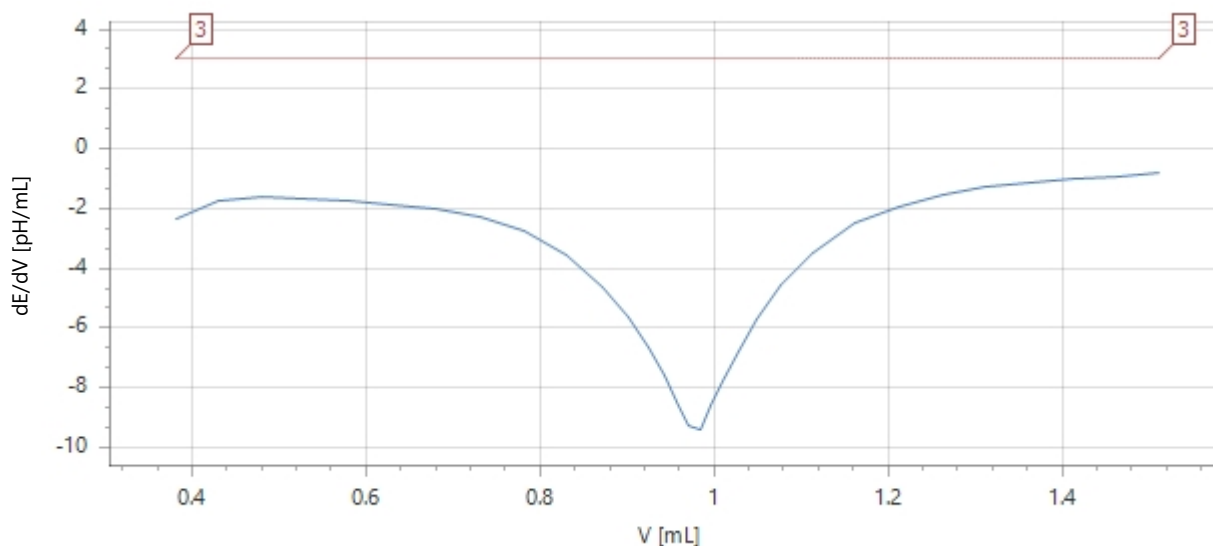

Measured Values Scope 1/1, Sample 10/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.473		25.2			
29	0.17150	7.023	0.17150	25.2	-0.450		
59	0.25725	6.731	0.08575	25.2	-0.292		
90	0.30000	6.566	0.04275	25.2	-0.165		
181	0.33125	6.475	0.03125	25.2	-0.091		
211	0.38125	6.396	0.05000	25.2	-0.079	-2.32	
241	0.43125	6.310	0.05000	25.2	-0.086	-1.74	
271	0.48125	6.221	0.05000	25.2	-0.089	-1.59	
302	0.53125	6.136	0.05000	25.2	-0.085	-1.68	
332	0.58125	6.048	0.05000	25.1	-0.088	-1.77	
362	0.63125	5.959	0.05000	25.1	-0.089	-1.85	
392	0.68125	5.862	0.05000	25.1	-0.097	-2.02	
423	0.73125	5.754	0.05000	25.1	-0.108	-2.31	
453	0.78125	5.628	0.05000	25.1	-0.126	-2.78	

Titration report

LabX

483	0.83125	5.480	0.05000	25.1	-0.148	-3.58	EQP
513	0.87300	5.321	0.04175	25.1	-0.159	-4.65	
544	0.90075	5.190	0.02775	25.1	-0.131	-5.65	
574	0.92375	5.054	0.02300	25.1	-0.136	-6.67	
604	0.94225	4.920	0.01850	25.1	-0.134	-7.57	
634	0.95750	4.790	0.01525	25.1	-0.130	-8.43	
664	0.97100	4.660	0.01350	25.1	-0.130	-9.28	
	0.979050 ⁽⁹⁾	4.579 ⁽⁹⁾					
695	0.98350	4.534	0.01250	25.1	-0.126	-9.37	
725	0.99625	4.413	0.01275	25.1	-0.121	-8.62	
755	1.01125	4.280	0.01500	25.1	-0.133	-7.64	
786	1.02775	4.157	0.01650	25.1	-0.123	-6.77	
816	1.04925	4.027	0.02150	25.1	-0.130	-5.71	
846	1.07700	3.895	0.02775	25.1	-0.132	-4.57	
876	1.11300	3.764	0.03600	25.1	-0.131	-3.48	
906	1.16200	3.629	0.04900	25.1	-0.135	-2.51	
937	1.21200	3.524	0.05000	25.1	-0.105	-1.91	
967	1.26200	3.438	0.05000	25.1	-0.086	-1.54	
997	1.31200	3.366	0.05000	25.1	-0.072	-1.30	
1027	1.36200	3.304	0.05000	25.1	-0.062	-1.14	
1058	1.41200	3.251	0.05000	25.1	-0.053	-1.01	
1088	1.46200	3.202	0.05000	25.1	-0.049	-0.91	
1118	1.51200	3.158	0.05000	25.1	-0.044	-0.82	
1148	1.56200	3.119	0.05000	25.1	-0.039		
1179	1.61200	3.083	0.05000	25.1	-0.036		
1209	1.66200	3.050	0.05000	25.1	-0.033		
1239	1.71200	3.019	0.05000	25.1	-0.031		
1269	1.76200	2.990	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 11/30

Measured Values Scope 1/1, Sample 11/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.496	23.6	
2	4.946	23.8	0.450
3	5.215	24.0	0.269
4	5.411	24.2	0.196
5	5.569	24.4	0.158
6	5.700	24.6	0.131
7	5.814	24.7	0.114
8	5.916	24.9	0.102
9	5.998	24.9	0.082
10	6.083	25.0	0.085
11	6.162	25.1	0.079
12	6.236	25.1	0.074
13	6.304	25.1	0.068
14	6.371	25.1	0.067
15	6.427	25.2	0.056
16	6.485	25.2	0.058

Titration report

LabX


17	6.540	25.2	0.055
18	6.592	25.2	0.052
19	6.642	25.2	0.050
20	6.689	25.2	0.047
21	6.734	25.2	0.045
22	6.773	25.2	0.039
23	6.814	25.2	0.041
24	6.855	25.2	0.041
25	6.892	25.2	0.037
26	6.924	25.2	0.032
27	6.959	25.2	0.035
28	6.990	25.2	0.031
29	7.016	25.2	0.026
30	7.039	25.2	0.023
31	7.072	25.2	0.033
32	7.099	25.2	0.027
33	7.124	25.2	0.025
34	7.148	25.2	0.024
35	7.171	25.2	0.023
36	7.190	25.2	0.019
37	7.211	25.2	0.021
38	7.231	25.2	0.020
39	7.249	25.2	0.018
40	7.266	25.2	0.017
41	7.282	25.2	0.016
42	7.298	25.2	0.016
43	7.312	25.2	0.014
44	7.326	25.2	0.014
45	7.339	25.2	0.013
46	7.352	25.2	0.013
47	7.364	25.2	0.012
48	7.376	25.2	0.012
49	7.387	25.2	0.011
50	7.397	25.2	0.010
51	7.407	25.2	0.010
52	7.416	25.2	0.009
53	7.425	25.2	0.009
54	7.433	25.2	0.008
55	7.442	25.2	0.009
56	7.450	25.2	0.008
57	7.456	25.2	0.006
58	7.463	25.2	0.007
59	7.470	25.2	0.007
60	7.477	25.2	0.007

Titration (EQP) [1] Scope 1/1, Sample 11/30


E - V Curve Scope 1/1, Sample 11/30

Titration report

LabX

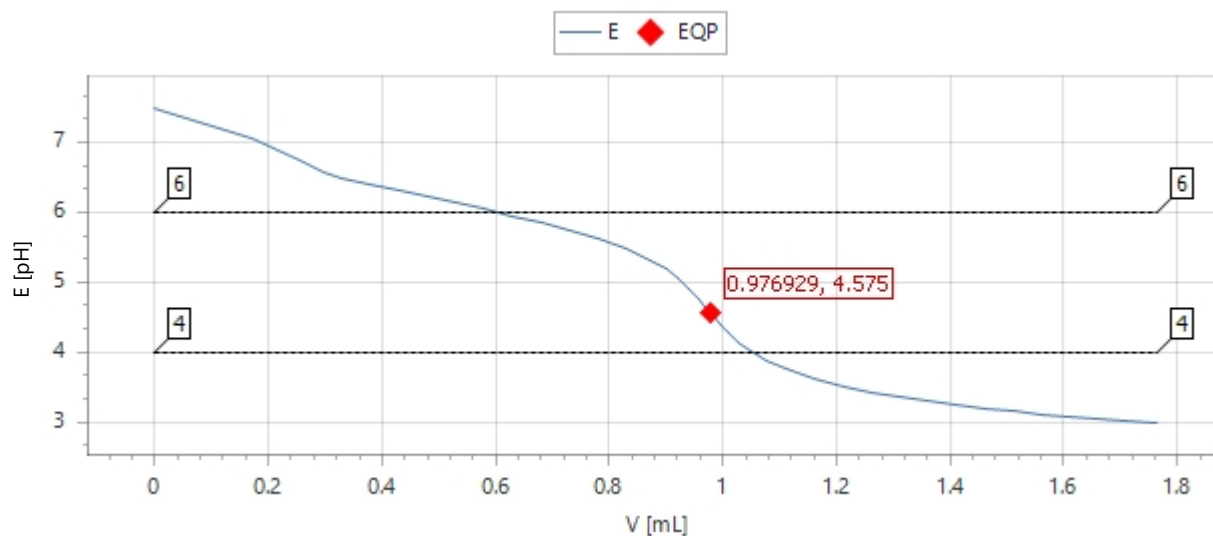

dE/dV - V Curve Scope 1/1, Sample 11/30

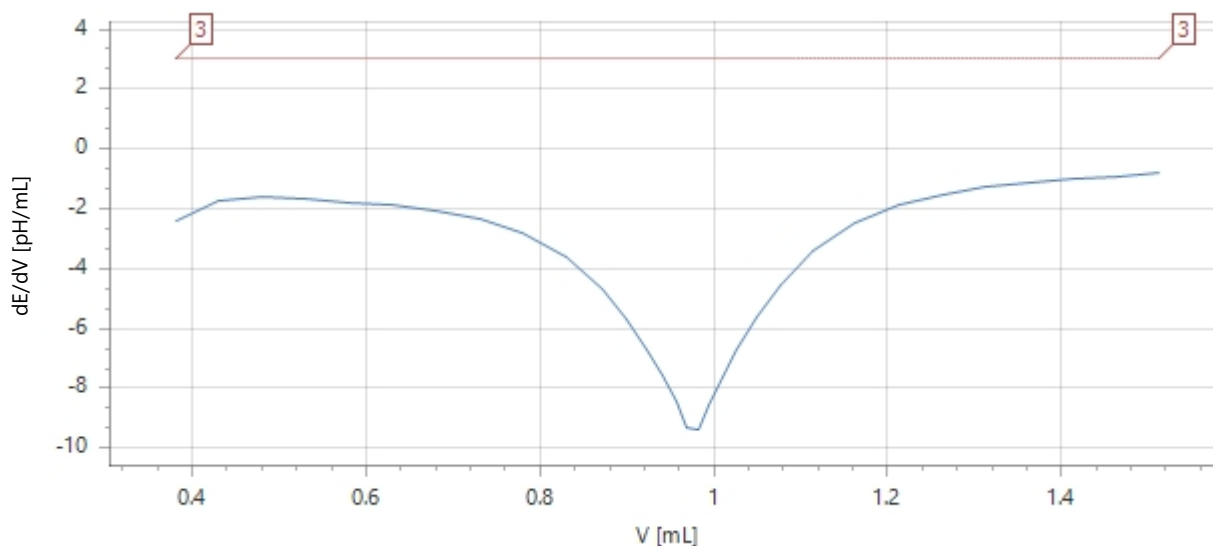

Measured Values Scope 1/1, Sample 11/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.490		25.2			
28	0.17150	7.054	0.17150	25.2	-0.436		
58	0.25725	6.742	0.08575	25.2	-0.312		
88	0.30000	6.571	0.04275	25.2	-0.171		
179	0.33075	6.472	0.03075	25.2	-0.099		
209	0.38075	6.389	0.05000	25.2	-0.083	-2.42	
239	0.43075	6.304	0.05000	25.2	-0.085	-1.75	
269	0.48075	6.219	0.05000	25.2	-0.085	-1.58	
300	0.53075	6.131	0.05000	25.2	-0.088	-1.68	
330	0.58075	6.043	0.05000	25.2	-0.088	-1.79	
360	0.63075	5.951	0.05000	25.1	-0.092	-1.87	
390	0.68075	5.854	0.05000	25.1	-0.097	-2.04	
421	0.73075	5.743	0.05000	25.1	-0.111	-2.32	
451	0.78075	5.620	0.05000	25.1	-0.123	-2.79	

Titration report

LabX

481	0.83075	5.471	0.05000	25.1	-0.149	-3.61	EQP
511	0.87200	5.311	0.04125	25.1	-0.160	-4.69	
542	0.89900	5.183	0.02700	25.1	-0.128	-5.66	
572	0.92225	5.045	0.02325	25.1	-0.138	-6.71	
602	0.94050	4.913	0.01825	25.1	-0.132	-7.59	
632	0.95600	4.780	0.01550	25.1	-0.133	-8.46	
662	0.96925	4.652	0.01325	25.1	-0.128	-9.30	
	0.976929 ⁽⁹⁾	4.575 ⁽⁹⁾					
693	0.98175	4.526	0.01250	25.1	-0.126	-9.35	
723	0.99475	4.401	0.01300	25.1	-0.125	-8.55	
753	1.00925	4.275	0.01450	25.1	-0.126	-7.62	
784	1.02650	4.149	0.01725	25.1	-0.126	-6.71	
814	1.04900	4.013	0.02250	25.1	-0.136	-5.61	
845	1.07575	3.888	0.02675	25.1	-0.125	-4.53	
875	1.11375	3.752	0.03800	25.0	-0.136	-3.40	
905	1.16275	3.620	0.04900	25.1	-0.132	-2.46	
935	1.21275	3.515	0.05000	25.0	-0.105	-1.89	
966	1.26275	3.431	0.05000	25.0	-0.084	-1.53	
996	1.31275	3.360	0.05000	25.0	-0.071	-1.29	
1026	1.36275	3.298	0.05000	25.0	-0.062	-1.13	
1056	1.41275	3.244	0.05000	25.0	-0.054	-1.00	
1087	1.46275	3.197	0.05000	25.0	-0.047	-0.90	
1117	1.51275	3.154	0.05000	25.0	-0.043	-0.82	
1147	1.56275	3.115	0.05000	25.0	-0.039		
1177	1.61275	3.079	0.05000	25.0	-0.036		
1208	1.66275	3.045	0.05000	25.0	-0.034		
1238	1.71275	3.014	0.05000	25.0	-0.031		
1268	1.76275	2.986	0.05000	25.0	-0.028		

Measure (normal) [1] Scope 1/1, Sample 12/30

Measured Values Scope 1/1, Sample 12/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.615	23.7	
2	4.961	23.9	0.346
3	5.225	24.1	0.264
4	5.429	24.3	0.204
5	5.597	24.5	0.168
6	5.740	24.7	0.143
7	5.866	24.8	0.126
8	5.979	24.9	0.113
9	6.073	25.0	0.094
10	6.171	25.0	0.098
11	6.261	25.1	0.090
12	6.345	25.1	0.084
13	6.424	25.1	0.079
14	6.497	25.2	0.073
15	6.568	25.2	0.071
16	6.629	25.2	0.061

Titration report

LabX


17	6.692	25.2	0.063
18	6.753	25.2	0.061
19	6.812	25.2	0.059
20	6.867	25.2	0.055
21	6.914	25.2	0.047
22	6.955	25.2	0.041
23	6.999	25.2	0.044
24	7.033	25.2	0.034
25	7.079	25.2	0.046
26	7.118	25.2	0.039
27	7.154	25.2	0.036
28	7.185	25.2	0.031
29	7.213	25.2	0.028
30	7.241	25.2	0.028
31	7.269	25.2	0.028
32	7.295	25.2	0.026
33	7.320	25.2	0.025
34	7.344	25.2	0.024
35	7.366	25.3	0.022
36	7.384	25.2	0.018
37	7.404	25.2	0.020
38	7.422	25.2	0.018
39	7.440	25.2	0.018
40	7.456	25.3	0.016
41	7.471	25.2	0.015
42	7.484	25.3	0.013
43	7.497	25.3	0.013
44	7.509	25.2	0.012
45	7.520	25.2	0.011
46	7.531	25.2	0.011
47	7.541	25.2	0.010
48	7.551	25.2	0.010
49	7.561	25.2	0.010
50	7.569	25.2	0.008
51	7.577	25.2	0.008
52	7.585	25.2	0.008
53	7.592	25.2	0.007
54	7.599	25.2	0.007
55	7.605	25.2	0.006
56	7.611	25.2	0.006
57	7.616	25.2	0.005
58	7.622	25.2	0.006
59	7.627	25.2	0.005
60	7.632	25.2	0.005

Titration (EQP) [1] Scope 1/1, Sample 12/30


E - V Curve Scope 1/1, Sample 12/30

Titration report

LabX

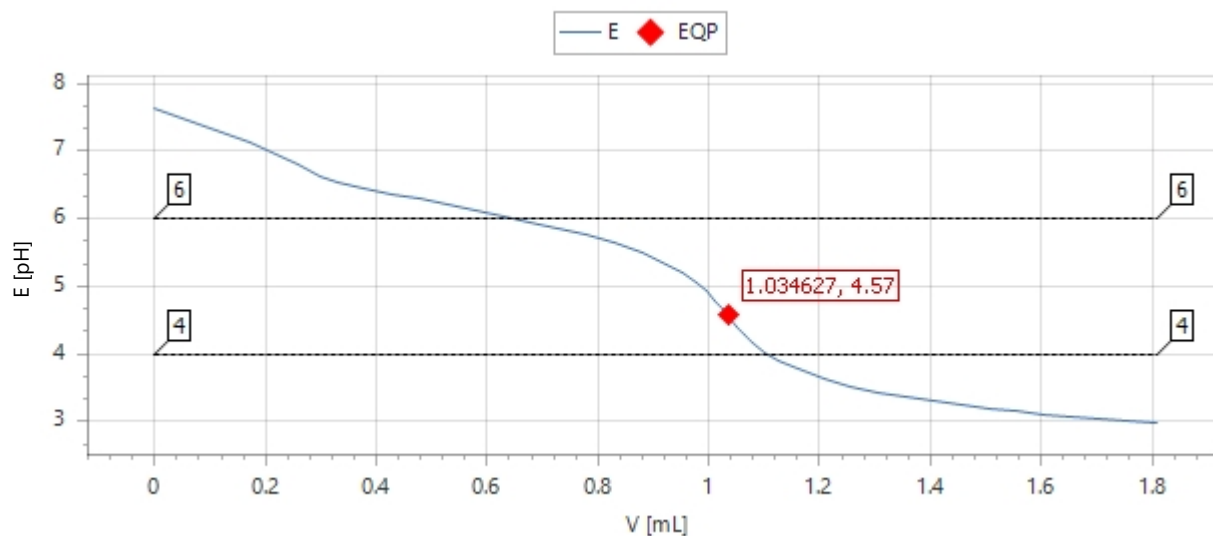

dE/dV - V Curve Scope 1/1, Sample 12/30

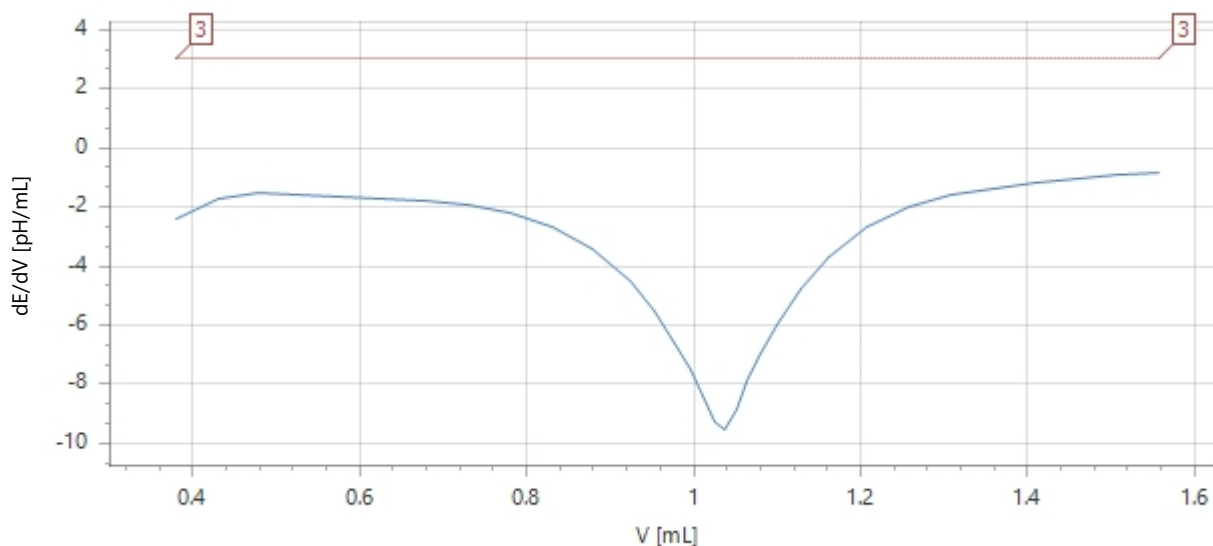

Measured Values Scope 1/1, Sample 12/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.643		25.2			
28	0.17150	7.131	0.17150	25.2	-0.512		
58	0.25725	6.806	0.08575	25.2	-0.325		
88	0.30000	6.628	0.04275	25.2	-0.178		
179	0.32975	6.534	0.02975	25.2	-0.094		
209	0.37975	6.456	0.05000	25.2	-0.078	-2.38	
239	0.42975	6.370	0.05000	25.2	-0.086	-1.69	
269	0.47975	6.288	0.05000	25.2	-0.082	-1.50	
300	0.52975	6.207	0.05000	25.2	-0.081	-1.59	
330	0.57975	6.122	0.05000	25.2	-0.085	-1.65	
360	0.62975	6.041	0.05000	25.1	-0.081	-1.69	
390	0.67975	5.953	0.05000	25.1	-0.088	-1.78	
421	0.72975	5.860	0.05000	25.1	-0.093	-1.94	
451	0.77975	5.756	0.05000	25.1	-0.104	-2.21	

Titration report

LabX

481	0.82975	5.637	0.05000	25.1	-0.119	-2.66	EQP
511	0.87975	5.493	0.05000	25.1	-0.144	-3.41	
542	0.92425	5.331	0.04450	25.1	-0.162	-4.51	
572	0.95350	5.197	0.02925	25.1	-0.134	-5.54	
602	0.97700	5.060	0.02350	25.1	-0.137	-6.58	
632	0.99525	4.931	0.01825	25.1	-0.129	-7.47	
662	1.01100	4.799	0.01575	25.1	-0.132	-8.35	
693	1.02450	4.671	0.01350	25.1	-0.128	-9.25	
	1.034627 ⁽⁹⁾	4.570 ⁽⁹⁾					
723	1.03725	4.544	0.01275	25.1	-0.127	-9.50	
753	1.05025	4.414	0.01300	25.1	-0.130	-8.85	
784	1.06375	4.293	0.01350	25.1	-0.121	-7.89	
814	1.08075	4.159	0.01700	25.1	-0.134	-6.96	
844	1.10050	4.035	0.01975	25.1	-0.124	-5.97	
874	1.12725	3.901	0.02675	25.1	-0.134	-4.80	
904	1.16150	3.770	0.03425	25.1	-0.131	-3.68	
935	1.20775	3.635	0.04625	25.1	-0.135	-2.68	
965	1.25775	3.524	0.05000	25.1	-0.111	-2.02	
995	1.30775	3.435	0.05000	25.1	-0.089	-1.61	
1026	1.35775	3.361	0.05000	25.1	-0.074	-1.35	
1056	1.40775	3.297	0.05000	25.1	-0.064	-1.17	
1086	1.45775	3.241	0.05000	25.1	-0.056	-1.04	
1116	1.50775	3.192	0.05000	25.0	-0.049	-0.93	
1146	1.55775	3.147	0.05000	25.0	-0.045	-0.84	
1177	1.60775	3.107	0.05000	25.1	-0.040		
1207	1.65775	3.071	0.05000	25.1	-0.036		
1237	1.70775	3.037	0.05000	25.1	-0.034		
1267	1.75775	3.005	0.05000	25.1	-0.032		
1298	1.80775	2.976	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 13/30

Measured Values Scope 1/1, Sample 13/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.320	23.5	
2	4.905	23.7	0.585
3	5.198	23.9	0.293
4	5.405	24.2	0.207
5	5.574	24.4	0.169
6	5.719	24.6	0.145
7	5.849	24.7	0.130
8	5.954	24.8	0.105
9	6.062	24.9	0.108
10	6.162	25.0	0.100
11	6.255	25.0	0.093
12	6.343	25.1	0.088
13	6.425	25.1	0.082
14	6.501	25.1	0.076
15	6.566	25.1	0.065

Titration report

LabX


16	6.634	25.2	0.068
17	6.698	25.2	0.064
18	6.759	25.2	0.061
19	6.816	25.2	0.057
20	6.873	25.2	0.057
21	6.920	25.2	0.047
22	6.962	25.2	0.042
23	7.007	25.2	0.045
24	7.044	25.2	0.037
25	7.088	25.2	0.044
26	7.127	25.2	0.039
27	7.162	25.2	0.035
28	7.195	25.2	0.033
29	7.223	25.2	0.028
30	7.253	25.2	0.030
31	7.281	25.2	0.028
32	7.308	25.2	0.027
33	7.333	25.2	0.025
34	7.358	25.2	0.025
35	7.380	25.2	0.022
36	7.398	25.2	0.018
37	7.420	25.2	0.022
38	7.431	25.2	0.011
39	7.450	25.2	0.019
40	7.470	25.2	0.020
41	7.494	25.2	0.024
42	7.500	25.2	0.006
43	7.511	25.2	0.011
44	7.524	25.2	0.013
45	7.536	25.2	0.012
46	7.547	25.2	0.011
47	7.558	25.2	0.011
48	7.568	25.2	0.010
49	7.577	25.2	0.009
50	7.585	25.2	0.008
51	7.593	25.2	0.008
52	7.601	25.2	0.008
53	7.609	25.2	0.008
54	7.616	25.2	0.007
55	7.622	25.2	0.006
56	7.628	25.2	0.006
57	7.634	25.2	0.006
58	7.639	25.2	0.005
59	7.644	25.2	0.005
60	7.649	25.2	0.005

Titration (EQP) [1] Scope 1/1, Sample 13/30


E - V Curve Scope 1/1, Sample 13/30

Titration report

LabX

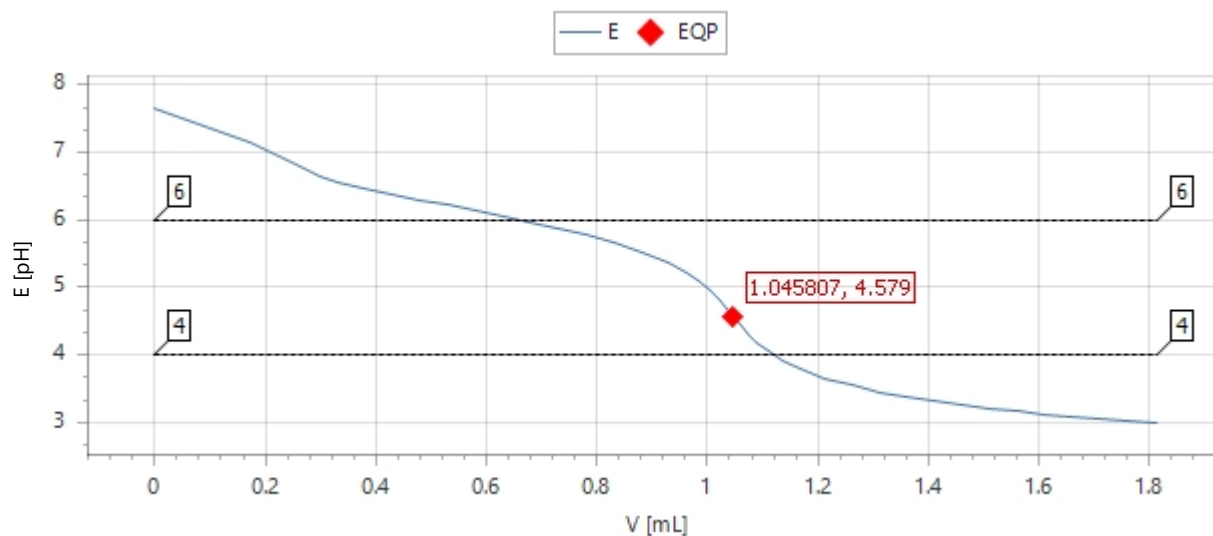

dE/dV - V Curve Scope 1/1, Sample 13/30

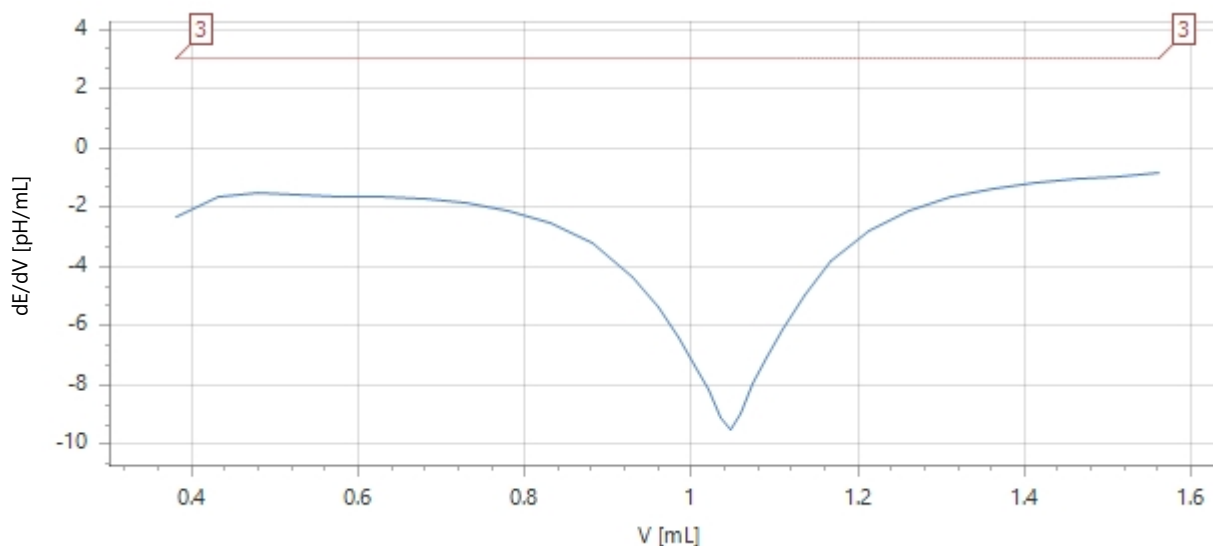

Measured Values Scope 1/1, Sample 13/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.660		25.2			
29	0.17150	7.140	0.17150	25.2	-0.520		
59	0.25725	6.814	0.08575	25.2	-0.326		
89	0.30000	6.640	0.04275	25.2	-0.174		
180	0.33075	6.543	0.03075	25.2	-0.097		
210	0.38075	6.467	0.05000	25.2	-0.076	-2.35	
240	0.43075	6.383	0.05000	25.2	-0.084	-1.67	
270	0.48075	6.301	0.05000	25.2	-0.082	-1.50	
301	0.53075	6.217	0.05000	25.1	-0.084	-1.58	
331	0.58075	6.137	0.05000	25.1	-0.080	-1.64	
361	0.63075	6.054	0.05000	25.1	-0.083	-1.66	
391	0.68075	5.968	0.05000	25.1	-0.086	-1.73	
422	0.73075	5.878	0.05000	25.1	-0.090	-1.88	
452	0.78075	5.778	0.05000	25.1	-0.100	-2.13	

Titration report

LabX

482	0.83075	5.663	0.05000	25.1	-0.115	-2.52	EQP
512	0.88075	5.526	0.05000	25.1	-0.137	-3.20	
543	0.93025	5.354	0.04950	25.1	-0.172	-4.33	
573	0.96075	5.219	0.03050	25.1	-0.135	-5.35	
603	0.98475	5.087	0.02400	25.0	-0.132	-6.36	
633	1.00475	4.949	0.02000	25.0	-0.138	-7.30	
664	1.02050	4.823	0.01575	25.0	-0.126	-8.13	
694	1.03500	4.687	0.01450	25.0	-0.136	-9.08	
	1.045807 ⁽⁹⁾	4.579 ⁽⁹⁾					
724	1.04725	4.565	0.01225	25.0	-0.122	-9.47	
754	1.06000	4.440	0.01275	25.0	-0.125	-8.94	
784	1.07400	4.311	0.01400	25.0	-0.129	-7.97	
815	1.08975	4.185	0.01575	25.0	-0.126	-7.08	
845	1.10900	4.060	0.01925	25.0	-0.125	-6.11	
875	1.13500	3.926	0.02600	25.0	-0.134	-4.95	
906	1.16800	3.794	0.03300	25.0	-0.132	-3.83	
936	1.21200	3.661	0.04400	25.0	-0.133	-2.81	
966	1.26200	3.546	0.05000	25.0	-0.115	-2.10	
996	1.31200	3.453	0.05000	25.0	-0.093	-1.66	
1026	1.36200	3.377	0.05000	25.0	-0.076	-1.38	
1057	1.41200	3.312	0.05000	25.0	-0.065	-1.19	
1087	1.46200	3.255	0.05000	25.0	-0.057	-1.05	
1117	1.51200	3.205	0.05000	25.0	-0.050	-0.94	
1147	1.56200	3.160	0.05000	25.0	-0.045	-0.86	
1178	1.61200	3.119	0.05000	25.0	-0.041		
1208	1.66200	3.081	0.05000	25.0	-0.038		
1238	1.71200	3.047	0.05000	25.0	-0.034		
1268	1.76200	3.015	0.05000	25.0	-0.032		
1298	1.81200	2.986	0.05000	25.0	-0.029		

Measure (normal) [1] Scope 1/1, Sample 14/30

Measured Values Scope 1/1, Sample 14/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.710	23.3	
2	5.015	23.5	0.305
3	5.270	23.7	0.255
4	5.470	24.0	0.200
5	5.632	24.2	0.162
6	5.767	24.4	0.135
7	5.885	24.6	0.118
8	5.982	24.7	0.097
9	6.080	24.8	0.098
10	6.170	24.9	0.090
11	6.253	24.9	0.083
12	6.331	25.0	0.078
13	6.402	25.0	0.071
14	6.468	25.1	0.066
15	6.526	25.1	0.058

Titration report

LabX


16	6.588	25.1	0.062
17	6.646	25.1	0.058
18	6.700	25.1	0.054
19	6.752	25.1	0.052
20	6.802	25.1	0.050
21	6.848	25.1	0.046
22	6.887	25.1	0.039
23	6.924	25.2	0.037
24	6.962	25.2	0.038
25	6.998	25.2	0.036
26	7.025	25.2	0.027
27	7.059	25.2	0.034
28	7.091	25.2	0.032
29	7.118	25.2	0.027
30	7.146	25.2	0.028
31	7.171	25.2	0.025
32	7.194	25.2	0.023
33	7.216	25.2	0.022
34	7.238	25.2	0.022
35	7.259	25.2	0.021
36	7.276	25.2	0.017
37	7.294	25.2	0.018
38	7.310	25.2	0.016
39	7.327	25.2	0.017
40	7.343	25.2	0.016
41	7.357	25.2	0.014
42	7.371	25.2	0.014
43	7.383	25.2	0.012
44	7.395	25.2	0.012
45	7.406	25.2	0.011
46	7.417	25.2	0.011
47	7.427	25.2	0.010
48	7.437	25.2	0.010
49	7.446	25.2	0.009
50	7.454	25.2	0.008
51	7.463	25.2	0.009
52	7.470	25.2	0.007
53	7.477	25.2	0.007
54	7.484	25.2	0.007
55	7.491	25.2	0.007
56	7.496	25.2	0.005
57	7.502	25.2	0.006
58	7.508	25.2	0.006
59	7.513	25.2	0.005
60	7.518	25.2	0.005

Titration (EQP) [1] Scope 1/1, Sample 14/30


E - V Curve Scope 1/1, Sample 14/30

Titration report

LabX

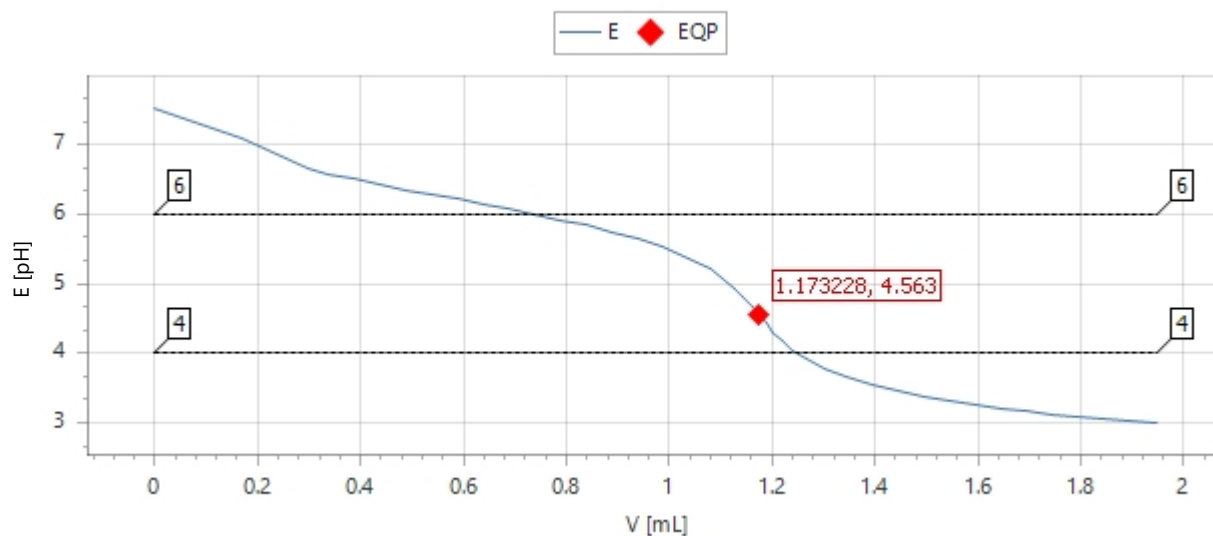

dE/dV - V Curve Scope 1/1, Sample 14/30

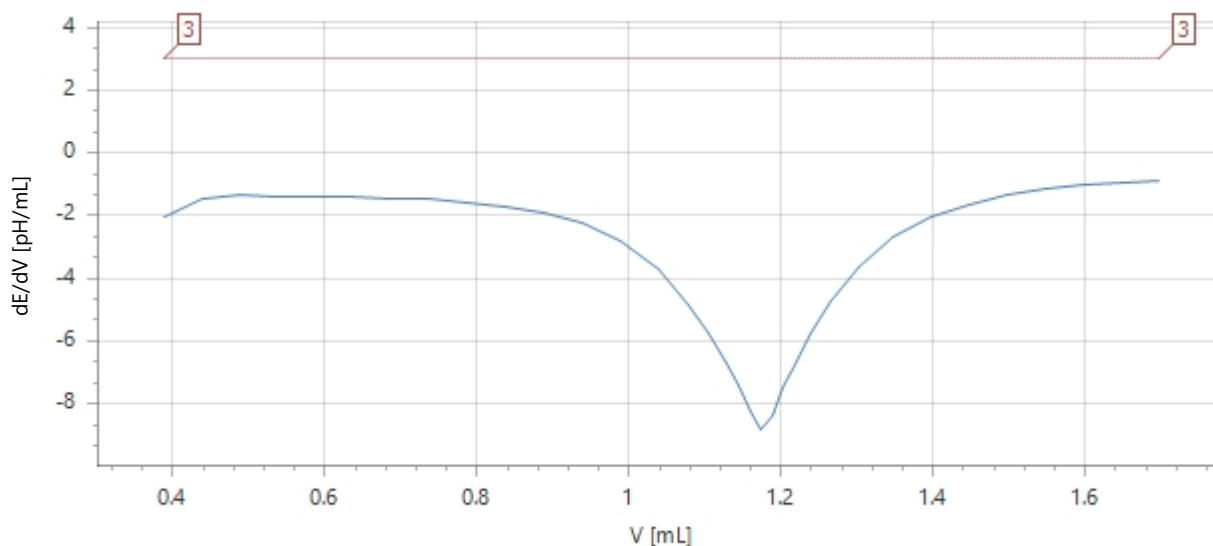

Measured Values Scope 1/1, Sample 14/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.529		25.2			
28	0.17150	7.091	0.17150	25.2	-0.438		
58	0.25725	6.806	0.08575	25.2	-0.285		
88	0.30000	6.661	0.04275	25.1	-0.145		
179	0.33875	6.574	0.03875	25.1	-0.087		
209	0.38875	6.499	0.05000	25.1	-0.075	-2.01	
239	0.43875	6.424	0.05000	25.1	-0.075	-1.45	
269	0.48875	6.353	0.05000	25.1	-0.071	-1.34	
300	0.53875	6.282	0.05000	25.1	-0.071	-1.38	
330	0.58875	6.212	0.05000	25.1	-0.070	-1.40	
360	0.63875	6.137	0.05000	25.1	-0.075	-1.41	
390	0.68875	6.070	0.05000	25.1	-0.067	-1.43	
421	0.73875	5.998	0.05000	25.1	-0.072	-1.48	
451	0.78875	5.920	0.05000	25.1	-0.078	-1.56	

Titration report

LabX

481	0.83875	5.838	0.05000	25.1	-0.082	-1.69	EQP
511	0.88875	5.747	0.05000	25.1	-0.091	-1.92	
542	0.93875	5.644	0.05000	25.1	-0.103	-2.25	
572	0.98875	5.521	0.05000	25.1	-0.123	-2.79	
602	1.03875	5.371	0.05000	25.1	-0.150	-3.70	
632	1.07925	5.209	0.04050	25.1	-0.162	-4.83	
662	1.10500	5.083	0.02575	25.1	-0.126	-5.77	
693	1.12775	4.940	0.02275	25.1	-0.143	-6.71	
723	1.14425	4.820	0.01650	25.1	-0.120	-7.44	
753	1.16050	4.680	0.01625	25.1	-0.140	-8.28	
	1.173228 ⁽⁹⁾	4.563 ⁽⁹⁾					
784	1.17375	4.558	0.01325	25.1	-0.122	-8.80	
814	1.18750	4.431	0.01375	25.1	-0.127	-8.37	
844	1.20200	4.307	0.01450	25.1	-0.124	-7.50	
874	1.21925	4.178	0.01725	25.1	-0.129	-6.66	
904	1.24000	4.049	0.02075	25.1	-0.129	-5.74	
935	1.26625	3.920	0.02625	25.1	-0.129	-4.71	
965	1.30100	3.787	0.03475	25.1	-0.133	-3.64	
995	1.34650	3.654	0.04550	25.1	-0.133	-2.69	
1025	1.39650	3.542	0.05000	25.1	-0.112	-2.03	
1056	1.44650	3.452	0.05000	25.1	-0.090	-1.62	
1086	1.49650	3.378	0.05000	25.1	-0.074	-1.35	
1116	1.54650	3.314	0.05000	25.1	-0.064	-1.17	
1146	1.59650	3.258	0.05000	25.1	-0.056	-1.04	
1177	1.64650	3.208	0.05000	25.1	-0.050	-0.93	
1207	1.69650	3.164	0.05000	25.1	-0.044	-0.85	
1237	1.74650	3.123	0.05000	25.1	-0.041		
1267	1.79650	3.086	0.05000	25.1	-0.037		
1298	1.84650	3.052	0.05000	25.1	-0.034		
1328	1.89650	3.020	0.05000	25.1	-0.032		
1358	1.94650	2.991	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 15/30

Measured Values Scope 1/1, Sample 15/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.743	23.5	
2	5.044	23.7	0.301
3	5.295	23.9	0.251
4	5.490	24.1	0.195
5	5.651	24.4	0.161
6	5.787	24.6	0.136
7	5.906	24.7	0.119
8	6.013	24.8	0.107
9	6.102	24.9	0.089
10	6.194	25.0	0.092
11	6.278	25.0	0.084
12	6.355	25.1	0.077
13	6.429	25.1	0.074

Titration report

LabX


14	6.499	25.2	0.070
15	6.564	25.2	0.065
16	6.618	25.2	0.054
17	6.677	25.2	0.059
18	6.733	25.2	0.056
19	6.787	25.2	0.054
20	6.840	25.2	0.053
21	6.887	25.2	0.047
22	6.925	25.2	0.038
23	6.967	25.2	0.042
24	7.005	25.2	0.038
25	7.035	25.2	0.030
26	7.076	25.2	0.041
27	7.109	25.2	0.033
28	7.141	25.2	0.032
29	7.166	25.2	0.025
30	7.194	25.2	0.028
31	7.220	25.2	0.026
32	7.245	25.2	0.025
33	7.269	25.3	0.024
34	7.290	25.2	0.021
35	7.311	25.3	0.021
36	7.329	25.3	0.018
37	7.348	25.3	0.019
38	7.365	25.3	0.017
39	7.382	25.3	0.017
40	7.397	25.3	0.015
41	7.411	25.3	0.014
42	7.424	25.3	0.013
43	7.435	25.3	0.011
44	7.448	25.3	0.013
45	7.459	25.3	0.011
46	7.469	25.3	0.010
47	7.479	25.3	0.010
48	7.489	25.3	0.010
49	7.498	25.3	0.009
50	7.506	25.3	0.008
51	7.514	25.3	0.008
52	7.521	25.3	0.007
53	7.528	25.3	0.007
54	7.535	25.2	0.007
55	7.541	25.2	0.006
56	7.548	25.3	0.007
57	7.553	25.3	0.005
58	7.558	25.3	0.005
59	7.564	25.3	0.006
60	7.568	25.3	0.004

Titration (EQP) [1] Scope 1/1, Sample 15/30


E - V Curve Scope 1/1, Sample 15/30

Titration report

LabX

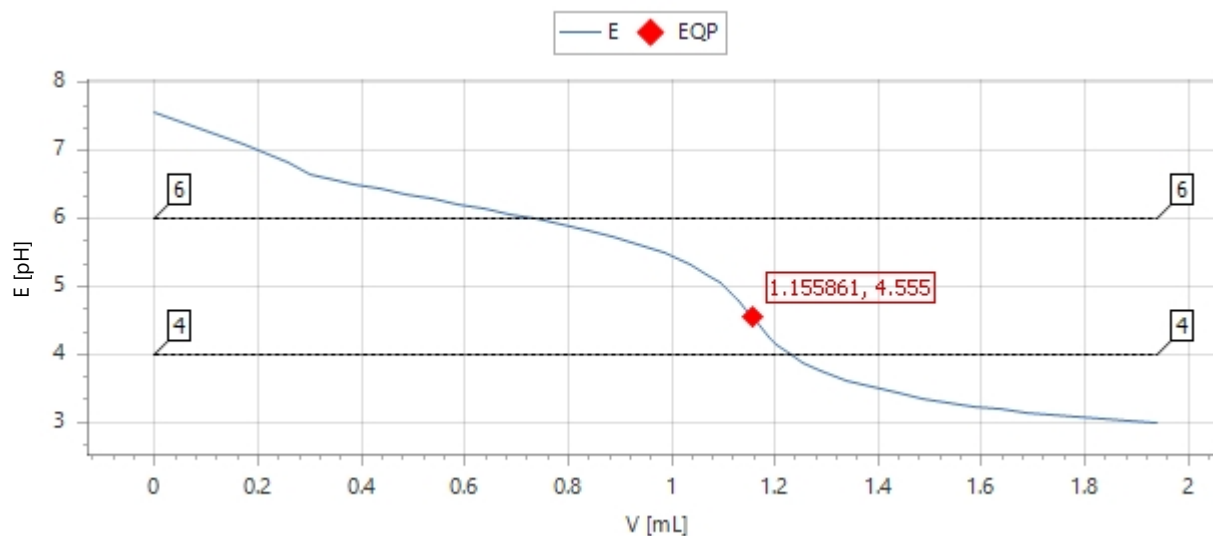

dE/dV - V Curve Scope 1/1, Sample 15/30

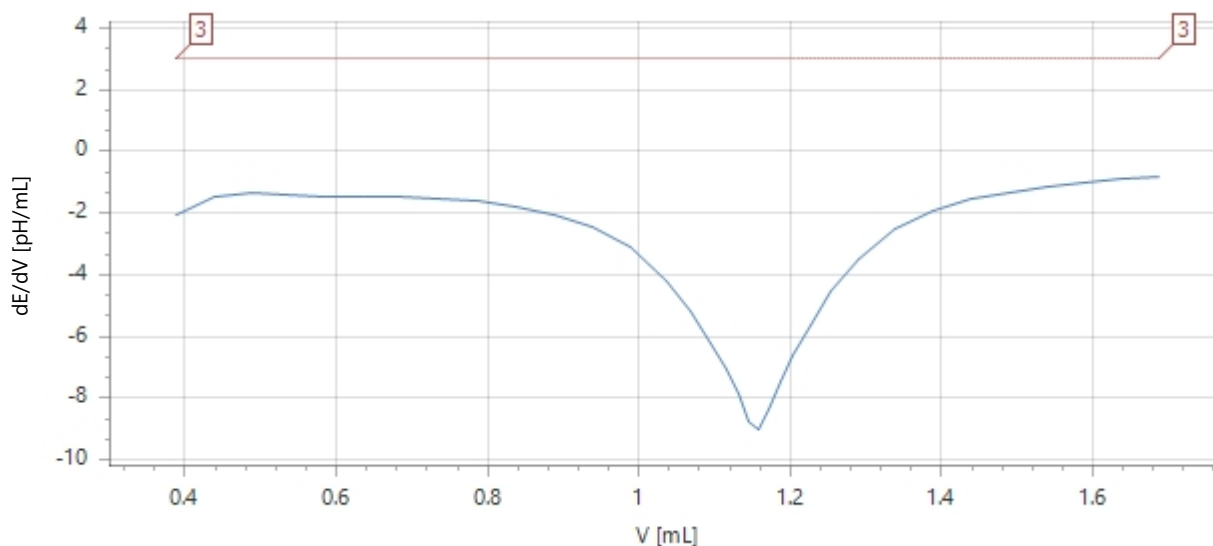

Measured Values Scope 1/1, Sample 15/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.579		25.3			
28	0.17150	7.109	0.17150	25.2	-0.470		
58	0.25725	6.814	0.08575	25.2	-0.295		
88	0.30000	6.663	0.04275	25.2	-0.151		
180	0.33750	6.579	0.03750	25.2	-0.084		
210	0.38750	6.506	0.05000	25.2	-0.073	-2.04	
240	0.43750	6.430	0.05000	25.2	-0.076	-1.47	
270	0.48750	6.354	0.05000	25.2	-0.076	-1.35	
301	0.53750	6.283	0.05000	25.2	-0.071	-1.41	
331	0.58750	6.211	0.05000	25.2	-0.072	-1.44	
361	0.63750	6.139	0.05000	25.2	-0.072	-1.44	
391	0.68750	6.068	0.05000	25.2	-0.071	-1.47	
422	0.73750	5.989	0.05000	25.2	-0.079	-1.53	
452	0.78750	5.910	0.05000	25.2	-0.079	-1.61	

Titration report

LabX

482	0.83750	5.829	0.05000	25.2	-0.081	-1.77	EQP
512	0.88750	5.731	0.05000	25.2	-0.098	-2.02	
542	0.93750	5.624	0.05000	25.2	-0.107	-2.41	
573	0.98750	5.491	0.05000	25.2	-0.133	-3.09	
603	1.03750	5.321	0.05000	25.2	-0.170	-4.19	
633	1.06925	5.185	0.03175	25.2	-0.136	-5.19	
664	1.09450	5.049	0.02525	25.2	-0.136	-6.20	
694	1.11450	4.914	0.02000	25.2	-0.135	-7.06	
724	1.13075	4.786	0.01625	25.2	-0.128	-7.84	
754	1.14550	4.655	0.01475	25.2	-0.131	-8.73	
	1.155861 ⁽⁹⁾	4.555 ⁽⁹⁾					
784	1.15900	4.525	0.01350	25.2	-0.130	-8.99	
815	1.17225	4.402	0.01325	25.2	-0.123	-8.38	
845	1.18725	4.273	0.01500	25.2	-0.129	-7.45	
875	1.20425	4.148	0.01700	25.2	-0.125	-6.61	
906	1.22575	4.019	0.02150	25.2	-0.129	-5.63	
936	1.25350	3.888	0.02775	25.2	-0.131	-4.54	
966	1.28975	3.755	0.03625	25.2	-0.133	-3.47	
996	1.33800	3.622	0.04825	25.2	-0.133	-2.53	
1026	1.38800	3.515	0.05000	25.2	-0.107	-1.93	
1057	1.43800	3.429	0.05000	25.2	-0.086	-1.56	
1087	1.48800	3.357	0.05000	25.2	-0.072	-1.31	
1117	1.53800	3.295	0.05000	25.2	-0.062	-1.14	
1147	1.58800	3.240	0.05000	25.2	-0.055	-1.02	
1178	1.63800	3.192	0.05000	25.2	-0.048	-0.91	
1208	1.68800	3.148	0.05000	25.2	-0.044	-0.83	
1238	1.73800	3.108	0.05000	25.1	-0.040		
1268	1.78800	3.072	0.05000	25.1	-0.036		
1299	1.83800	3.038	0.05000	25.1	-0.034		
1329	1.88800	3.007	0.05000	25.1	-0.031		
1359	1.93800	2.978	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 16/30

Measured Values Scope 1/1, Sample 16/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.695	23.3	
2	5.004	23.5	0.309
3	5.277	23.7	0.273
4	5.482	24.0	0.205
5	5.650	24.3	0.168
6	5.794	24.5	0.144
7	5.923	24.7	0.129
8	6.029	24.8	0.106
9	6.135	24.9	0.106
10	6.232	25.0	0.097
11	6.324	25.0	0.092
12	6.410	25.1	0.086
13	6.489	25.1	0.079

Titration report

LabX


14	6.563	25.2	0.074
15	6.625	25.2	0.062
16	6.692	25.2	0.067
17	6.754	25.2	0.062
18	6.813	25.2	0.059
19	6.870	25.2	0.057
20	6.918	25.2	0.048
21	6.964	25.2	0.046
22	7.004	25.2	0.040
23	7.041	25.2	0.037
24	7.086	25.2	0.045
25	7.124	25.2	0.038
26	7.160	25.2	0.036
27	7.192	25.2	0.032
28	7.223	25.2	0.031
29	7.249	25.3	0.026
30	7.277	25.2	0.028
31	7.302	25.2	0.025
32	7.326	25.2	0.024
33	7.348	25.3	0.022
34	7.369	25.3	0.021
35	7.389	25.3	0.020
36	7.406	25.2	0.017
37	7.422	25.2	0.016
38	7.438	25.2	0.016
39	7.453	25.3	0.015
40	7.467	25.2	0.014
41	7.481	25.2	0.014
42	7.494	25.2	0.013
43	7.505	25.2	0.011
44	7.516	25.2	0.011
45	7.526	25.2	0.010
46	7.536	25.2	0.010
47	7.545	25.2	0.009
48	7.554	25.2	0.009
49	7.562	25.2	0.008
50	7.569	25.2	0.007
51	7.576	25.2	0.007
52	7.583	25.2	0.007
53	7.589	25.2	0.006
54	7.595	25.2	0.006
55	7.601	25.2	0.006
56	7.606	25.2	0.005
57	7.610	25.2	0.004
58	7.615	25.2	0.005
59	7.619	25.2	0.004
60	7.623	25.2	0.004

Titration (EQP) [1] Scope 1/1, Sample 16/30


E - V Curve Scope 1/1, Sample 16/30

Titration report

LabX

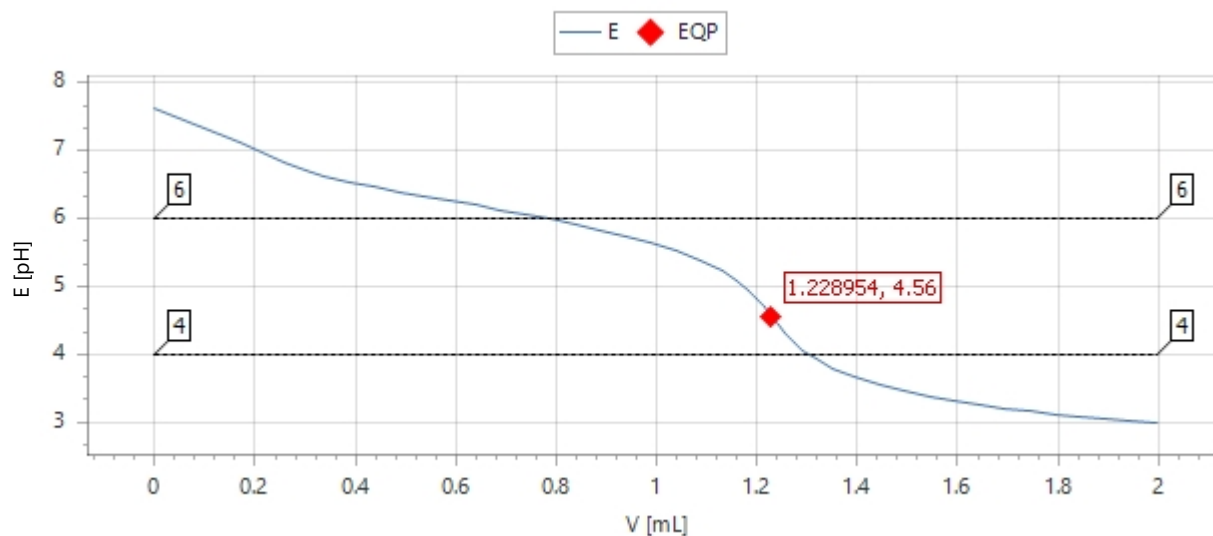

dE/dV - V Curve Scope 1/1, Sample 16/30

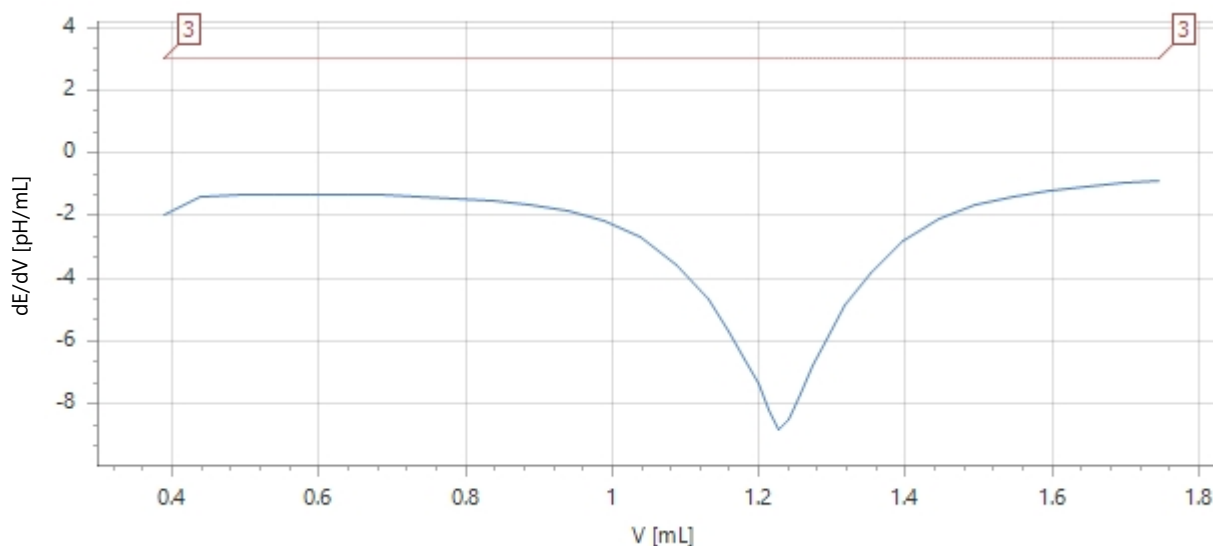

Measured Values Scope 1/1, Sample 16/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.633		25.2			
28	0.17150	7.137	0.17150	25.2	-0.496		
58	0.25725	6.843	0.08575	25.2	-0.294		
88	0.30000	6.696	0.04275	25.2	-0.147		
179	0.33900	6.616	0.03900	25.2	-0.080		
209	0.38900	6.544	0.05000	25.2	-0.072	-1.95	
239	0.43900	6.470	0.05000	25.2	-0.074	-1.40	
270	0.48900	6.398	0.05000	25.2	-0.072	-1.30	
300	0.53900	6.330	0.05000	25.2	-0.068	-1.34	
330	0.58900	6.262	0.05000	25.2	-0.068	-1.35	
360	0.63900	6.196	0.05000	25.2	-0.066	-1.33	
390	0.68900	6.128	0.05000	25.2	-0.068	-1.34	
421	0.73900	6.060	0.05000	25.2	-0.068	-1.38	
451	0.78900	5.990	0.05000	25.2	-0.070	-1.43	

Titration report

LabX

481	0.83900	5.916	0.05000	25.2	-0.074	-1.52	
511	0.88900	5.836	0.05000	25.2	-0.080	-1.64	
542	0.93900	5.747	0.05000	25.2	-0.089	-1.86	
572	0.98900	5.647	0.05000	25.2	-0.100	-2.18	
602	1.03900	5.530	0.05000	25.2	-0.117	-2.70	
632	1.08900	5.382	0.05000	25.2	-0.148	-3.56	
663	1.13100	5.222	0.04200	25.2	-0.160	-4.68	
693	1.15875	5.089	0.02775	25.2	-0.133	-5.68	
723	1.18125	4.951	0.02250	25.2	-0.138	-6.60	
753	1.19850	4.826	0.01725	25.2	-0.125	-7.35	
784	1.21425	4.695	0.01575	25.2	-0.131	-8.18	
814	1.22825	4.567	0.01400	25.2	-0.128	-8.80	
	1.228954 ⁽⁹⁾	4.560 ⁽⁹⁾					EQP
844	1.24200	4.438	0.01375	25.2	-0.129	-8.51	
874	1.25600	4.317	0.01400	25.2	-0.121	-7.66	
905	1.27300	4.185	0.01700	25.2	-0.132	-6.79	
935	1.29250	4.060	0.01950	25.2	-0.125	-5.92	
965	1.31800	3.929	0.02550	25.2	-0.131	-4.88	
995	1.35125	3.798	0.03325	25.2	-0.131	-3.80	
1026	1.39550	3.662	0.04425	25.2	-0.136	-2.81	
1056	1.44550	3.547	0.05000	25.1	-0.115	-2.11	
1086	1.49550	3.454	0.05000	25.1	-0.093	-1.67	
1116	1.54550	3.378	0.05000	25.1	-0.076	-1.39	
1146	1.59550	3.312	0.05000	25.1	-0.066	-1.20	
1177	1.64550	3.255	0.05000	25.1	-0.057	-1.06	
1207	1.69550	3.204	0.05000	25.1	-0.051	-0.95	
1237	1.74550	3.159	0.05000	25.1	-0.045	-0.86	
1268	1.79550	3.118	0.05000	25.1	-0.041		
1298	1.84550	3.081	0.05000	25.1	-0.037		
1328	1.89550	3.046	0.05000	25.1	-0.035		
1358	1.94550	3.014	0.05000	25.1	-0.032		
1388	1.99550	2.984	0.05000	25.1	-0.030		

Measure (normal) [1] Scope 1/1, Sample 17/30

Measured Values Scope 1/1, Sample 17/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.737	23.1	
2	5.037	23.3	0.300
3	5.299	23.6	0.262
4	5.503	23.9	0.204
5	5.672	24.2	0.169
6	5.817	24.4	0.145
7	5.945	24.6	0.128
8	6.052	24.7	0.107
9	6.164	24.9	0.112
10	6.266	24.9	0.102
11	6.359	25.0	0.093
12	6.448	25.1	0.089

Titration report

LabX

13	6.531	25.1	0.083
14	6.608	25.1	0.077
15	6.675	25.2	0.067
16	6.743	25.2	0.068
17	6.807	25.2	0.064
18	6.868	25.2	0.061
19	6.921	25.2	0.053
20	6.973	25.2	0.052
21	7.016	25.2	0.043
22	7.058	25.2	0.042
23	7.105	25.2	0.047
24	7.145	25.2	0.040
25	7.182	25.2	0.037
26	7.216	25.2	0.034
27	7.248	25.2	0.032
28	7.279	25.3	0.031
29	7.304	25.3	0.025
30	7.332	25.3	0.028
31	7.357	25.3	0.025
32	7.380	25.3	0.023
33	7.402	25.3	0.022
34	7.422	25.3	0.020
35	7.441	25.3	0.019
36	7.457	25.3	0.016
37	7.473	25.3	0.016
38	7.489	25.3	0.016
39	7.503	25.3	0.014
40	7.517	25.3	0.014
41	7.530	25.3	0.013
42	7.542	25.3	0.012
43	7.553	25.3	0.011
44	7.563	25.3	0.010
45	7.573	25.3	0.010
46	7.582	25.3	0.009
47	7.590	25.3	0.008
48	7.598	25.3	0.008
49	7.606	25.3	0.008
50	7.613	25.3	0.007
51	7.619	25.3	0.006
52	7.625	25.3	0.006
53	7.630	25.3	0.005
54	7.636	25.3	0.006
55	7.641	25.3	0.005
56	7.645	25.3	0.004
57	7.649	25.3	0.004
58	7.654	25.3	0.005
59	7.657	25.3	0.003
60	7.661	25.3	0.004


Titration (EQP) [1]

Scope 1/1, Sample 17/30


Titration report

LabX

E - V Curve Scope 1/1, Sample 17/30

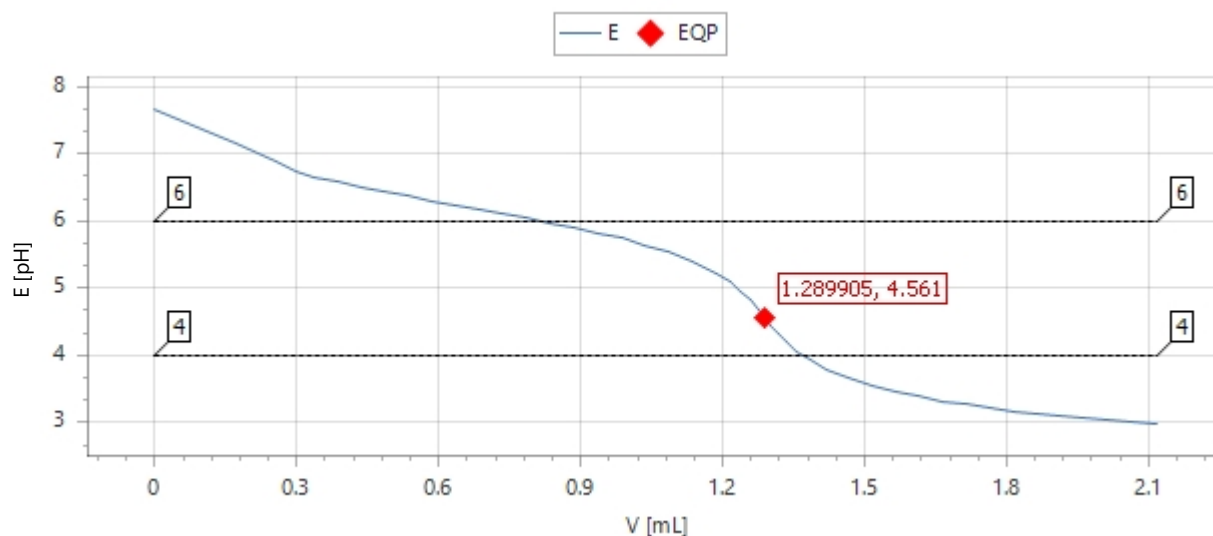

dE/dV - V Curve Scope 1/1, Sample 17/30

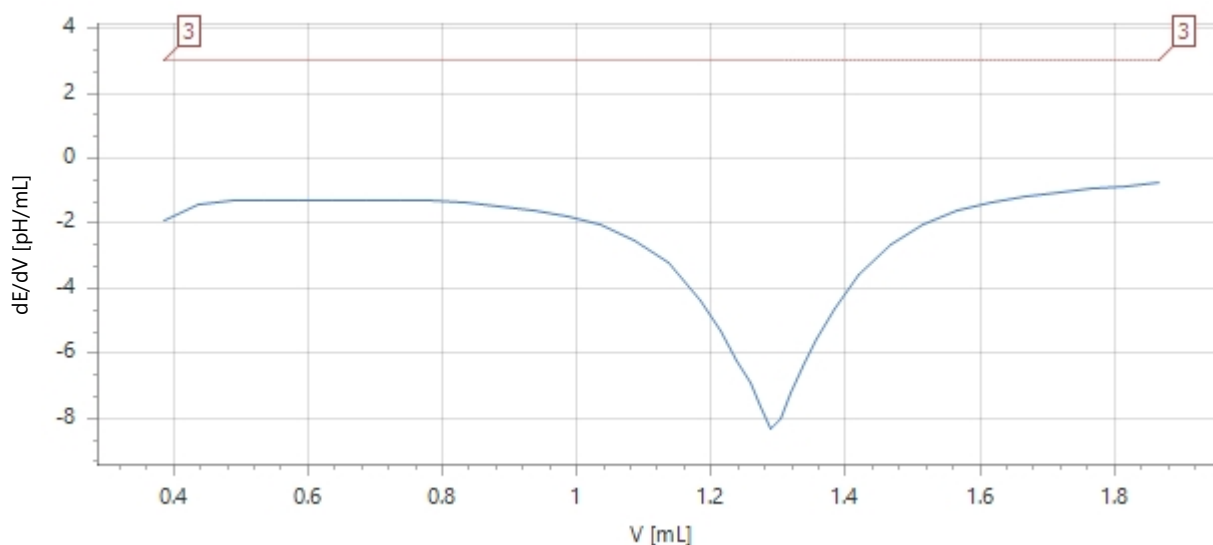

Measured Values Scope 1/1, Sample 17/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.670		25.3			
28	0.17150	7.164	0.17150	25.2	-0.506		
58	0.25725	6.877	0.08575	25.2	-0.287		
88	0.30000	6.726	0.04275	25.2	-0.151		
179	0.33650	6.649	0.03650	25.2	-0.077		
209	0.38650	6.578	0.05000	25.2	-0.071	-1.94	
239	0.43650	6.506	0.05000	25.2	-0.072	-1.39	
269	0.48650	6.435	0.05000	25.2	-0.071	-1.27	
300	0.53650	6.367	0.05000	25.2	-0.068	-1.32	
330	0.58650	6.301	0.05000	25.2	-0.066	-1.31	
360	0.63650	6.237	0.05000	25.2	-0.064	-1.28	
390	0.68650	6.171	0.05000	25.2	-0.066	-1.27	

Titration report

LabX

421	0.73650	6.109	0.05000	25.2	-0.062	-1.28	EQP
451	0.78650	6.045	0.05000	25.2	-0.064	-1.32	
481	0.83650	5.976	0.05000	25.2	-0.069	-1.38	
511	0.88650	5.905	0.05000	25.2	-0.071	-1.47	
542	0.93650	5.826	0.05000	25.2	-0.079	-1.59	
572	0.98650	5.742	0.05000	25.2	-0.084	-1.78	
602	1.03650	5.647	0.05000	25.2	-0.095	-2.06	
632	1.08650	5.535	0.05000	25.2	-0.112	-2.51	
662	1.13650	5.398	0.05000	25.1	-0.137	-3.23	
693	1.18550	5.224	0.04900	25.1	-0.174	-4.38	
723	1.21475	5.093	0.02925	25.1	-0.131	-5.31	
753	1.23900	4.956	0.02425	25.1	-0.137	-6.20	
784	1.25800	4.826	0.01900	25.1	-0.130	-6.93	
814	1.27450	4.695	0.01650	25.1	-0.131	-7.67	
844	1.28900	4.569	0.01450	25.1	-0.126	-8.30	
	1.289905	4.561					
874	1.30325	4.445	0.01425	25.1	-0.124	-8.03	
904	1.31875	4.317	0.01550	25.1	-0.128	-7.23	
935	1.33625	4.190	0.01750	25.1	-0.127	-6.44	
965	1.35775	4.059	0.02150	25.1	-0.131	-5.57	
995	1.38425	3.931	0.02650	25.1	-0.128	-4.60	
1026	1.41950	3.798	0.03525	25.1	-0.133	-3.58	
1056	1.46600	3.664	0.04650	25.1	-0.134	-2.65	
1086	1.51600	3.553	0.05000	25.1	-0.111	-2.02	
1116	1.56600	3.463	0.05000	25.1	-0.090	-1.62	
1146	1.61600	3.388	0.05000	25.1	-0.075	-1.36	
1177	1.66600	3.324	0.05000	25.1	-0.064	-1.17	
1207	1.71600	3.269	0.05000	25.1	-0.055	-1.04	
1237	1.76600	3.219	0.05000	25.1	-0.050	-0.93	
1267	1.81600	3.175	0.05000	25.1	-0.044	-0.84	
1298	1.86600	3.134	0.05000	25.1	-0.041	-0.77	
1328	1.91600	3.097	0.05000	25.1	-0.037		
1358	1.96600	3.063	0.05000	25.1	-0.034		
1388	2.01600	3.032	0.05000	25.1	-0.031		
1419	2.06600	3.002	0.05000	25.1	-0.030		
1449	2.11600	2.974	0.05000	25.1	-0.028		

Measure (normal) [1] Scope 1/1, Sample 18/30

Measured Values Scope 1/1, Sample 18/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.718	23.0	
2	5.008	23.2	0.290
3	5.285	23.5	0.277
4	5.502	23.9	0.217
5	5.682	24.2	0.180
6	5.840	24.4	0.158
7	5.983	24.6	0.143
8	6.100	24.7	0.117

Titration report

LabX

9	6.220	24.9	0.120
10	6.330	24.9	0.110
11	6.434	25.0	0.104
12	6.531	25.1	0.097
13	6.623	25.1	0.092
14	6.709	25.2	0.086
15	6.788	25.2	0.079
16	6.856	25.2	0.068
17	6.923	25.2	0.067
18	6.987	25.2	0.064
19	7.040	25.2	0.053
20	7.098	25.2	0.058
21	7.149	25.3	0.051
22	7.191	25.3	0.042
23	7.236	25.3	0.045
24	7.277	25.3	0.041
25	7.315	25.3	0.038
26	7.350	25.3	0.035
27	7.382	25.3	0.032
28	7.412	25.3	0.030
29	7.437	25.3	0.025
30	7.463	25.3	0.026
31	7.487	25.3	0.024
32	7.508	25.3	0.021
33	7.527	25.3	0.019
34	7.546	25.3	0.019
35	7.562	25.3	0.016
36	7.576	25.3	0.014
37	7.591	25.3	0.015
38	7.604	25.3	0.013
39	7.617	25.3	0.013
40	7.628	25.3	0.011
41	7.639	25.3	0.011
42	7.648	25.3	0.009
43	7.656	25.3	0.008
44	7.664	25.3	0.008
45	7.673	25.3	0.009
46	7.681	25.3	0.008
47	7.687	25.3	0.006
48	7.694	25.3	0.007
49	7.699	25.3	0.005
50	7.704	25.3	0.005
51	7.709	25.3	0.005
52	7.713	25.3	0.004
53	7.718	25.3	0.005
54	7.722	25.3	0.004
55	7.726	25.3	0.004
56	7.730	25.3	0.004
57	7.733	25.3	0.003
58	7.736	25.3	0.003
59	7.739	25.3	0.003

Titration report

LabX

60

7.742

25.3


0.003

Titration (EQP) [1]

Scope 1/1, Sample 18/30


E - V Curve

Scope 1/1, Sample 18/30

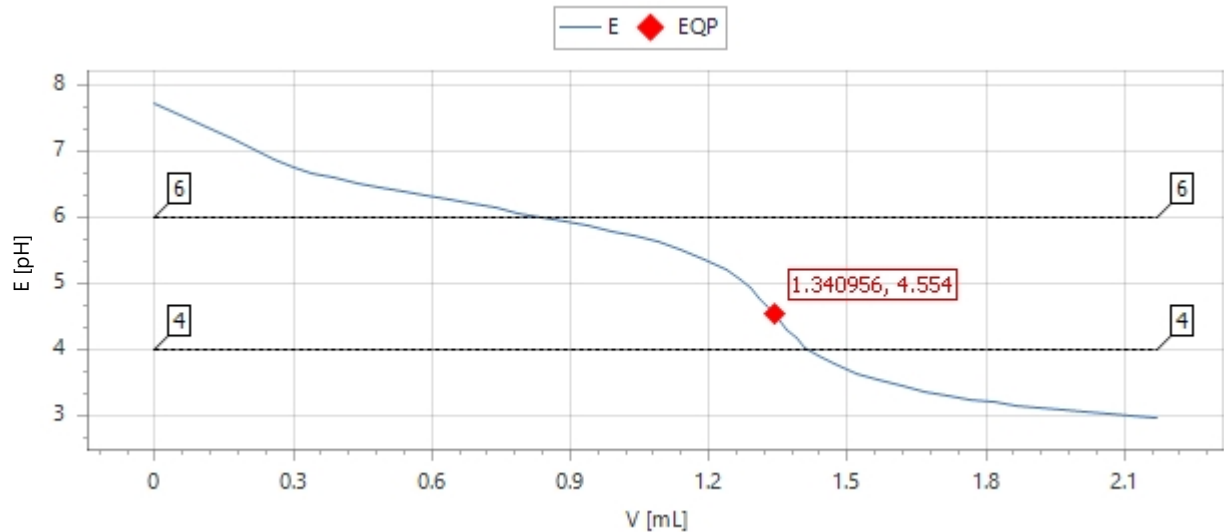

dE/dV - V Curve

Scope 1/1, Sample 18/30

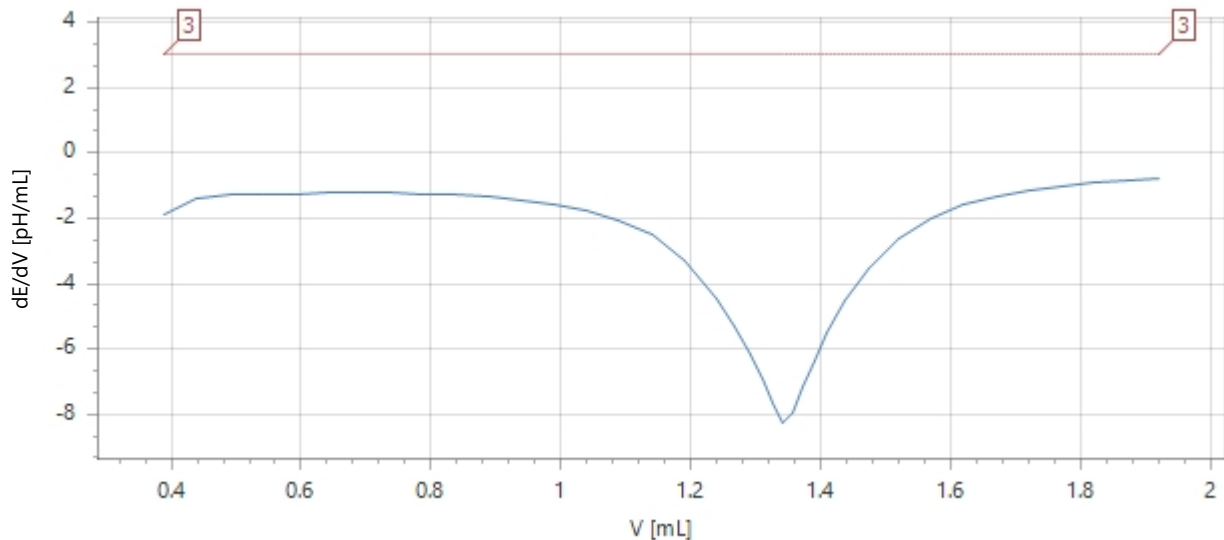

Measured Values

Scope 1/1, Sample 18/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.747		25.3			
28	0.17150	7.194	0.17150	25.3	-0.553		
58	0.25725	6.904	0.08575	25.3	-0.290		
88	0.30000	6.759	0.04275	25.3	-0.145		
180	0.33950	6.679	0.03950	25.2	-0.080		
210	0.38950	6.606	0.05000	25.2	-0.073	-1.88	
240	0.43950	6.535	0.05000	25.2	-0.071	-1.37	
270	0.48950	6.464	0.05000	25.2	-0.071	-1.25	

Titration report

LabX

301	0.53950	6.398	0.05000	25.2	-0.066	-1.27	EQP
331	0.58950	6.336	0.05000	25.2	-0.062	-1.26	
361	0.63950	6.274	0.05000	25.2	-0.062	-1.23	
391	0.68950	6.213	0.05000	25.2	-0.061	-1.22	
422	0.73950	6.151	0.05000	25.2	-0.062	-1.23	
452	0.78950	6.087	0.05000	25.2	-0.064	-1.26	
482	0.83950	6.025	0.05000	25.2	-0.062	-1.30	
512	0.88950	5.959	0.05000	25.2	-0.066	-1.36	
543	0.93950	5.888	0.05000	25.2	-0.071	-1.44	
573	0.98950	5.811	0.05000	25.2	-0.077	-1.56	
603	1.03950	5.728	0.05000	25.2	-0.083	-1.76	
633	1.08950	5.633	0.05000	25.2	-0.095	-2.05	
664	1.13950	5.522	0.05000	25.2	-0.111	-2.51	
694	1.18950	5.384	0.05000	25.2	-0.138	-3.27	
724	1.23800	5.210	0.04850	25.2	-0.174	-4.43	
754	1.26675	5.079	0.02875	25.2	-0.131	-5.34	
784	1.29025	4.946	0.02350	25.2	-0.133	-6.20	
815	1.30950	4.812	0.01925	25.2	-0.134	-6.93	
845	1.32550	4.685	0.01600	25.2	-0.127	-7.66	
875	1.34050	4.558	0.01500	25.2	-0.127	-8.24	
	1.340956	4.554					
905	1.35525	4.430	0.01475	25.2	-0.128	-7.94	
936	1.37050	4.306	0.01525	25.2	-0.124	-7.14	
966	1.38825	4.178	0.01775	25.2	-0.128	-6.35	
996	1.40925	4.053	0.02100	25.2	-0.125	-5.52	
1026	1.43700	3.920	0.02775	25.2	-0.133	-4.52	
1057	1.47225	3.789	0.03525	25.2	-0.131	-3.52	
1087	1.51900	3.656	0.04675	25.2	-0.133	-2.60	
1117	1.56900	3.547	0.05000	25.2	-0.109	-1.99	
1147	1.61900	3.459	0.05000	25.2	-0.088	-1.60	
1178	1.66900	3.385	0.05000	25.2	-0.074	-1.34	
1208	1.71900	3.321	0.05000	25.2	-0.064	-1.16	
1238	1.76900	3.266	0.05000	25.2	-0.055	-1.03	
1268	1.81900	3.217	0.05000	25.2	-0.049	-0.93	
1298	1.86900	3.173	0.05000	25.2	-0.044	-0.84	
1329	1.91900	3.132	0.05000	25.2	-0.041	-0.77	
1359	1.96900	3.095	0.05000	25.2	-0.037		
1389	2.01900	3.062	0.05000	25.2	-0.033		
1419	2.06900	3.030	0.05000	25.2	-0.032		
1450	2.11900	3.001	0.05000	25.2	-0.029		
1480	2.16900	2.973	0.05000	25.2	-0.028		

Measure (normal) [1] Scope 1/1, Sample 19/30

Measured Values Scope 1/1, Sample 19/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.731	22.7	
2	5.045	23.0	0.314
3	5.321	23.3	0.276

Titration report

LabX

4	5.538	23.7	0.217
5	5.722	24.0	0.184
6	5.887	24.3	0.165
7	6.034	24.5	0.147
8	6.154	24.7	0.120
9	6.275	24.8	0.121
10	6.387	24.9	0.112
11	6.491	25.0	0.104
12	6.592	25.0	0.101
13	6.684	25.1	0.092
14	6.771	25.1	0.087
15	6.847	25.1	0.076
16	6.922	25.2	0.075
17	6.991	25.2	0.069
18	7.049	25.2	0.058
19	7.113	25.2	0.064
20	7.167	25.2	0.054
21	7.216	25.2	0.049
22	7.257	25.2	0.041
23	7.301	25.2	0.044
24	7.341	25.2	0.040
25	7.378	25.2	0.037
26	7.413	25.2	0.035
27	7.445	25.2	0.032
28	7.472	25.2	0.027
29	7.496	25.2	0.024
30	7.521	25.2	0.025
31	7.543	25.2	0.022
32	7.564	25.2	0.021
33	7.584	25.2	0.020
34	7.603	25.2	0.019
35	7.620	25.2	0.017
36	7.635	25.2	0.015
37	7.649	25.2	0.014
38	7.661	25.2	0.012
39	7.673	25.2	0.012
40	7.684	25.2	0.011
41	7.695	25.2	0.011
42	7.704	25.2	0.009
43	7.712	25.2	0.008
44	7.721	25.2	0.009
45	7.728	25.2	0.007
46	7.735	25.2	0.007
47	7.741	25.2	0.006
48	7.747	25.2	0.006
49	7.752	25.2	0.005
50	7.757	25.2	0.005
51	7.761	25.2	0.004
52	7.766	25.2	0.005
53	7.769	25.2	0.003
54	7.773	25.2	0.004


Titration report

LabX


55	7.777	25.2	0.004
56	7.780	25.2	0.003
57	7.783	25.2	0.003
58	7.786	25.2	0.003
59	7.788	25.2	0.002
60	7.791	25.2	0.003

Titration (EQP) [1] Scope 1/1, Sample 19/30

E - V Curve Scope 1/1, Sample 19/30

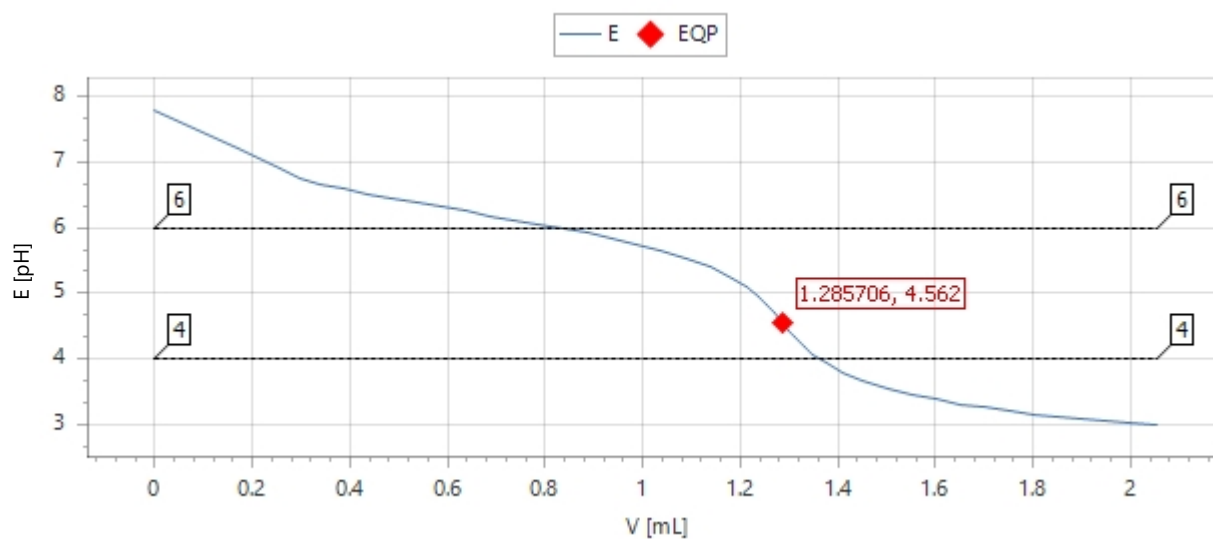

dE/dV - V Curve Scope 1/1, Sample 19/30

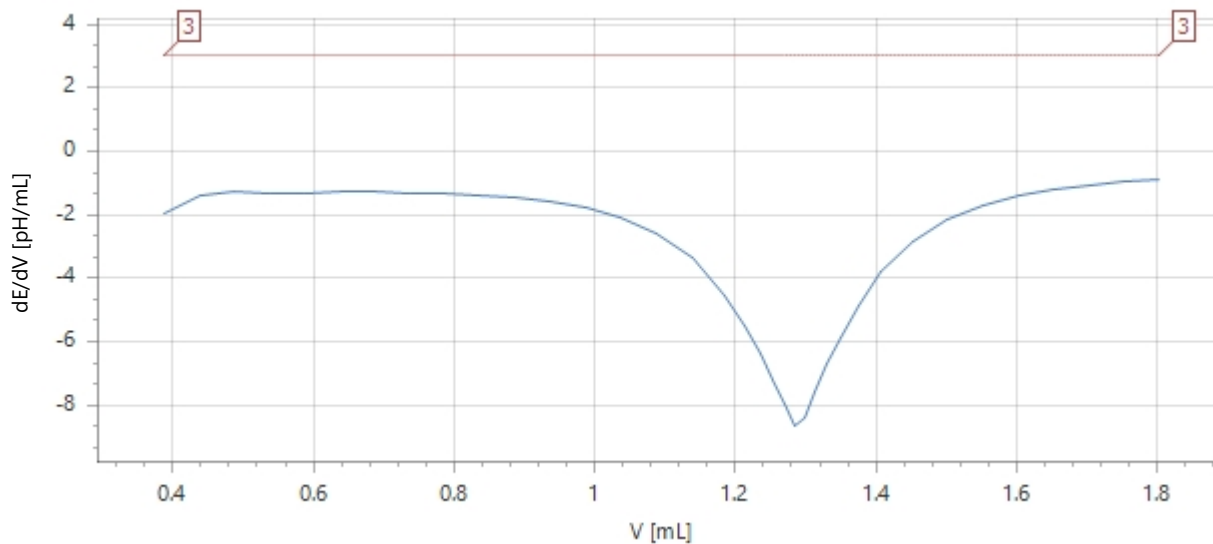

Measured Values Scope 1/1, Sample 19/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.796		25.2			
28	0.17150	7.214	0.17150	25.2	-0.582		
58	0.25725	6.907	0.08575	25.2	-0.307		

Titration report

LabX

88	0.30000	6.754	0.04275	25.2	-0.153	
179	0.33800	6.668	0.03800	25.2	-0.086	
209	0.38800	6.595	0.05000	25.2	-0.073	-1.98
239	0.43800	6.521	0.05000	25.2	-0.074	-1.41
269	0.48800	6.449	0.05000	25.2	-0.072	-1.29
300	0.53800	6.382	0.05000	25.2	-0.067	-1.32
330	0.58800	6.314	0.05000	25.2	-0.068	-1.31
360	0.63800	6.250	0.05000	25.2	-0.064	-1.29
390	0.68800	6.186	0.05000	25.2	-0.064	-1.29
421	0.73800	6.121	0.05000	25.2	-0.065	-1.30
451	0.78800	6.055	0.05000	25.2	-0.066	-1.34
481	0.83800	5.985	0.05000	25.2	-0.070	-1.39
511	0.88800	5.914	0.05000	25.2	-0.071	-1.48
542	0.93800	5.837	0.05000	25.2	-0.077	-1.60
572	0.98800	5.751	0.05000	25.2	-0.086	-1.80
602	1.03800	5.654	0.05000	25.2	-0.097	-2.10
632	1.08800	5.539	0.05000	25.2	-0.115	-2.59
663	1.13800	5.398	0.05000	25.2	-0.141	-3.37
693	1.18400	5.230	0.04600	25.2	-0.168	-4.52
723	1.21250	5.098	0.02850	25.2	-0.132	-5.49
753	1.23550	4.965	0.02300	25.2	-0.133	-6.39
784	1.25425	4.833	0.01875	25.2	-0.132	-7.17
814	1.27025	4.700	0.01600	25.2	-0.133	-7.96
844	1.28400	4.578	0.01375	25.2	-0.122	-8.62
	1.285706 ⁽⁹⁾	4.562 ⁽⁹⁾				EQP
874	1.29825	4.446	0.01425	25.2	-0.132	-8.39
904	1.31250	4.325	0.01425	25.2	-0.121	-7.58
935	1.33000	4.192	0.01750	25.2	-0.133	-6.72
965	1.34975	4.066	0.01975	25.2	-0.126	-5.86
995	1.37525	3.935	0.02550	25.2	-0.131	-4.85
1025	1.40800	3.804	0.03275	25.2	-0.131	-3.81
1056	1.45200	3.670	0.04400	25.2	-0.134	-2.83
1086	1.50200	3.553	0.05000	25.2	-0.117	-2.12
1116	1.55200	3.460	0.05000	25.2	-0.093	-1.68
1146	1.60200	3.383	0.05000	25.2	-0.077	-1.40
1177	1.65200	3.317	0.05000	25.2	-0.066	-1.20
1207	1.70200	3.260	0.05000	25.2	-0.057	-1.06
1237	1.75200	3.209	0.05000	25.2	-0.051	-0.95
1267	1.80200	3.164	0.05000	25.2	-0.045	-0.86
1298	1.85200	3.123	0.05000	25.2	-0.041	
1328	1.90200	3.085	0.05000	25.2	-0.038	
1358	1.95200	3.051	0.05000	25.2	-0.034	
1388	2.00200	3.018	0.05000	25.1	-0.033	
1419	2.05200	2.989	0.05000	25.1	-0.029	

Measure (normal) [1] Scope 1/1, Sample 20/30

Measured Values Scope 1/1, Sample 20/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
-------------	---------------------	---------------	------------

Titration report

LabX

1	4.727	22.9	
2	5.030	23.1	0.303
3	5.329	23.5	0.299
4	5.565	23.8	0.236
5	5.766	24.1	0.201
6	5.942	24.4	0.176
7	6.109	24.6	0.167
8	6.248	24.7	0.139
9	6.393	24.9	0.145
10	6.526	24.9	0.133
11	6.650	25.0	0.124
12	6.765	25.1	0.115
13	6.873	25.1	0.108
14	6.966	25.2	0.093
15	7.043	25.2	0.077
16	7.128	25.2	0.085
17	7.201	25.2	0.073
18	7.268	25.2	0.067
19	7.330	25.3	0.062
20	7.387	25.3	0.057
21	7.436	25.3	0.049
22	7.476	25.3	0.040
23	7.516	25.3	0.040
24	7.556	25.3	0.040
25	7.590	25.3	0.034
26	7.622	25.3	0.032
27	7.649	25.3	0.027
28	7.675	25.3	0.026
29	7.696	25.3	0.021
30	7.718	25.3	0.022
31	7.738	25.3	0.020
32	7.757	25.3	0.019
33	7.772	25.3	0.015
34	7.787	25.3	0.015
35	7.800	25.3	0.013
36	7.811	25.3	0.011
37	7.821	25.3	0.010
38	7.832	25.3	0.011
39	7.841	25.3	0.009
40	7.849	25.3	0.008
41	7.856	25.3	0.007
42	7.863	25.3	0.007
43	7.869	25.3	0.006
44	7.875	25.3	0.006
45	7.880	25.3	0.005
46	7.885	25.3	0.005
47	7.889	25.3	0.004
48	7.893	25.3	0.004
49	7.897	25.3	0.004
50	7.900	25.3	0.003
51	7.903	25.3	0.003


Titration report

LabX


52	7.906	25.3	0.003
53	7.909	25.3	0.003
54	7.912	25.3	0.003
55	7.914	25.3	0.002
56	7.916	25.3	0.002
57	7.918	25.3	0.002
58	7.920	25.3	0.002
59	7.921	25.3	0.001
60	7.922	25.3	0.001

Titration (EQP) [1] Scope 1/1, Sample 20/30

E - V Curve Scope 1/1, Sample 20/30

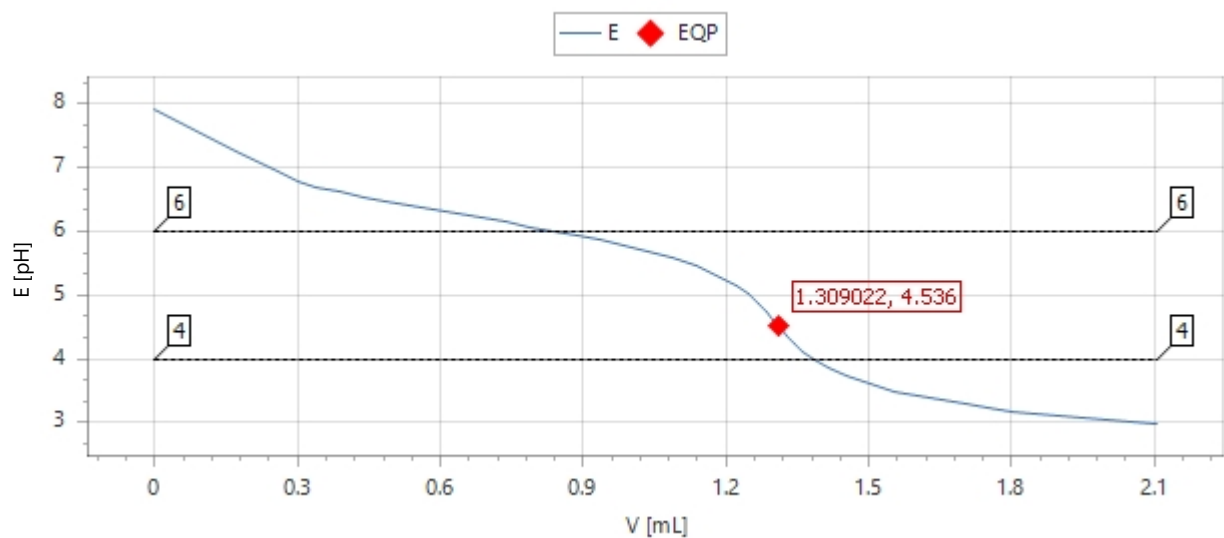

dE/dV - V Curve Scope 1/1, Sample 20/30

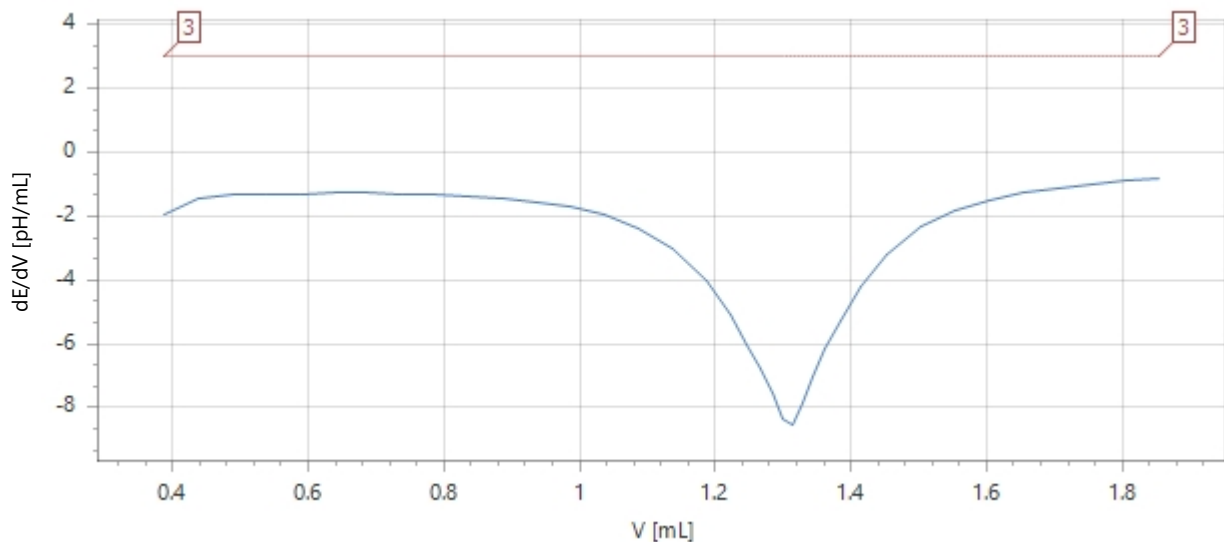

Measured Values Scope 1/1, Sample 20/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
----------	--------	------------------	---------	------------	---------	---------------	-----

Titration report

LabX

0	0.00000	7.926		25.3		
30	0.17150	7.244	0.17150	25.3	-0.682	
60	0.25725	6.939	0.08575	25.3	-0.305	
91	0.30000	6.788	0.04275	25.3	-0.151	
182	0.33800	6.699	0.03800	25.3	-0.089	
212	0.38800	6.623	0.05000	25.3	-0.076	-1.96
242	0.43800	6.546	0.05000	25.3	-0.077	-1.46
272	0.48800	6.472	0.05000	25.3	-0.074	-1.32
303	0.53800	6.404	0.05000	25.3	-0.068	-1.33
333	0.58800	6.336	0.05000	25.3	-0.068	-1.32
363	0.63800	6.272	0.05000	25.3	-0.064	-1.28
393	0.68800	6.208	0.05000	25.3	-0.064	-1.28
424	0.73800	6.144	0.05000	25.3	-0.064	-1.29
454	0.78800	6.078	0.05000	25.3	-0.066	-1.32
484	0.83800	6.011	0.05000	25.3	-0.067	-1.37
514	0.88800	5.940	0.05000	25.3	-0.071	-1.45
545	0.93800	5.866	0.05000	25.3	-0.074	-1.55
575	0.98800	5.781	0.05000	25.3	-0.085	-1.71
605	1.03800	5.690	0.05000	25.3	-0.091	-1.97
635	1.08800	5.584	0.05000	25.3	-0.106	-2.36
666	1.13800	5.456	0.05000	25.3	-0.128	-3.00
696	1.18800	5.290	0.05000	25.3	-0.166	-4.04
726	1.22150	5.151	0.03350	25.3	-0.139	-5.06
756	1.24750	5.013	0.02600	25.3	-0.138	-6.02
786	1.26725	4.883	0.01975	25.3	-0.130	-6.81
817	1.28400	4.754	0.01675	25.3	-0.129	-7.55
847	1.29900	4.625	0.01500	25.3	-0.129	-8.37
	1.309022 ⁽⁹⁾	4.536 ⁽⁹⁾				EQP
877	1.31300	4.500	0.01400	25.3	-0.125	-8.52
908	1.32750	4.371	0.01450	25.3	-0.129	-7.87
938	1.34300	4.248	0.01550	25.3	-0.123	-7.02
968	1.36200	4.116	0.01900	25.3	-0.132	-6.18
998	1.38450	3.990	0.02250	25.3	-0.126	-5.26
1028	1.41450	3.859	0.03000	25.3	-0.131	-4.21
1059	1.45375	3.725	0.03925	25.3	-0.134	-3.20
1089	1.50375	3.596	0.05000	25.3	-0.129	-2.34
1119	1.55375	3.495	0.05000	25.3	-0.101	-1.82
1149	1.60375	3.413	0.05000	25.3	-0.082	-1.49
1180	1.65375	3.344	0.05000	25.3	-0.069	-1.26
1210	1.70375	3.284	0.05000	25.3	-0.060	-1.11
1240	1.75375	3.231	0.05000	25.3	-0.053	-0.99
1270	1.80375	3.183	0.05000	25.3	-0.048	-0.89
1301	1.85375	3.141	0.05000	25.3	-0.042	-0.81
1331	1.90375	3.102	0.05000	25.3	-0.039	
1361	1.95375	3.067	0.05000	25.3	-0.035	
1391	2.00375	3.034	0.05000	25.3	-0.033	
1422	2.05375	3.003	0.05000	25.3	-0.031	
1452	2.10375	2.974	0.05000	25.3	-0.029	

Measure (normal) [1] Scope 1/1, Sample 21/30

Titration report

LabX

Measured Values Scope 1/1, Sample 21/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.736	23.2	
2	5.028	23.4	0.292
3	5.302	23.7	0.274
4	5.532	24.0	0.230
5	5.733	24.2	0.201
6	5.912	24.5	0.179
7	6.073	24.7	0.161
8	6.224	24.8	0.151
9	6.351	24.9	0.127
10	6.482	25.0	0.131
11	6.602	25.1	0.120
12	6.714	25.2	0.112
13	6.817	25.2	0.103
14	6.912	25.2	0.095
15	6.994	25.3	0.082
16	7.062	25.3	0.068
17	7.141	25.3	0.079
18	7.209	25.3	0.068
19	7.271	25.3	0.062
20	7.329	25.3	0.058
21	7.382	25.3	0.053
22	7.425	25.3	0.043
23	7.468	25.4	0.043
24	7.509	25.4	0.041
25	7.545	25.4	0.036
26	7.579	25.4	0.034
27	7.609	25.4	0.030
28	7.638	25.4	0.029
29	7.660	25.4	0.022
30	7.685	25.4	0.025
31	7.707	25.4	0.022
32	7.728	25.4	0.021
33	7.746	25.4	0.018
34	7.763	25.4	0.017
35	7.778	25.4	0.015
36	7.791	25.4	0.013
37	7.804	25.4	0.013
38	7.817	25.4	0.013
39	7.827	25.4	0.010
40	7.837	25.4	0.010
41	7.846	25.4	0.009
42	7.855	25.4	0.009
43	7.862	25.4	0.007
44	7.869	25.4	0.007
45	7.876	25.4	0.007
46	7.882	25.4	0.006
47	7.887	25.4	0.005


Titration report

LabX


48	7.892	25.4	0.005
49	7.896	25.4	0.004
50	7.900	25.4	0.004
51	7.905	25.4	0.005
52	7.908	25.4	0.003
53	7.911	25.4	0.003
54	7.915	25.4	0.004
55	7.917	25.4	0.002
56	7.920	25.4	0.003
57	7.922	25.4	0.002
58	7.924	25.4	0.002
59	7.927	25.4	0.003
60	7.929	25.4	0.002

Titration (EQP) [1] Scope 1/1, Sample 21/30

E - V Curve Scope 1/1, Sample 21/30

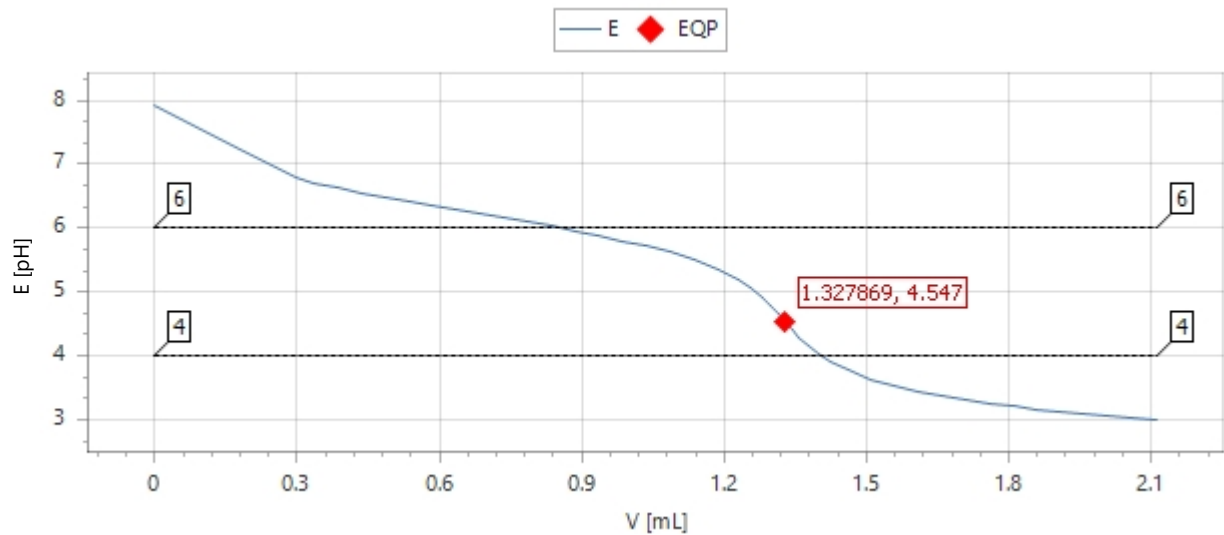

dE/dV - V Curve Scope 1/1, Sample 21/30

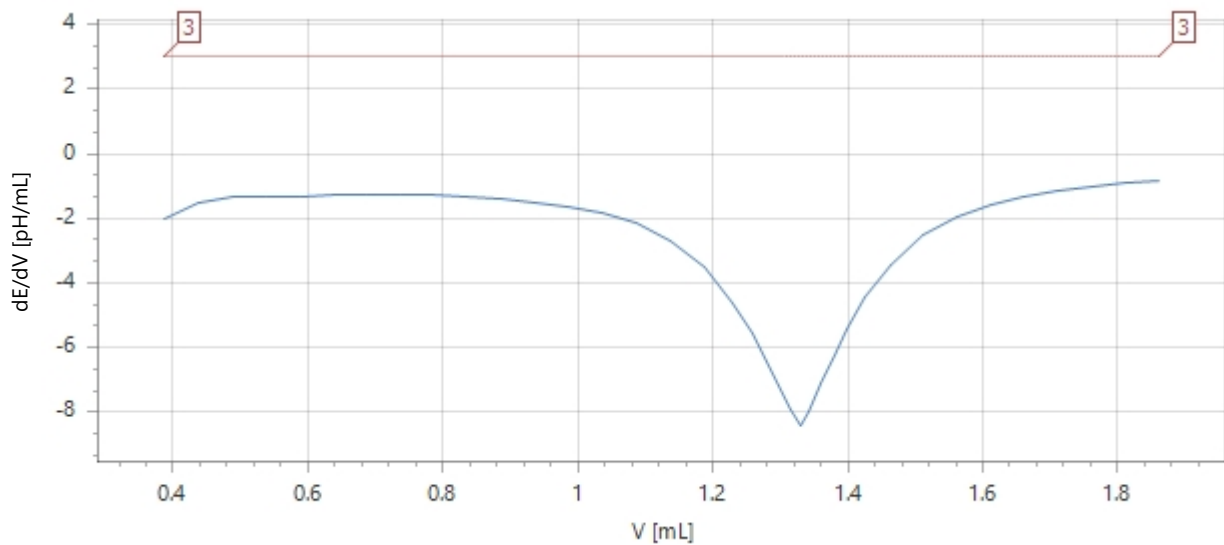

Titration report

LabX

Measured Values Scope 1/1, Sample 21/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.933		25.4			
30	0.17150	7.276	0.17150	25.4	-0.657		
60	0.25725	6.965	0.08575	25.3	-0.311		
90	0.30000	6.809	0.04275	25.3	-0.156		
181	0.33700	6.710	0.03700	25.4	-0.099		
211	0.38700	6.636	0.05000	25.3	-0.074	-2.04	
241	0.43700	6.559	0.05000	25.3	-0.077	-1.49	
271	0.48700	6.485	0.05000	25.3	-0.074	-1.31	
302	0.53700	6.416	0.05000	25.3	-0.069	-1.33	
332	0.58700	6.349	0.05000	25.3	-0.067	-1.32	
362	0.63700	6.284	0.05000	25.3	-0.065	-1.28	
392	0.68700	6.220	0.05000	25.3	-0.064	-1.27	
423	0.73700	6.157	0.05000	25.3	-0.063	-1.27	
453	0.78700	6.093	0.05000	25.3	-0.064	-1.30	
483	0.83700	6.027	0.05000	25.3	-0.066	-1.34	
513	0.88700	5.957	0.05000	25.3	-0.070	-1.41	
544	0.93700	5.884	0.05000	25.3	-0.073	-1.50	
574	0.98700	5.805	0.05000	25.3	-0.079	-1.63	
604	1.03700	5.717	0.05000	25.3	-0.088	-1.84	
634	1.08700	5.619	0.05000	25.3	-0.098	-2.15	
664	1.13700	5.502	0.05000	25.3	-0.117	-2.67	
695	1.18700	5.356	0.05000	25.3	-0.146	-3.52	
725	1.23025	5.191	0.04325	25.3	-0.165	-4.64	
755	1.25750	5.061	0.02725	25.3	-0.130	-5.56	
786	1.28025	4.925	0.02275	25.3	-0.136	-6.42	
816	1.29850	4.794	0.01825	25.3	-0.131	-7.15	
846	1.31400	4.668	0.01550	25.3	-0.126	-7.91	
	1.327869 ⁽⁹⁾	4.547 ⁽⁹⁾					EQP
876	1.32875	4.539	0.01475	25.3	-0.129	-8.40	
906	1.34300	4.414	0.01425	25.3	-0.125	-7.98	
937	1.35875	4.284	0.01575	25.3	-0.130	-7.13	
967	1.37625	4.160	0.01750	25.3	-0.124	-6.35	
997	1.39825	4.030	0.02200	25.3	-0.130	-5.46	
1027	1.42600	3.901	0.02775	25.3	-0.129	-4.46	
1058	1.46300	3.766	0.03700	25.3	-0.135	-3.43	
1088	1.51075	3.634	0.04775	25.3	-0.132	-2.53	
1118	1.56075	3.527	0.05000	25.3	-0.107	-1.94	
1148	1.61075	3.441	0.05000	25.3	-0.086	-1.57	
1179	1.66075	3.368	0.05000	25.3	-0.073	-1.32	
1209	1.71075	3.306	0.05000	25.3	-0.062	-1.15	
1239	1.76075	3.251	0.05000	25.3	-0.055	-1.02	
1269	1.81075	3.202	0.05000	25.3	-0.049	-0.92	
1300	1.86075	3.159	0.05000	25.3	-0.043	-0.83	
1330	1.91075	3.119	0.05000	25.3	-0.040		
1360	1.96075	3.082	0.05000	25.3	-0.037		
1390	2.01075	3.049	0.05000	25.3	-0.033		
1421	2.06075	3.018	0.05000	25.3	-0.031		

Titration report

LabX

1451 2.11075 2.989 0.05000 25.3 -0.029

Measure (normal) [1] Scope 1/1, Sample 22/30

Measured Values Scope 1/1, Sample 22/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.416	23.5	
2	4.796	23.7	0.380
3	5.223	23.9	0.427
4	5.601	24.2	0.378
5	5.940	24.4	0.339
6	6.240	24.6	0.300
7	6.493	24.8	0.253
8	6.687	24.9	0.194
9	6.875	25.0	0.188
10	7.033	25.1	0.158
11	7.172	25.1	0.139
12	7.287	25.2	0.115
13	7.382	25.2	0.095
14	7.462	25.3	0.080
15	7.521	25.3	0.059
16	7.578	25.3	0.057
17	7.626	25.3	0.048
18	7.666	25.3	0.040
19	7.701	25.3	0.035
20	7.730	25.3	0.029
21	7.754	25.3	0.024
22	7.773	25.3	0.019
23	7.790	25.4	0.017
24	7.804	25.4	0.014
25	7.817	25.4	0.013
26	7.827	25.4	0.010
27	7.836	25.4	0.009
28	7.844	25.4	0.008
29	7.851	25.4	0.007
30	7.856	25.4	0.005
31	7.860	25.4	0.004
32	7.864	25.4	0.004
33	7.868	25.4	0.004
34	7.872	25.4	0.004
35	7.875	25.4	0.003
36	7.877	25.4	0.002
37	7.879	25.4	0.002
38	7.881	25.4	0.002
39	7.882	25.4	0.001
40	7.884	25.4	0.002
41	7.885	25.4	0.001
42	7.886	25.4	0.001
43	7.888	25.4	0.002


Titration report

LabX

44	7.889	25.4	0.001
45	7.890	25.4	0.001
46	7.891	25.4	0.001
47	7.891	25.4	0.000
48	7.892	25.4	0.001
49	7.893	25.4	0.001
50	7.894	25.4	0.001
51	7.894	25.4	0.000
52	7.895	25.4	0.001
53	7.895	25.4	0.000
54	7.896	25.4	0.001
55	7.897	25.4	0.001
56	7.897	25.4	0.000
57	7.898	25.4	0.001
58	7.898	25.4	0.000
59	7.898	25.4	0.000
60	7.899	25.4	0.001

Titration (EQP) [1] Scope 1/1, Sample 22/30


E - V Curve Scope 1/1, Sample 22/30

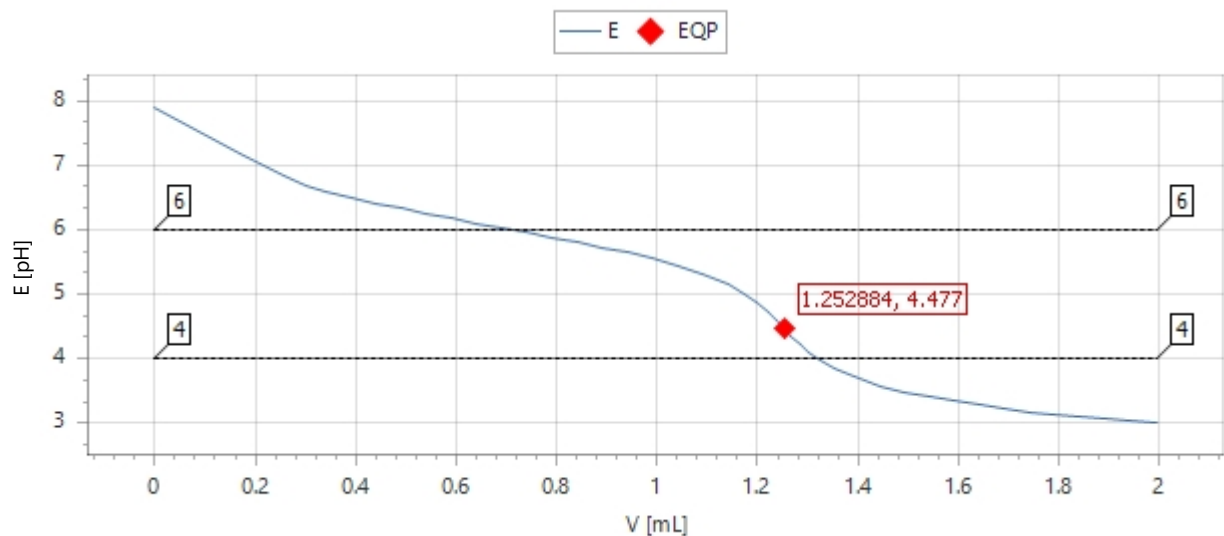

dE/dV - V Curve Scope 1/1, Sample 22/30

Titration report

LabX

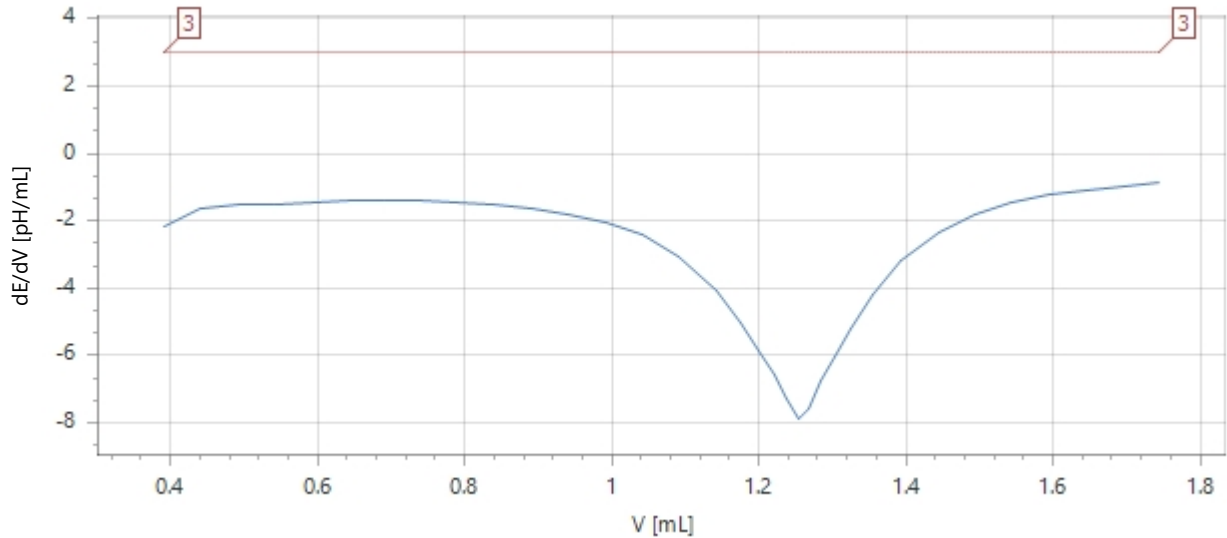

Measured Values Scope 1/1, Sample 22/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.900		25.4			
28	0.17150	7.176	0.17150	25.4	-0.724		
58	0.25725	6.840	0.08575	25.3	-0.336		
88	0.30000	6.688	0.04275	25.3	-0.152		
179	0.34175	6.586	0.04175	25.3	-0.102		
209	0.39175	6.496	0.05000	25.3	-0.090	-2.16	
239	0.44175	6.406	0.05000	25.3	-0.090	-1.66	
269	0.49175	6.321	0.05000	25.3	-0.085	-1.54	
300	0.54175	6.244	0.05000	25.3	-0.077	-1.52	
330	0.59175	6.169	0.05000	25.3	-0.075	-1.47	
360	0.64175	6.096	0.05000	25.3	-0.073	-1.42	
390	0.69175	6.025	0.05000	25.3	-0.071	-1.41	
421	0.74175	5.956	0.05000	25.3	-0.069	-1.42	
451	0.79175	5.884	0.05000	25.3	-0.072	-1.46	
481	0.84175	5.808	0.05000	25.3	-0.076	-1.52	
511	0.89175	5.728	0.05000	25.3	-0.080	-1.63	
542	0.94175	5.643	0.05000	25.3	-0.085	-1.79	
572	0.99175	5.547	0.05000	25.3	-0.096	-2.03	
602	1.04175	5.438	0.05000	25.3	-0.109	-2.42	
632	1.09175	5.306	0.05000	25.3	-0.132	-3.05	
662	1.14175	5.140	0.05000	25.3	-0.166	-4.08	
693	1.17475	4.999	0.03300	25.3	-0.141	-5.02	
723	1.19975	4.867	0.02500	25.3	-0.132	-5.86	
753	1.22025	4.733	0.02050	25.3	-0.134	-6.58	
784	1.23700	4.609	0.01675	25.3	-0.124	-7.28	
	1.252884 ⁽⁹⁾	4.477 ⁽⁹⁾					EQP
814	1.25325	4.474	0.01625	25.3	-0.135	-7.87	
844	1.26775	4.353	0.01450	25.3	-0.121	-7.59	
874	1.28400	4.224	0.01625	25.3	-0.129	-6.78	
904	1.30200	4.101	0.01800	25.3	-0.123	-6.04	
935	1.32500	3.972	0.02300	25.3	-0.129	-5.18	

Titration report

LabX

965	1.35475	3.841	0.02975	25.3	-0.131	-4.18
995	1.39375	3.710	0.03900	25.3	-0.131	-3.18
1025	1.44375	3.582	0.05000	25.3	-0.128	-2.32
1056	1.49375	3.483	0.05000	25.3	-0.099	-1.79
1086	1.54375	3.402	0.05000	25.3	-0.081	-1.46
1116	1.59375	3.334	0.05000	25.3	-0.068	-1.24
1146	1.64375	3.275	0.05000	25.3	-0.059	-1.08
1177	1.69375	3.223	0.05000	25.3	-0.052	-0.97
1207	1.74375	3.177	0.05000	25.3	-0.046	-0.88
1237	1.79375	3.135	0.05000	25.3	-0.042	
1267	1.84375	3.097	0.05000	25.3	-0.038	
1298	1.89375	3.061	0.05000	25.3	-0.036	
1328	1.94375	3.029	0.05000	25.3	-0.032	
1358	1.99375	2.999	0.05000	25.3	-0.030	

Measure (normal) [1] Scope 1/1, Sample 23/30

Measured Values Scope 1/1, Sample 23/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.945	23.5	
2	5.214	23.7	0.269
3	5.491	23.9	0.277
4	5.724	24.2	0.233
5	5.926	24.4	0.202
6	6.105	24.6	0.179
7	6.269	24.8	0.164
8	6.409	24.9	0.140
9	6.553	25.0	0.144
10	6.689	25.0	0.136
11	6.813	25.1	0.124
12	6.924	25.1	0.111
13	7.021	25.2	0.097
14	7.117	25.2	0.096
15	7.193	25.2	0.076
16	7.269	25.2	0.076
17	7.343	25.2	0.074
18	7.409	25.3	0.066
19	7.466	25.3	0.057
20	7.518	25.3	0.052
21	7.567	25.3	0.049
22	7.607	25.3	0.040
23	7.646	25.3	0.039
24	7.682	25.3	0.036
25	7.714	25.3	0.032
26	7.743	25.3	0.029
27	7.768	25.3	0.025
28	7.792	25.3	0.024
29	7.811	25.3	0.019
30	7.831	25.3	0.020


Titration report

LabX

31	7.848	25.3	0.017
32	7.864	25.3	0.016
33	7.879	25.3	0.015
34	7.892	25.3	0.013
35	7.903	25.3	0.011
36	7.913	25.3	0.010
37	7.923	25.3	0.010
38	7.932	25.3	0.009
39	7.940	25.3	0.008
40	7.948	25.3	0.008
41	7.954	25.3	0.006
42	7.961	25.3	0.007
43	7.966	25.3	0.005
44	7.971	25.3	0.005
45	7.975	25.3	0.004
46	7.979	25.3	0.004
47	7.983	25.3	0.004
48	7.987	25.3	0.004
49	7.990	25.3	0.003
50	7.992	25.3	0.002
51	7.995	25.3	0.003
52	7.997	25.3	0.002
53	7.999	25.3	0.002
54	8.001	25.3	0.002
55	8.003	25.3	0.002
56	8.005	25.3	0.002
57	8.006	25.3	0.001
58	8.007	25.3	0.001
59	8.009	25.3	0.002
60	8.010	25.3	0.001

Titration (EQP) [1] Scope 1/1, Sample 23/30


E - V Curve Scope 1/1, Sample 23/30

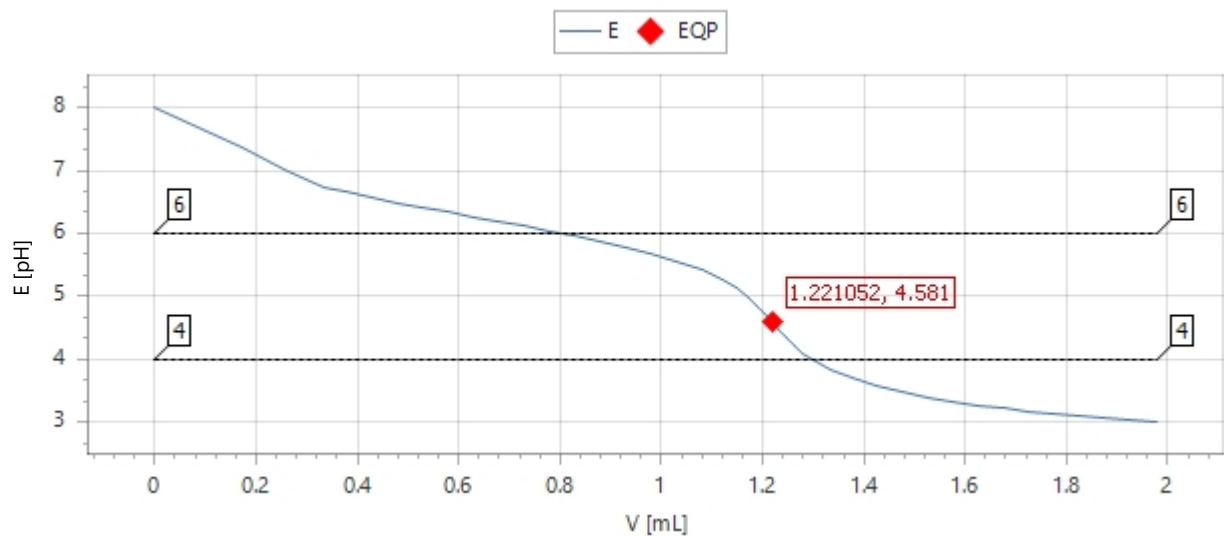

Titration report

LabX

dE/dV - V Curve Scope 1/1, Sample 23/30

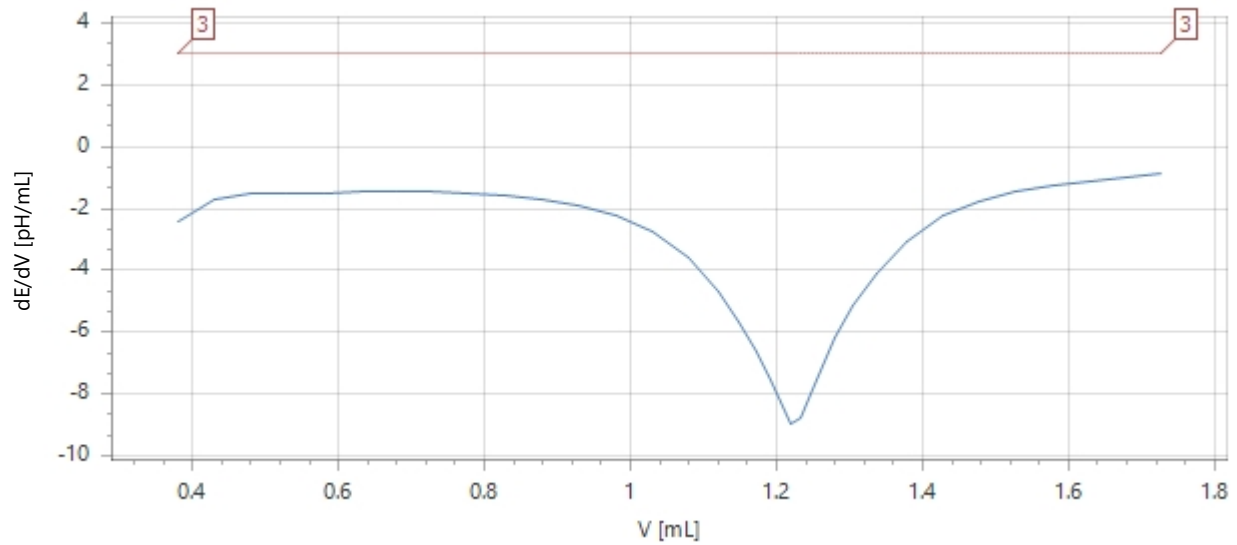

Measured Values Scope 1/1, Sample 23/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	8.012		25.3			
28	0.17150	7.364	0.17150	25.3	-0.648		
58	0.25725	7.019	0.08575	25.3	-0.345		
88	0.30000	6.850	0.04275	25.3	-0.169		
179	0.33450	6.731	0.03450	25.3	-0.119		
209	0.37975	6.656	0.04525	25.3	-0.075	-2.41	
239	0.42975	6.570	0.05000	25.3	-0.086	-1.73	
270	0.47975	6.489	0.05000	25.3	-0.081	-1.49	
300	0.52975	6.409	0.05000	25.2	-0.080	-1.48	
330	0.57975	6.334	0.05000	25.2	-0.075	-1.51	
360	0.62975	6.259	0.05000	25.3	-0.075	-1.47	
390	0.67975	6.186	0.05000	25.3	-0.073	-1.46	
421	0.72975	6.114	0.05000	25.3	-0.072	-1.47	
451	0.77975	6.038	0.05000	25.2	-0.076	-1.52	
481	0.82975	5.961	0.05000	25.2	-0.077	-1.58	
511	0.87975	5.877	0.05000	25.2	-0.084	-1.70	
542	0.92975	5.787	0.05000	25.2	-0.090	-1.90	
572	0.97975	5.684	0.05000	25.2	-0.103	-2.21	
602	1.02975	5.566	0.05000	25.2	-0.118	-2.72	
632	1.07975	5.417	0.05000	25.2	-0.149	-3.57	
663	1.12125	5.258	0.04150	25.2	-0.159	-4.67	
693	1.14875	5.127	0.02750	25.2	-0.131	-5.65	
723	1.17175	4.990	0.02300	25.2	-0.137	-6.61	
753	1.18975	4.858	0.01800	25.2	-0.132	-7.41	
784	1.20500	4.730	0.01525	25.2	-0.128	-8.23	
814	1.21900	4.601	0.01400	25.2	-0.129	-8.94	
	1.221052	4.581					EQP
844	1.23225	4.475	0.01325	25.2	-0.126	-8.79	
874	1.24625	4.350	0.01400	25.2	-0.125	-7.97	

Titration report

LabX

905	1.26250	4.218	0.01625	25.2	-0.132	-7.06
935	1.28075	4.094	0.01825	25.2	-0.124	-6.20
965	1.30475	3.964	0.02400	25.2	-0.130	-5.16
995	1.33575	3.832	0.03100	25.2	-0.132	-4.08
1026	1.37700	3.699	0.04125	25.2	-0.133	-3.05
1056	1.42700	3.574	0.05000	25.2	-0.125	-2.25
1086	1.47700	3.477	0.05000	25.2	-0.097	-1.76
1116	1.52700	3.397	0.05000	25.2	-0.080	-1.45
1146	1.57700	3.330	0.05000	25.2	-0.067	-1.24
1177	1.62700	3.270	0.05000	25.2	-0.060	-1.09
1207	1.67700	3.219	0.05000	25.2	-0.051	-0.97
1237	1.72700	3.172	0.05000	25.2	-0.047	-0.88
1268	1.77700	3.130	0.05000	25.2	-0.042	
1298	1.82700	3.092	0.05000	25.2	-0.038	
1328	1.87700	3.056	0.05000	25.2	-0.036	
1358	1.92700	3.024	0.05000	25.2	-0.032	
1388	1.97700	2.994	0.05000	25.2	-0.030	

Measure (normal) [1] Scope 1/1, Sample 24/30

Measured Values Scope 1/1, Sample 24/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.776	23.2	
2	5.068	23.4	0.292
3	5.349	23.7	0.281
4	5.584	24.0	0.235
5	5.788	24.3	0.204
6	5.970	24.5	0.182
7	6.139	24.6	0.169
8	6.279	24.8	0.140
9	6.425	24.9	0.146
10	6.558	24.9	0.133
11	6.679	25.0	0.121
12	6.793	25.1	0.114
13	6.902	25.1	0.109
14	6.995	25.1	0.093
15	7.082	25.1	0.087
16	7.155	25.2	0.073
17	7.227	25.2	0.072
18	7.294	25.2	0.067
19	7.356	25.2	0.062
20	7.414	25.2	0.058
21	7.465	25.2	0.051
22	7.507	25.2	0.042
23	7.550	25.2	0.043
24	7.590	25.2	0.040
25	7.625	25.2	0.035
26	7.658	25.2	0.033
27	7.689	25.2	0.031

Titration report

LabX


28	7.715	25.2	0.026
29	7.736	25.2	0.021
30	7.758	25.2	0.022
31	7.778	25.2	0.020
32	7.797	25.2	0.019
33	7.814	25.2	0.017
34	7.829	25.2	0.015
35	7.842	25.2	0.013
36	7.853	25.2	0.011
37	7.864	25.2	0.011
38	7.875	25.2	0.011
39	7.886	25.2	0.011
40	7.895	25.2	0.009
41	7.902	25.2	0.007
42	7.910	25.2	0.008
43	7.916	25.2	0.006
44	7.921	25.2	0.005
45	7.927	25.2	0.006
46	7.932	25.2	0.005
47	7.937	25.2	0.005
48	7.941	25.2	0.004
49	7.945	25.2	0.004
50	7.948	25.2	0.003
51	7.951	25.2	0.003
52	7.954	25.2	0.003
53	7.957	25.2	0.003
54	7.960	25.2	0.003
55	7.962	25.2	0.002
56	7.964	25.2	0.002
57	7.966	25.2	0.002
58	7.968	25.2	0.002
59	7.969	25.2	0.001
60	7.971	25.2	0.002

Titration (EQP) [1] Scope 1/1, Sample 24/30


E - V Curve Scope 1/1, Sample 24/30

Titration report

LabX

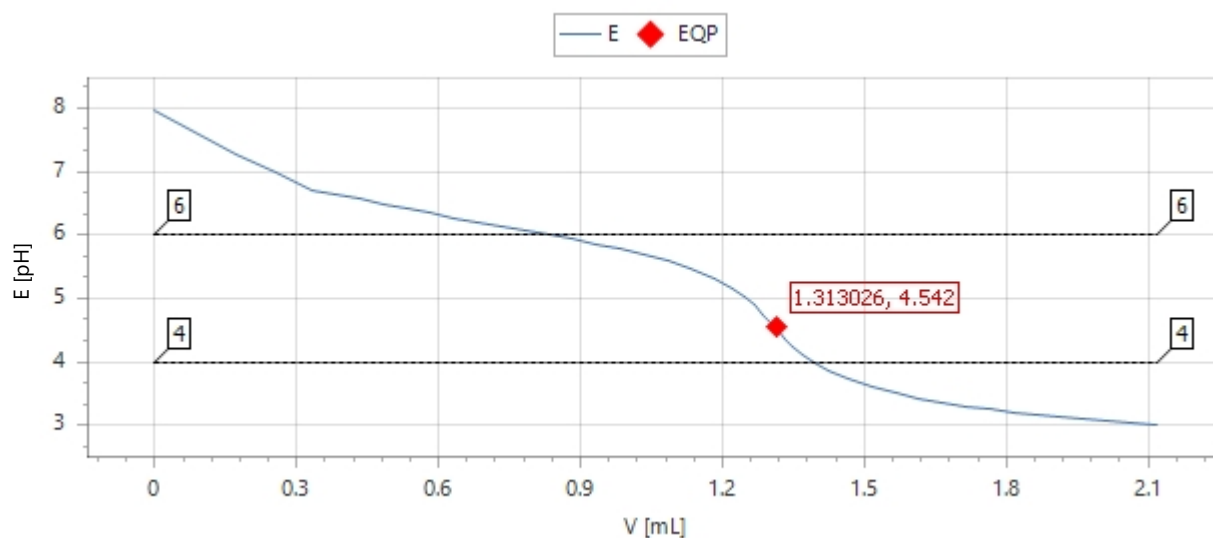

dE/dV - V Curve Scope 1/1, Sample 24/30

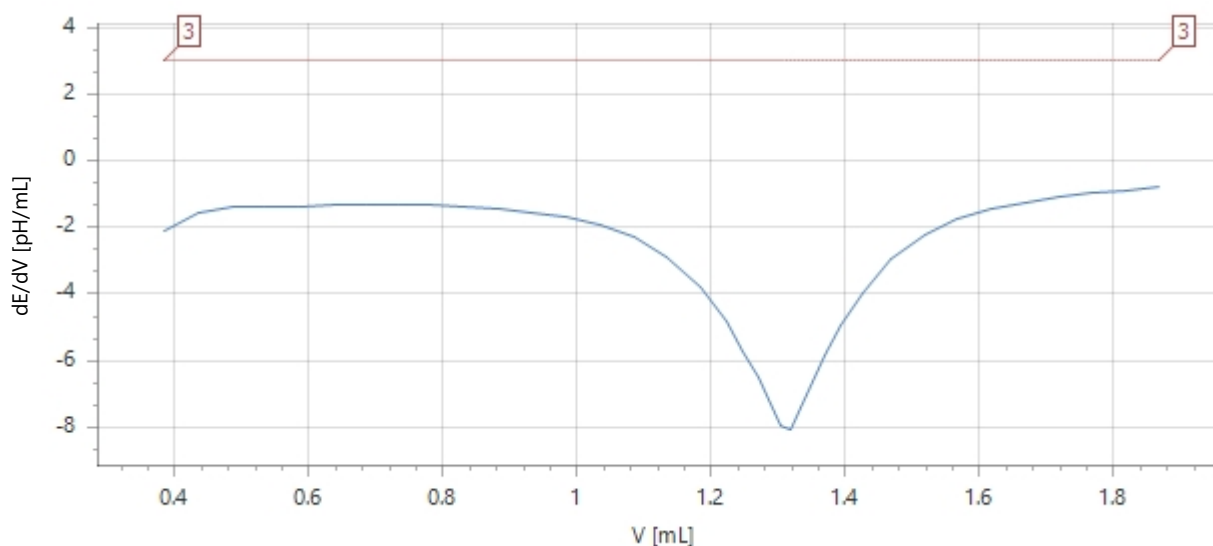

Measured Values Scope 1/1, Sample 24/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.974		25.2			
31	0.17150	7.285	0.17150	25.2	-0.689		
62	0.25725	6.980	0.08575	25.2	-0.305		
92	0.30000	6.821	0.04275	25.2	-0.159		
183	0.33450	6.714	0.03450	25.2	-0.107		
213	0.38450	6.639	0.05000	25.2	-0.075	-2.11	
243	0.43450	6.561	0.05000	25.2	-0.078	-1.58	
273	0.48450	6.485	0.05000	25.2	-0.076	-1.36	
304	0.53450	6.411	0.05000	25.2	-0.074	-1.36	
334	0.58450	6.340	0.05000	25.2	-0.071	-1.38	
364	0.63450	6.276	0.05000	25.2	-0.064	-1.33	
394	0.68450	6.210	0.05000	25.2	-0.066	-1.31	
424	0.73450	6.144	0.05000	25.1	-0.066	-1.31	
455	0.78450	6.074	0.05000	25.1	-0.070	-1.34	

Titration report

LabX

485	0.83450	6.010	0.05000	25.1	-0.064	-1.38	EQP
515	0.88450	5.939	0.05000	25.1	-0.071	-1.44	
546	0.93450	5.862	0.05000	25.1	-0.077	-1.54	
576	0.98450	5.781	0.05000	25.1	-0.081	-1.69	
606	1.03450	5.690	0.05000	25.1	-0.091	-1.94	
636	1.08450	5.586	0.05000	25.1	-0.104	-2.29	
666	1.13450	5.461	0.05000	25.1	-0.125	-2.88	
697	1.18450	5.303	0.05000	25.1	-0.158	-3.82	
727	1.22150	5.155	0.03700	25.1	-0.148	-4.82	
757	1.24825	5.019	0.02675	25.1	-0.136	-5.72	
787	1.26925	4.888	0.02100	25.1	-0.131	-6.48	
818	1.28700	4.759	0.01775	25.1	-0.129	-7.19	
848	1.30300	4.628	0.01600	25.1	-0.131	-7.92	
	1.313026 ⁽⁹⁾	4.542 ⁽⁹⁾					
878	1.31775	4.502	0.01475	25.1	-0.126	-8.05	
908	1.33275	4.377	0.01500	25.1	-0.125	-7.43	
939	1.34975	4.249	0.01700	25.1	-0.128	-6.63	
969	1.36950	4.121	0.01975	25.1	-0.128	-5.84	
999	1.39375	3.993	0.02425	25.1	-0.128	-4.94	
1029	1.42525	3.862	0.03150	25.1	-0.131	-3.95	
1060	1.46725	3.728	0.04200	25.1	-0.134	-2.97	
1090	1.51725	3.606	0.05000	25.1	-0.122	-2.21	
1120	1.56725	3.510	0.05000	25.1	-0.096	-1.74	
1150	1.61725	3.431	0.05000	25.1	-0.079	-1.44	
1180	1.66725	3.363	0.05000	25.1	-0.068	-1.23	
1211	1.71725	3.304	0.05000	25.1	-0.059	-1.08	
1241	1.76725	3.253	0.05000	25.1	-0.051	-0.97	
1271	1.81725	3.207	0.05000	25.1	-0.046	-0.87	
1302	1.86725	3.165	0.05000	25.1	-0.042	-0.79	
1332	1.91725	3.127	0.05000	25.1	-0.038		
1362	1.96725	3.092	0.05000	25.1	-0.035		
1392	2.01725	3.059	0.05000	25.1	-0.033		
1422	2.06725	3.029	0.05000	25.1	-0.030		
1453	2.11725	3.000	0.05000	25.1	-0.029		

Measure (normal) [1] Scope 1/1, Sample 25/30

Measured Values Scope 1/1, Sample 25/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.735	23.1	
2	5.042	23.3	0.307
3	5.315	23.6	0.273
4	5.542	23.9	0.227
5	5.743	24.2	0.201
6	5.926	24.4	0.183
7	6.095	24.6	0.169
8	6.254	24.7	0.159
9	6.386	24.8	0.132
10	6.523	24.9	0.137

Titration report

LabX

11	6.652	25.0	0.129
12	6.769	25.1	0.117
13	6.878	25.1	0.109
14	6.975	25.1	0.097
15	7.053	25.2	0.078
16	7.142	25.2	0.089
17	7.219	25.2	0.077
18	7.288	25.2	0.069
19	7.351	25.2	0.063
20	7.408	25.2	0.057
21	7.460	25.2	0.052
22	7.504	25.2	0.044
23	7.548	25.2	0.044
24	7.587	25.2	0.039
25	7.623	25.3	0.036
26	7.657	25.3	0.034
27	7.688	25.3	0.031
28	7.716	25.3	0.028
29	7.740	25.3	0.024
30	7.763	25.3	0.023
31	7.783	25.3	0.020
32	7.802	25.3	0.019
33	7.818	25.3	0.016
34	7.833	25.3	0.015
35	7.847	25.3	0.014
36	7.859	25.3	0.012
37	7.871	25.3	0.012
38	7.881	25.3	0.010
39	7.891	25.3	0.010
40	7.900	25.3	0.009
41	7.908	25.3	0.008
42	7.915	25.3	0.007
43	7.921	25.3	0.006
44	7.928	25.3	0.007
45	7.933	25.3	0.005
46	7.938	25.3	0.005
47	7.942	25.3	0.004
48	7.946	25.3	0.004
49	7.950	25.3	0.004
50	7.954	25.3	0.004
51	7.957	25.3	0.003
52	7.960	25.3	0.003
53	7.963	25.3	0.003
54	7.965	25.3	0.002
55	7.968	25.3	0.003
56	7.970	25.3	0.002
57	7.972	25.3	0.002
58	7.973	25.3	0.001
59	7.975	25.3	0.002
60	7.977	25.3	0.002

Titration report


LabX

Titration (EQP) [1]

Scope 1/1, Sample 25/30


E - V Curve

Scope 1/1, Sample 25/30

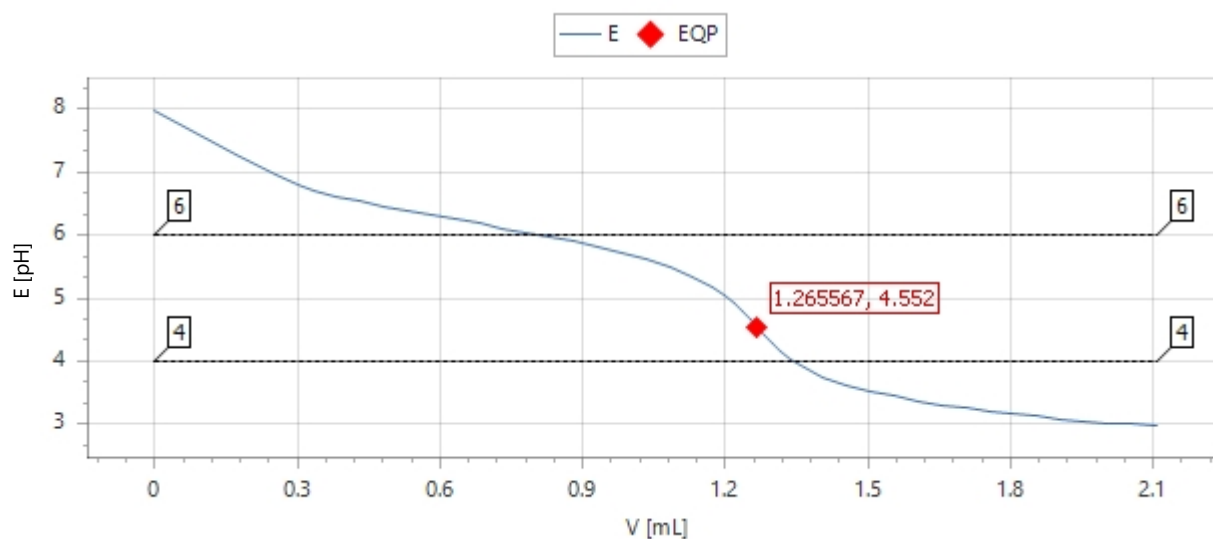

dE/dV - V Curve

Scope 1/1, Sample 25/30

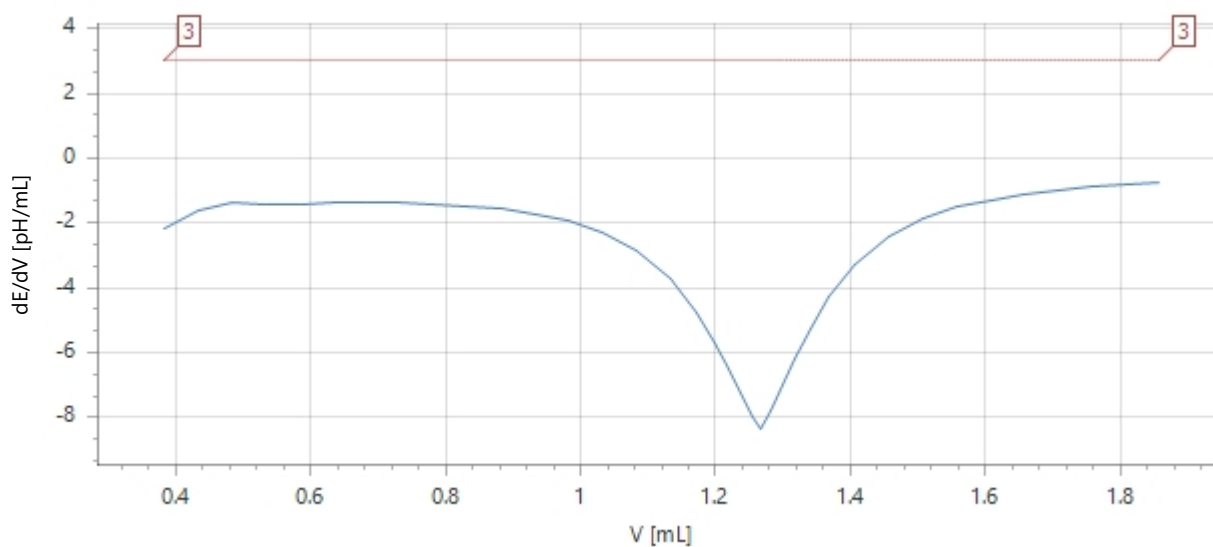

Measured Values

Scope 1/1, Sample 25/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.980		25.3			
31	0.17150	7.281	0.17150	25.2	-0.699		
61	0.25725	6.968	0.08575	25.2	-0.313		
92	0.30000	6.802	0.04275	25.2	-0.166		
183	0.33350	6.694	0.03350	25.2	-0.108		
213	0.38325	6.618	0.04975	25.2	-0.076	-2.17	
243	0.43325	6.539	0.05000	25.2	-0.079	-1.61	
273	0.48325	6.461	0.05000	25.2	-0.078	-1.38	
304	0.53325	6.386	0.05000	25.2	-0.075	-1.40	
334	0.58325	6.316	0.05000	25.2	-0.070	-1.42	

Titration report

LabX

364	0.63325	6.246	0.05000	25.2	-0.070	-1.38	EQP
394	0.68325	6.177	0.05000	25.2	-0.069	-1.36	
424	0.73325	6.109	0.05000	25.2	-0.068	-1.37	
455	0.78325	6.040	0.05000	25.1	-0.069	-1.41	
485	0.83325	5.967	0.05000	25.1	-0.073	-1.47	
515	0.88325	5.892	0.05000	25.1	-0.075	-1.57	
546	0.93325	5.808	0.05000	25.1	-0.084	-1.71	
576	0.98325	5.717	0.05000	25.1	-0.091	-1.93	
606	1.03325	5.612	0.05000	25.1	-0.105	-2.27	
636	1.08325	5.490	0.05000	25.1	-0.122	-2.83	
666	1.13325	5.336	0.05000	25.1	-0.154	-3.74	
697	1.17150	5.183	0.03825	25.1	-0.153	-4.77	
727	1.19800	5.052	0.02650	25.1	-0.131	-5.68	
757	1.22000	4.916	0.02200	25.1	-0.136	-6.52	
787	1.23750	4.790	0.01750	25.1	-0.126	-7.23	
818	1.25350	4.659	0.01600	25.1	-0.131	-8.01	
	1.265567 ⁽⁹⁾	4.552 ⁽⁹⁾					
848	1.26800	4.531	0.01450	25.1	-0.128	-8.33	
878	1.28225	4.406	0.01425	25.1	-0.125	-7.81	
908	1.29800	4.281	0.01575	25.1	-0.125	-6.98	
939	1.31675	4.152	0.01875	25.1	-0.129	-6.18	
969	1.33950	4.023	0.02275	25.1	-0.129	-5.28	
999	1.36825	3.894	0.02875	25.1	-0.129	-4.28	
1029	1.40675	3.760	0.03850	25.1	-0.134	-3.27	
1060	1.45675	3.629	0.05000	25.1	-0.131	-2.39	
1090	1.50675	3.526	0.05000	25.1	-0.103	-1.85	
1120	1.55675	3.443	0.05000	25.1	-0.083	-1.51	
1150	1.60675	3.373	0.05000	25.1	-0.070	-1.28	
1181	1.65675	3.312	0.05000	25.1	-0.061	-1.12	
1211	1.70675	3.258	0.05000	25.1	-0.054	-1.00	
1241	1.75675	3.211	0.05000	25.1	-0.047	-0.90	
1271	1.80675	3.168	0.05000	25.1	-0.043	-0.82	
1302	1.85675	3.129	0.05000	25.1	-0.039	-0.74	
1332	1.90675	3.093	0.05000	25.1	-0.036		
1362	1.95675	3.060	0.05000	25.1	-0.033		
1392	2.00675	3.029	0.05000	25.1	-0.031		
1422	2.05675	3.001	0.05000	25.1	-0.028		
1453	2.10675	2.974	0.05000	25.1	-0.027		

Measure (normal) [1] Scope 1/1, Sample 26/30

Measured Values Scope 1/1, Sample 26/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.667	23.3	
2	4.974	23.5	0.307
3	5.270	23.7	0.296
4	5.514	24.0	0.244
5	5.720	24.2	0.206
6	5.909	24.4	0.189

Titration report

LabX

7	6.083	24.6	0.174
8	6.228	24.7	0.145
9	6.379	24.8	0.151
10	6.519	24.9	0.140
11	6.651	24.9	0.132
12	6.771	25.0	0.120
13	6.887	25.0	0.116
14	6.989	25.0	0.102
15	7.073	25.1	0.084
16	7.161	25.1	0.088
17	7.238	25.1	0.077
18	7.309	25.1	0.071
19	7.372	25.1	0.063
20	7.431	25.1	0.059
21	7.484	25.1	0.053
22	7.526	25.1	0.042
23	7.571	25.1	0.045
24	7.612	25.1	0.041
25	7.649	25.1	0.037
26	7.683	25.1	0.034
27	7.713	25.1	0.030
28	7.740	25.1	0.027
29	7.764	25.1	0.024
30	7.784	25.1	0.020
31	7.804	25.1	0.020
32	7.822	25.1	0.018
33	7.840	25.1	0.018
34	7.856	25.1	0.016
35	7.870	25.1	0.014
36	7.880	25.1	0.010
37	7.892	25.1	0.012
38	7.903	25.1	0.011
39	7.912	25.1	0.009
40	7.920	25.1	0.008
41	7.928	25.1	0.008
42	7.935	25.1	0.007
43	7.941	25.1	0.006
44	7.946	25.1	0.005
45	7.952	25.1	0.006
46	7.957	25.1	0.005
47	7.962	25.1	0.005
48	7.966	25.1	0.004
49	7.970	25.1	0.004
50	7.973	25.1	0.003
51	7.976	25.1	0.003
52	7.979	25.1	0.003
53	7.981	25.1	0.002
54	7.984	25.1	0.003
55	7.986	25.1	0.002
56	7.988	25.1	0.002
57	7.990	25.1	0.002


Titration report

LabX


58	7.991	25.1	0.001
59	7.993	25.1	0.002
60	7.995	25.1	0.002

Titration (EQP) [1] Scope 1/1, Sample 26/30

E - V Curve Scope 1/1, Sample 26/30

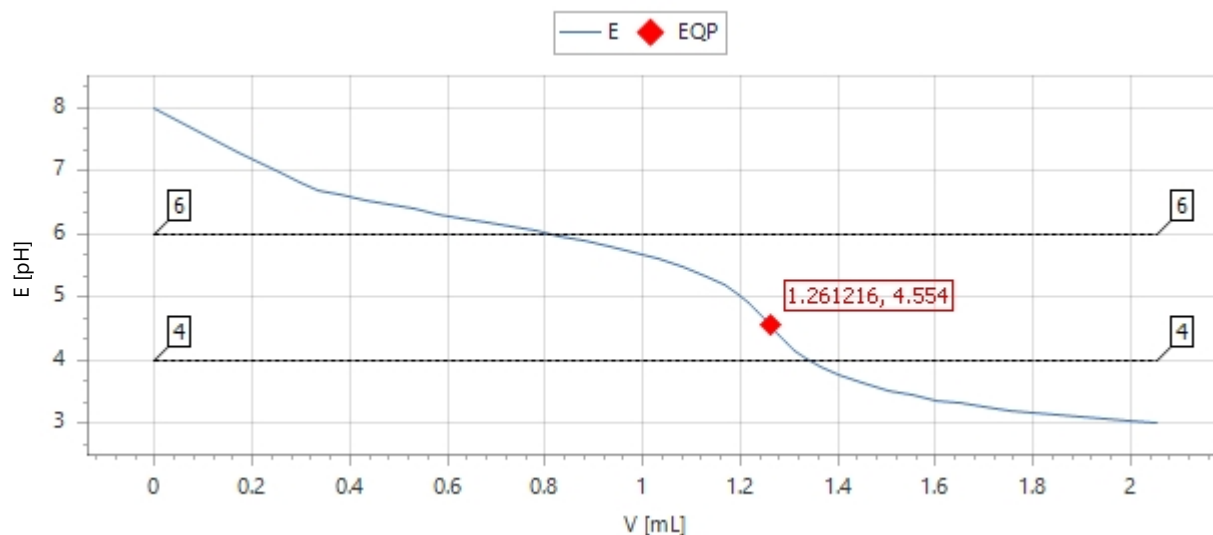

dE/dV - V Curve Scope 1/1, Sample 26/30

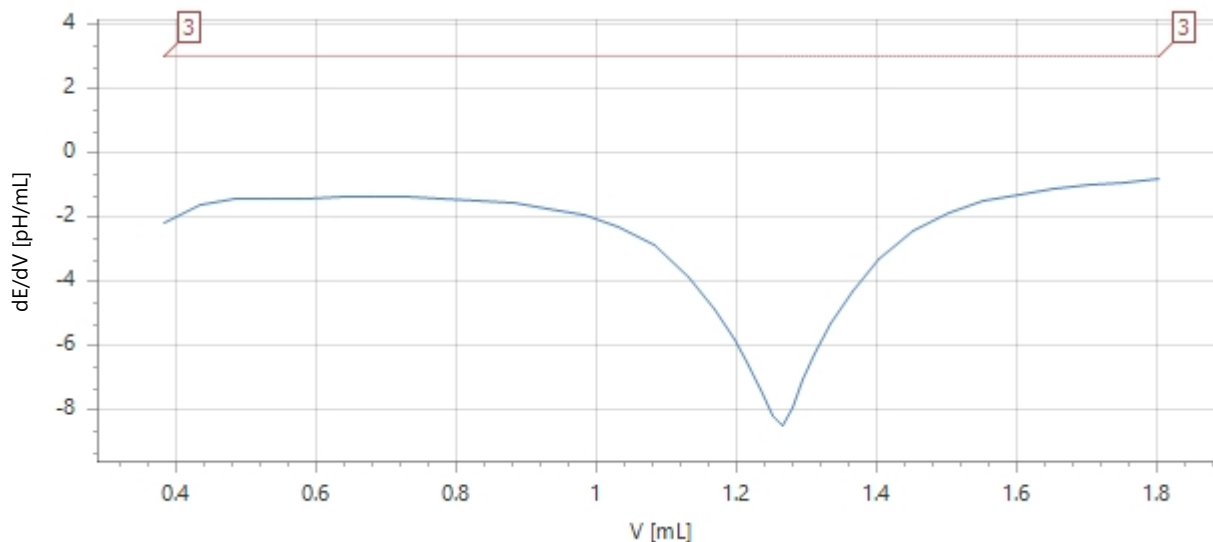

Measured Values Scope 1/1, Sample 26/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.997		25.1			
31	0.17150	7.288	0.17150	25.1	-0.709		
61	0.25725	6.971	0.08575	25.1	-0.317		
91	0.30000	6.806	0.04275	25.1	-0.165		
182	0.33300	6.698	0.03300	25.1	-0.108		
212	0.38200	6.622	0.04900	25.1	-0.076	-2.19	

Titration report

LabX

242	0.43200	6.541	0.05000	25.1	-0.081	-1.62	
272	0.48200	6.464	0.05000	25.1	-0.077	-1.39	
303	0.53200	6.389	0.05000	25.1	-0.075	-1.40	
333	0.58200	6.319	0.05000	25.1	-0.070	-1.41	
363	0.63200	6.248	0.05000	25.1	-0.071	-1.37	
394	0.68200	6.182	0.05000	25.1	-0.066	-1.36	
424	0.73200	6.113	0.05000	25.1	-0.069	-1.37	
454	0.78200	6.042	0.05000	25.1	-0.071	-1.41	
484	0.83200	5.970	0.05000	25.1	-0.072	-1.47	
514	0.88200	5.894	0.05000	25.1	-0.076	-1.57	
545	0.93200	5.811	0.05000	25.1	-0.083	-1.71	
575	0.98200	5.719	0.05000	25.1	-0.092	-1.94	
605	1.03200	5.615	0.05000	25.1	-0.104	-2.30	
636	1.08200	5.490	0.05000	25.1	-0.125	-2.88	
666	1.13200	5.333	0.05000	25.1	-0.157	-3.83	
696	1.16925	5.182	0.03725	25.1	-0.151	-4.87	
726	1.19525	5.050	0.02600	25.1	-0.132	-5.79	
756	1.21650	4.919	0.02125	25.1	-0.131	-6.61	
787	1.23450	4.785	0.01800	25.1	-0.134	-7.37	
817	1.24975	4.656	0.01525	25.1	-0.129	-8.15	
	1.261216 ⁽⁹⁾	4.554 ⁽⁹⁾					EQP
847	1.26375	4.531	0.01400	25.1	-0.125	-8.46	
877	1.27825	4.404	0.01450	25.1	-0.127	-7.90	
908	1.29400	4.277	0.01575	25.1	-0.127	-7.03	
938	1.31225	4.150	0.01825	25.1	-0.127	-6.23	
968	1.33475	4.022	0.02250	25.1	-0.128	-5.32	
998	1.36400	3.891	0.02925	25.1	-0.131	-4.29	
1029	1.40225	3.759	0.03825	25.1	-0.132	-3.27	
1059	1.45225	3.626	0.05000	25.1	-0.133	-2.39	
1089	1.50225	3.524	0.05000	25.1	-0.102	-1.85	
1119	1.55225	3.441	0.05000	25.1	-0.083	-1.51	
1150	1.60225	3.371	0.05000	25.1	-0.070	-1.28	
1180	1.65225	3.310	0.05000	25.1	-0.061	-1.12	
1210	1.70225	3.257	0.05000	25.1	-0.053	-1.00	
1240	1.75225	3.209	0.05000	25.1	-0.048	-0.90	
1270	1.80225	3.166	0.05000	25.1	-0.043	-0.82	
1301	1.85225	3.127	0.05000	25.1	-0.039		
1331	1.90225	3.091	0.05000	25.1	-0.036		
1361	1.95225	3.058	0.05000	25.1	-0.033		
1392	2.00225	3.027	0.05000	25.1	-0.031		
1422	2.05225	2.999	0.05000	25.1	-0.028		

Measure (normal) [1] Scope 1/1, Sample 27/30

Measured Values Scope 1/1, Sample 27/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.713	21.7	
2	5.015	22.0	0.302
3	5.280	22.4	0.265

Titration report

LabX

4	5.501	22.9	0.221
5	5.689	23.3	0.188
6	5.858	23.7	0.169
7	6.014	24.0	0.156
8	6.141	24.2	0.127
9	6.274	24.3	0.133
10	6.397	24.5	0.123
11	6.511	24.6	0.114
12	6.616	24.7	0.105
13	6.714	24.8	0.098
14	6.809	24.8	0.095
15	6.888	24.8	0.079
16	6.964	24.9	0.076
17	7.034	24.9	0.070
18	7.107	24.9	0.073
19	7.173	24.9	0.066
20	7.234	25.0	0.061
21	7.291	25.0	0.057
22	7.336	25.0	0.045
23	7.383	25.0	0.047
24	7.426	25.0	0.043
25	7.467	25.0	0.041
26	7.504	25.0	0.037
27	7.537	25.0	0.033
28	7.569	25.0	0.032
29	7.595	25.0	0.026
30	7.622	25.0	0.027
31	7.647	25.0	0.025
32	7.670	25.0	0.023
33	7.690	25.0	0.020
34	7.707	25.0	0.017
35	7.725	25.0	0.018
36	7.739	25.0	0.014
37	7.753	25.0	0.014
38	7.767	25.0	0.014
39	7.779	25.0	0.012
40	7.790	25.0	0.011
41	7.800	25.0	0.010
42	7.809	25.0	0.009
43	7.817	25.0	0.008
44	7.826	25.0	0.009
45	7.833	25.0	0.007
46	7.840	25.0	0.007
47	7.847	25.0	0.007
48	7.853	25.0	0.006
49	7.859	25.0	0.006
50	7.863	25.0	0.004
51	7.868	25.0	0.005
52	7.872	25.0	0.004
53	7.877	25.0	0.005
54	7.880	25.0	0.003


Titration report

LabX


55	7.884	25.0	0.004
56	7.887	25.0	0.003
57	7.889	25.0	0.002
58	7.892	25.0	0.003
59	7.894	25.0	0.002
60	7.897	25.0	0.003

Titration (EQP) [1] Scope 1/1, Sample 27/30

E - V Curve Scope 1/1, Sample 27/30

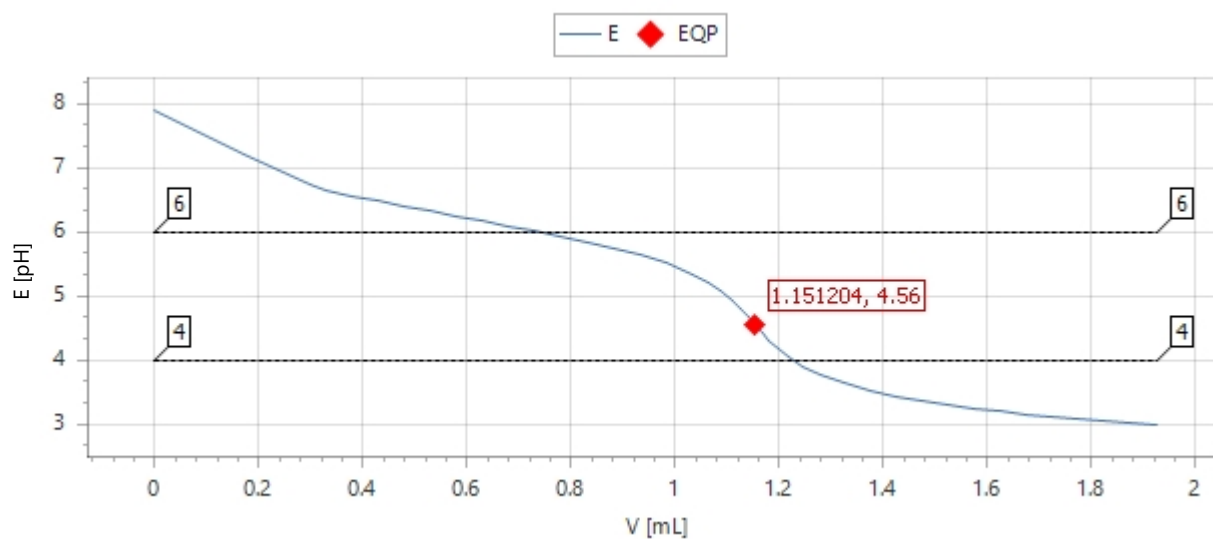

dE/dV - V Curve Scope 1/1, Sample 27/30

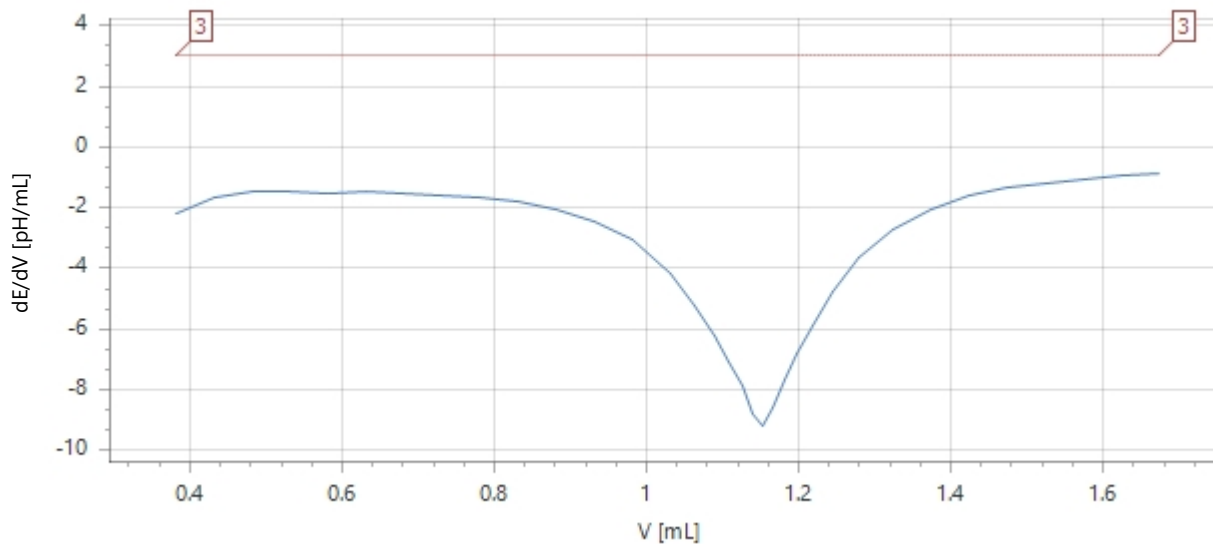

Measured Values Scope 1/1, Sample 27/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.901		25.0			
32	0.17150	7.216	0.17150	25.0	-0.685		
63	0.25725	6.906	0.08575	25.0	-0.310		

Titration report

LabX

93	0.30000	6.737	0.04275	25.0	-0.169	
184	0.33125	6.634	0.03125	25.0	-0.103	
214	0.38125	6.557	0.05000	25.0	-0.077	-2.22
244	0.43125	6.473	0.05000	25.1	-0.084	-1.66
274	0.48125	6.394	0.05000	25.1	-0.079	-1.44
305	0.53125	6.317	0.05000	25.1	-0.077	-1.48
335	0.58125	6.239	0.05000	25.1	-0.078	-1.51
365	0.63125	6.167	0.05000	25.0	-0.072	-1.49
395	0.68125	6.091	0.05000	25.0	-0.076	-1.52
426	0.73125	6.013	0.05000	25.0	-0.078	-1.57
456	0.78125	5.932	0.05000	25.0	-0.081	-1.66
486	0.83125	5.844	0.05000	25.0	-0.088	-1.81
516	0.88125	5.746	0.05000	25.0	-0.098	-2.04
546	0.93125	5.638	0.05000	25.0	-0.108	-2.43
577	0.98125	5.506	0.05000	25.0	-0.132	-3.07
607	1.03125	5.339	0.05000	25.1	-0.167	-4.15
637	1.06425	5.196	0.03300	25.1	-0.143	-5.21
668	1.08850	5.066	0.02425	25.1	-0.130	-6.19
698	1.10875	4.930	0.02025	25.1	-0.136	-7.09
728	1.12500	4.804	0.01625	25.1	-0.126	-7.89
758	1.14000	4.666	0.01500	25.1	-0.138	-8.79
	1.151204 ⁽⁹⁾	4.560 ⁽⁹⁾				
788	1.15250	4.548	0.01250	25.1	-0.118	-9.16
819	1.16650	4.413	0.01400	25.1	-0.135	-8.61
849	1.18025	4.292	0.01375	25.1	-0.121	-7.70
879	1.19725	4.162	0.01700	25.1	-0.130	-6.81
909	1.21750	4.035	0.02025	25.1	-0.127	-5.86
940	1.24400	3.904	0.02650	25.1	-0.131	-4.75
970	1.27850	3.771	0.03450	25.1	-0.133	-3.66
1000	1.32400	3.639	0.04550	25.1	-0.132	-2.69
1030	1.37400	3.527	0.05000	25.1	-0.112	-2.03
1060	1.42400	3.438	0.05000	25.1	-0.089	-1.62
1091	1.47400	3.363	0.05000	25.1	-0.075	-1.35
1121	1.52400	3.299	0.05000	25.1	-0.064	-1.17
1151	1.57400	3.243	0.05000	25.1	-0.056	-1.04
1181	1.62400	3.194	0.05000	25.1	-0.049	-0.93
1212	1.67400	3.149	0.05000	25.0	-0.045	-0.84
1242	1.72400	3.109	0.05000	25.0	-0.040	
1272	1.77400	3.072	0.05000	25.0	-0.037	
1302	1.82400	3.038	0.05000	25.0	-0.034	
1333	1.87400	3.007	0.05000	25.0	-0.031	
1363	1.92400	2.977	0.05000	25.0	-0.030	

EQP

Measure (normal) [1] Scope 1/1, Sample 28/30

Measured Values Scope 1/1, Sample 28/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.712	22.0	
2	5.007	22.3	0.295

Titration report

LabX

3	5.287	22.7	0.280
4	5.516	23.1	0.229
5	5.705	23.5	0.189
6	5.870	23.8	0.165
7	6.024	24.1	0.154
8	6.167	24.3	0.143
9	6.287	24.4	0.120
10	6.411	24.6	0.124
11	6.528	24.7	0.117
12	6.641	24.8	0.113
13	6.744	24.8	0.103
14	6.840	24.9	0.096
15	6.925	24.9	0.085
16	6.995	24.9	0.070
17	7.067	25.0	0.072
18	7.139	25.0	0.072
19	7.203	25.0	0.064
20	7.262	25.0	0.059
21	7.316	25.0	0.054
22	7.362	25.0	0.046
23	7.409	25.0	0.047
24	7.452	25.0	0.043
25	7.491	25.0	0.039
26	7.527	25.0	0.036
27	7.561	25.0	0.034
28	7.592	25.0	0.031
29	7.617	25.0	0.025
30	7.643	25.0	0.026
31	7.666	25.0	0.023
32	7.688	25.0	0.022
33	7.709	25.0	0.021
34	7.728	25.0	0.019
35	7.745	25.0	0.017
36	7.759	25.0	0.014
37	7.774	25.0	0.015
38	7.787	25.0	0.013
39	7.800	25.0	0.013
40	7.812	25.0	0.012
41	7.823	25.0	0.011
42	7.832	25.0	0.009
43	7.840	25.0	0.008
44	7.847	25.0	0.007
45	7.855	25.0	0.008
46	7.861	25.0	0.006
47	7.868	25.0	0.007
48	7.873	25.0	0.005
49	7.879	25.0	0.006
50	7.883	25.0	0.004
51	7.887	25.0	0.004
52	7.891	25.0	0.004
53	7.895	25.0	0.004


Titration report

LabX


54	7.899	25.0	0.004
55	7.902	25.0	0.003
56	7.905	25.0	0.003
57	7.907	25.0	0.002
58	7.910	25.0	0.003
59	7.912	25.0	0.002
60	7.914	25.0	0.002

Titration (EQP) [1] Scope 1/1, Sample 28/30

E - V Curve Scope 1/1, Sample 28/30

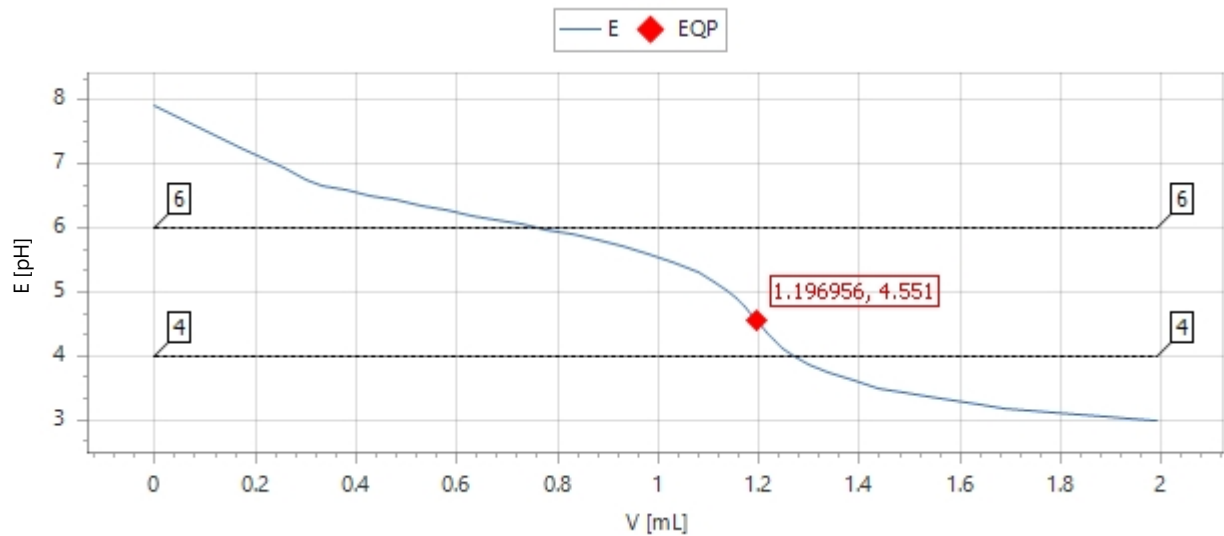

dE/dV - V Curve Scope 1/1, Sample 28/30

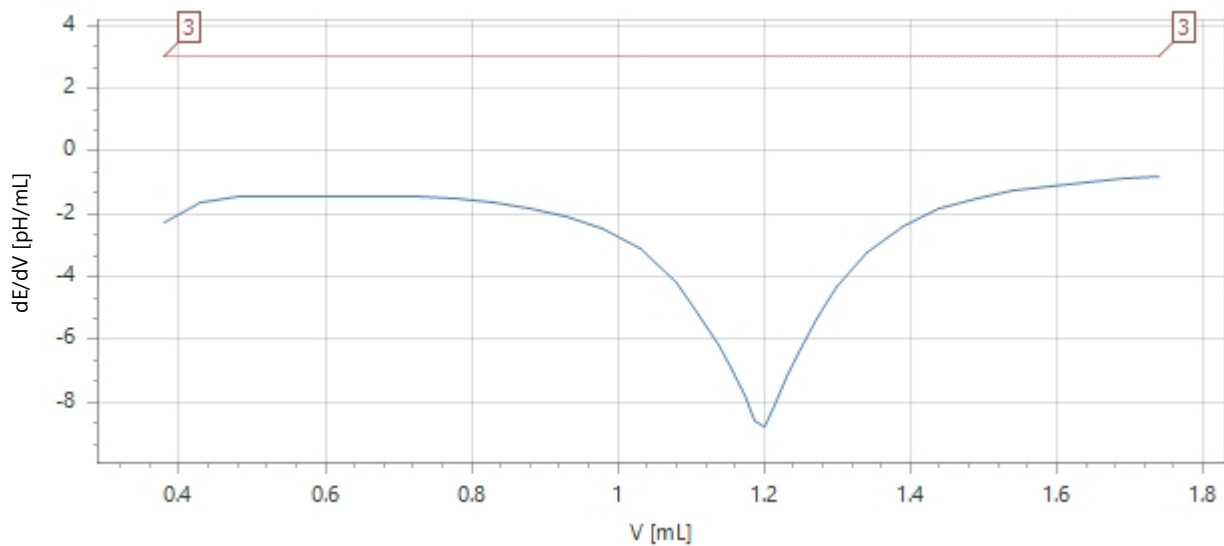

Measured Values Scope 1/1, Sample 28/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.919		25.0			
31	0.17150	7.248	0.17150	25.0	-0.671		

Titration report

LabX

61	0.25725	6.934	0.08575	25.0	-0.314	
91	0.30000	6.763	0.04275	25.0	-0.171	
182	0.33100	6.659	0.03100	25.1	-0.104	
212	0.38000	6.584	0.04900	25.1	-0.075	-2.24
242	0.43000	6.504	0.05000	25.1	-0.080	-1.65
272	0.48000	6.424	0.05000	25.1	-0.080	-1.42
303	0.53000	6.348	0.05000	25.1	-0.076	-1.45
333	0.58000	6.273	0.05000	25.1	-0.075	-1.47
363	0.63000	6.199	0.05000	25.1	-0.074	-1.44
393	0.68000	6.131	0.05000	25.1	-0.068	-1.44
424	0.73000	6.056	0.05000	25.1	-0.075	-1.47
454	0.78000	5.981	0.05000	25.1	-0.075	-1.53
484	0.83000	5.902	0.05000	25.1	-0.079	-1.63
514	0.88000	5.814	0.05000	25.1	-0.088	-1.80
545	0.93000	5.719	0.05000	25.1	-0.095	-2.05
575	0.98000	5.607	0.05000	25.1	-0.112	-2.45
605	1.03000	5.475	0.05000	25.1	-0.132	-3.11
635	1.08000	5.304	0.05000	25.1	-0.171	-4.21
666	1.11150	5.168	0.03150	25.1	-0.136	-5.20
696	1.13650	5.033	0.02500	25.1	-0.135	-6.17
726	1.15650	4.899	0.02000	25.1	-0.134	-7.01
756	1.17300	4.769	0.01650	25.1	-0.130	-7.77
786	1.18775	4.636	0.01475	25.1	-0.133	-8.61
	1.196956 ⁽⁹⁾	4.551 ⁽⁹⁾				EQP
817	1.20075	4.516	0.01300	25.1	-0.120	-8.75
847	1.21525	4.384	0.01450	25.1	-0.132	-8.06
877	1.23050	4.260	0.01525	25.1	-0.124	-7.18
908	1.24925	4.129	0.01875	25.1	-0.131	-6.31
938	1.27150	4.001	0.02225	25.1	-0.128	-5.35
968	1.30050	3.872	0.02900	25.1	-0.129	-4.30
998	1.33950	3.738	0.03900	25.1	-0.134	-3.25
1028	1.38950	3.607	0.05000	25.1	-0.131	-2.37
1059	1.43950	3.506	0.05000	25.1	-0.101	-1.84
1089	1.48950	3.423	0.05000	25.1	-0.083	-1.50
1119	1.53950	3.353	0.05000	25.1	-0.070	-1.27
1149	1.58950	3.293	0.05000	25.1	-0.060	-1.11
1180	1.63950	3.240	0.05000	25.1	-0.053	-0.99
1210	1.68950	3.193	0.05000	25.1	-0.047	-0.89
1240	1.73950	3.150	0.05000	25.1	-0.043	-0.81
1270	1.78950	3.111	0.05000	25.1	-0.039	
1301	1.83950	3.075	0.05000	25.1	-0.036	
1331	1.88950	3.042	0.05000	25.1	-0.033	
1361	1.93950	3.012	0.05000	25.1	-0.030	
1391	1.98950	2.983	0.05000	25.1	-0.029	

Measure (normal) [1] Scope 1/1, Sample 29/30

Measured Values Scope 1/1, Sample 29/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
-------------	---------------------	---------------	------------

Titration report

LabX

1	4.701	23.3	
2	4.994	23.5	0.293
3	5.269	23.7	0.275
4	5.488	24.0	0.219
5	5.672	24.2	0.184
6	5.837	24.4	0.165
7	5.988	24.6	0.151
8	6.113	24.7	0.125
9	6.244	24.8	0.131
10	6.367	24.9	0.123
11	6.482	25.0	0.115
12	6.588	25.0	0.106
13	6.684	25.0	0.096
14	6.776	25.1	0.092
15	6.857	25.1	0.081
16	6.934	25.1	0.077
17	7.005	25.1	0.071
18	7.076	25.1	0.071
19	7.142	25.2	0.066
20	7.199	25.2	0.057
21	7.253	25.2	0.054
22	7.298	25.2	0.045
23	7.345	25.2	0.047
24	7.386	25.2	0.041
25	7.425	25.2	0.039
26	7.462	25.2	0.037
27	7.495	25.2	0.033
28	7.526	25.2	0.031
29	7.550	25.2	0.024
30	7.577	25.2	0.027
31	7.601	25.2	0.024
32	7.623	25.2	0.022
33	7.644	25.2	0.021
34	7.663	25.2	0.019
35	7.680	25.2	0.017
36	7.696	25.2	0.016
37	7.711	25.2	0.015
38	7.725	25.2	0.014
39	7.738	25.2	0.013
40	7.750	25.2	0.012
41	7.761	25.2	0.011
42	7.772	25.2	0.011
43	7.780	25.2	0.008
44	7.789	25.2	0.009
45	7.797	25.2	0.008
46	7.804	25.2	0.007
47	7.811	25.2	0.007
48	7.817	25.2	0.006
49	7.824	25.2	0.007
50	7.829	25.2	0.005
51	7.834	25.2	0.005


Titration report

LabX


52	7.839	25.2	0.005
53	7.844	25.2	0.005
54	7.848	25.2	0.004
55	7.852	25.2	0.004
56	7.855	25.2	0.003
57	7.858	25.2	0.003
58	7.862	25.2	0.004
59	7.864	25.2	0.002
60	7.867	25.2	0.003

Titration (EQP) [1] Scope 1/1, Sample 29/30

E - V Curve Scope 1/1, Sample 29/30

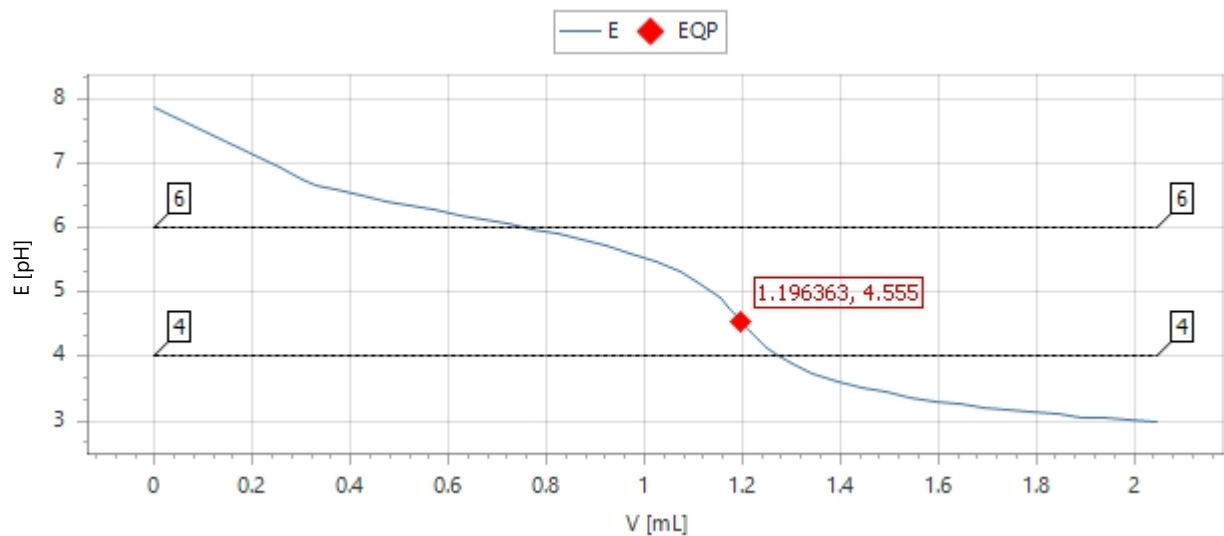

dE/dV - V Curve Scope 1/1, Sample 29/30

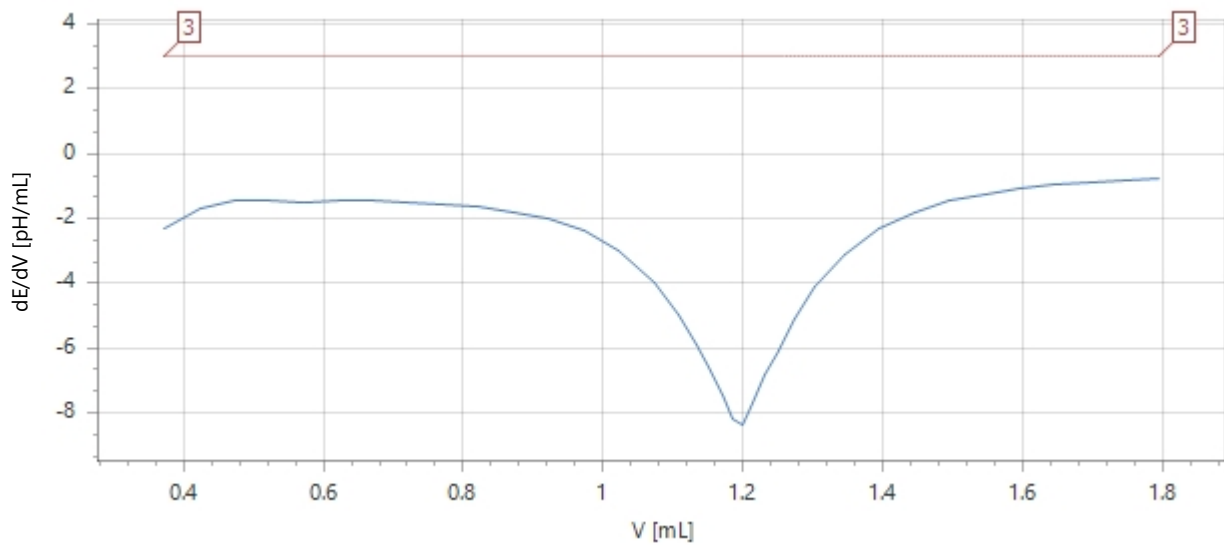

Measured Values Scope 1/1, Sample 29/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
----------	--------	------------------	---------	------------	---------	---------------	-----

Titration report

LabX

0	0.00000	7.873		25.2		
31	0.17150	7.234	0.17150	25.2	-0.639	
61	0.25725	6.924	0.08575	25.2	-0.310	
91	0.30000	6.754	0.04275	25.2	-0.170	
183	0.33125	6.641	0.03125	25.2	-0.113	
213	0.37250	6.577	0.04125	25.2	-0.064	-2.33
243	0.42250	6.500	0.05000	25.2	-0.077	-1.69
273	0.47250	6.421	0.05000	25.2	-0.079	-1.41
304	0.52250	6.344	0.05000	25.2	-0.077	-1.43
334	0.57250	6.270	0.05000	25.2	-0.074	-1.47
364	0.62250	6.198	0.05000	25.1	-0.072	-1.45
394	0.67250	6.125	0.05000	25.1	-0.073	-1.45
425	0.72250	6.052	0.05000	25.1	-0.073	-1.49
455	0.77250	5.976	0.05000	25.1	-0.076	-1.55
485	0.82250	5.896	0.05000	25.1	-0.080	-1.64
515	0.87250	5.808	0.05000	25.1	-0.088	-1.78
546	0.92250	5.714	0.05000	25.1	-0.094	-2.01
576	0.97250	5.606	0.05000	25.1	-0.108	-2.37
606	1.02250	5.478	0.05000	25.1	-0.128	-2.97
636	1.07250	5.316	0.05000	25.1	-0.162	-3.97
667	1.10750	5.169	0.03500	25.1	-0.147	-4.98
697	1.13275	5.039	0.02525	25.1	-0.130	-5.88
727	1.15425	4.901	0.02150	25.1	-0.138	-6.70
757	1.17125	4.773	0.01700	25.1	-0.128	-7.42
788	1.18650	4.644	0.01525	25.1	-0.129	-8.18
	1.196363 ⁽⁹⁾	4.555 ⁽⁹⁾				EQP
818	1.20075	4.516	0.01425	25.1	-0.128	-8.35
848	1.21500	4.396	0.01425	25.1	-0.120	-7.71
878	1.23175	4.263	0.01675	25.1	-0.133	-6.85
908	1.25025	4.140	0.01850	25.1	-0.123	-6.06
939	1.27425	4.007	0.02400	25.1	-0.133	-5.10
969	1.30375	3.881	0.02950	25.1	-0.126	-4.11
999	1.34500	3.746	0.04125	25.1	-0.135	-3.09
1030	1.39500	3.620	0.05000	25.1	-0.126	-2.28
1060	1.44500	3.522	0.05000	25.1	-0.098	-1.78
1090	1.49500	3.441	0.05000	25.1	-0.081	-1.46
1120	1.54500	3.373	0.05000	25.1	-0.068	-1.24
1150	1.59500	3.314	0.05000	25.1	-0.059	-1.09
1181	1.64500	3.262	0.05000	25.1	-0.052	-0.97
1211	1.69500	3.216	0.05000	25.1	-0.046	-0.88
1241	1.74500	3.174	0.05000	25.1	-0.042	-0.80
1271	1.79500	3.135	0.05000	25.1	-0.039	-0.73
1302	1.84500	3.100	0.05000	25.1	-0.035	
1332	1.89500	3.067	0.05000	25.1	-0.033	
1362	1.94500	3.037	0.05000	25.1	-0.030	
1392	1.99500	3.009	0.05000	25.1	-0.028	
1423	2.04500	2.982	0.05000	25.1	-0.027	

Measure (normal) [1] Scope 1/1, Sample 30/30

Titration report

LabX

Measured Values Scope 1/1, Sample 30/30

Time [s]	Meas. value [pH]	Temp. [°C]	dE [pH]
1	4.678	23.0	
2	4.998	23.2	0.320
3	5.286	23.5	0.288
4	5.512	23.8	0.226
5	5.703	24.1	0.191
6	5.875	24.3	0.172
7	6.030	24.5	0.155
8	6.157	24.7	0.127
9	6.290	24.8	0.133
10	6.415	24.9	0.125
11	6.530	24.9	0.115
12	6.638	25.0	0.108
13	6.738	25.0	0.100
14	6.834	25.1	0.096
15	6.911	25.1	0.077
16	6.989	25.1	0.078
17	7.059	25.1	0.070
18	7.131	25.1	0.072
19	7.194	25.2	0.063
20	7.252	25.2	0.058
21	7.305	25.2	0.053
22	7.352	25.2	0.047
23	7.399	25.2	0.047
24	7.441	25.2	0.042
25	7.480	25.2	0.039
26	7.516	25.2	0.036
27	7.548	25.2	0.032
28	7.578	25.2	0.030
29	7.604	25.2	0.026
30	7.630	25.2	0.026
31	7.653	25.2	0.023
32	7.674	25.2	0.021
33	7.695	25.2	0.021
34	7.714	25.2	0.019
35	7.730	25.2	0.016
36	7.744	25.2	0.014
37	7.759	25.2	0.015
38	7.772	25.2	0.013
39	7.785	25.2	0.013
40	7.796	25.2	0.011
41	7.807	25.2	0.011
42	7.816	25.2	0.009
43	7.824	25.2	0.008
44	7.833	25.2	0.009
45	7.840	25.2	0.007
46	7.846	25.2	0.006
47	7.853	25.2	0.007


Titration report

LabX


48	7.859	25.2	0.006
49	7.864	25.2	0.005
50	7.869	25.2	0.005
51	7.873	25.2	0.004
52	7.877	25.2	0.004
53	7.880	25.2	0.003
54	7.884	25.2	0.004
55	7.887	25.2	0.003
56	7.890	25.2	0.003
57	7.893	25.2	0.003
58	7.895	25.2	0.002
59	7.898	25.2	0.003
60	7.900	25.2	0.002

Titration (EQP) [1] Scope 1/1, Sample 30/30

E - V Curve Scope 1/1, Sample 30/30

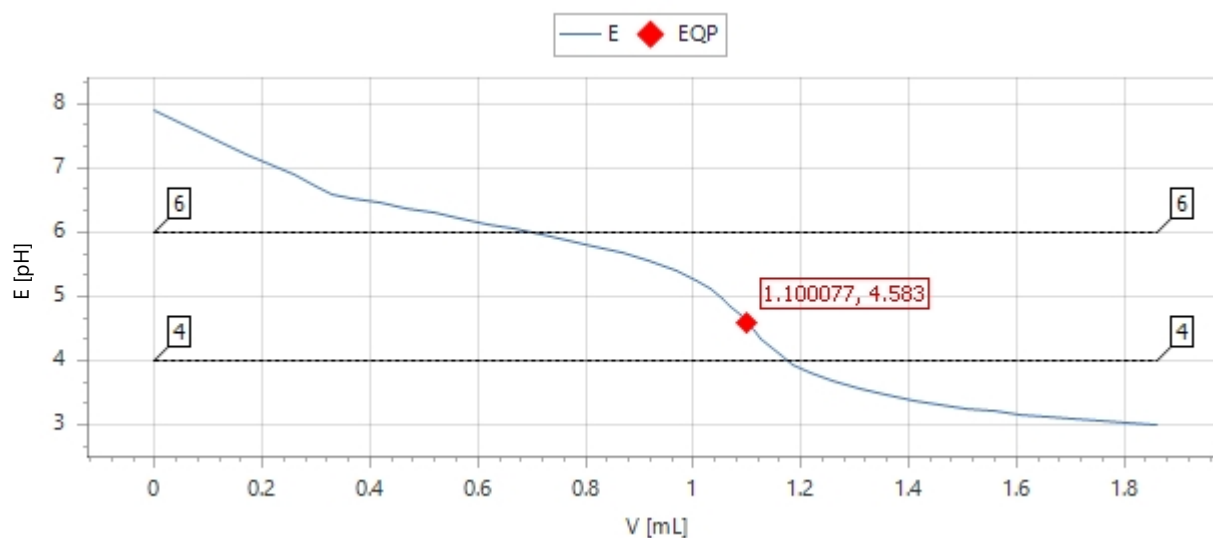

dE/dV - V Curve Scope 1/1, Sample 30/30

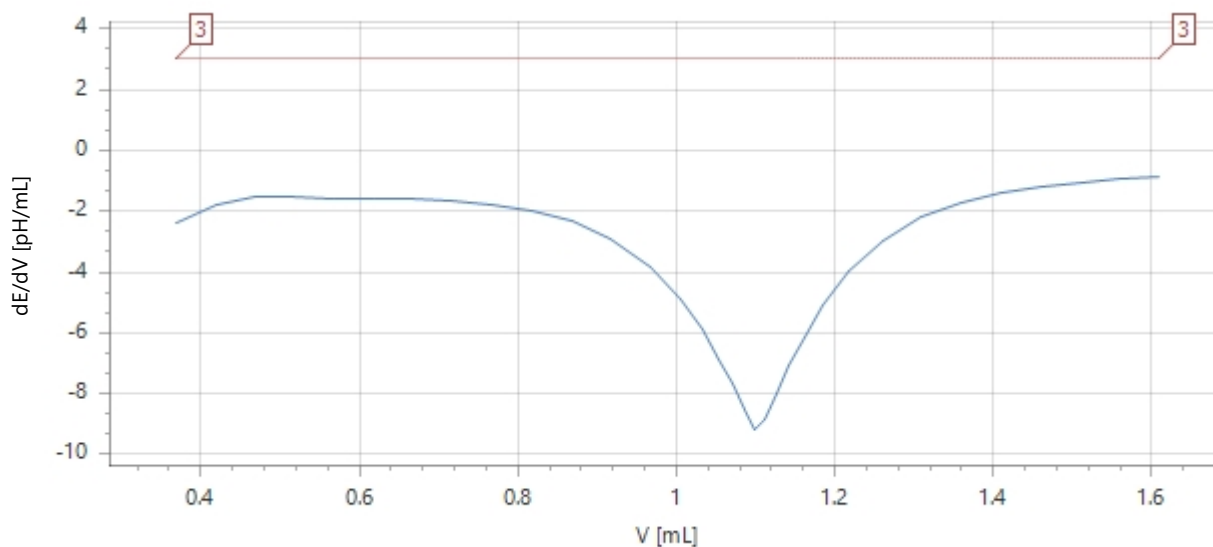

Titration report

LabX

Measured Values Scope 1/1, Sample 30/30

Time [s]	V [mL]	Meas. value [pH]	dV [mL]	Temp. [°C]	dE [pH]	dE/dV [pH/mL]	EQP
0	0.00000	7.906		25.2			
33	0.17150	7.204	0.17150	25.2	-0.702		
63	0.25725	6.891	0.08575	25.2	-0.313		
94	0.30000	6.714	0.04275	25.2	-0.177		
185	0.32925	6.604	0.02925	25.2	-0.110		
215	0.36825	6.540	0.03900	25.2	-0.064	-2.41	
245	0.41825	6.459	0.05000	25.2	-0.081	-1.80	
276	0.46825	6.377	0.05000	25.2	-0.082	-1.51	
306	0.51825	6.294	0.05000	25.2	-0.083	-1.54	
336	0.56825	6.214	0.05000	25.2	-0.080	-1.60	
366	0.61825	6.133	0.05000	25.2	-0.081	-1.60	
396	0.66825	6.055	0.05000	25.2	-0.078	-1.62	
427	0.71825	5.972	0.05000	25.2	-0.083	-1.68	
457	0.76825	5.882	0.05000	25.2	-0.090	-1.80	
487	0.81825	5.788	0.05000	25.2	-0.094	-2.01	
517	0.86825	5.680	0.05000	25.2	-0.108	-2.34	
548	0.91825	5.554	0.05000	25.2	-0.126	-2.89	
578	0.96825	5.396	0.05000	25.2	-0.158	-3.82	
608	1.00525	5.249	0.03700	25.2	-0.147	-4.88	
638	1.03200	5.115	0.02675	25.2	-0.134	-5.89	
669	1.05375	4.978	0.02175	25.2	-0.137	-6.84	
699	1.07100	4.849	0.01725	25.2	-0.129	-7.65	
729	1.08625	4.714	0.01525	25.2	-0.135	-8.55	
759	1.09925	4.591	0.01300	25.2	-0.123	-9.15	
	1.100077	4.583					EQP
790	1.11275	4.461	0.01350	25.2	-0.130	-8.87	
820	1.12675	4.333	0.01400	25.2	-0.128	-7.97	
850	1.14225	4.209	0.01550	25.2	-0.124	-7.09	
880	1.16150	4.080	0.01925	25.2	-0.129	-6.16	
911	1.18575	3.950	0.02425	25.2	-0.130	-5.09	
941	1.21725	3.820	0.03150	25.2	-0.130	-3.99	
971	1.25975	3.685	0.04250	25.2	-0.135	-2.95	
1001	1.30975	3.565	0.05000	25.2	-0.120	-2.18	
1032	1.35975	3.470	0.05000	25.2	-0.095	-1.72	
1062	1.40975	3.392	0.05000	25.2	-0.078	-1.42	
1092	1.45975	3.325	0.05000	25.2	-0.067	-1.22	
1122	1.50975	3.267	0.05000	25.2	-0.058	-1.07	
1152	1.55975	3.216	0.05000	25.2	-0.051	-0.96	
1183	1.60975	3.171	0.05000	25.2	-0.045	-0.86	
1213	1.65975	3.129	0.05000	25.2	-0.042		
1243	1.70975	3.092	0.05000	25.2	-0.037		
1274	1.75975	3.056	0.05000	25.2	-0.036		
1304	1.80975	3.025	0.05000	25.1	-0.031		
1334	1.85975	2.995	0.05000	25.1	-0.030		

Used Instruments

Titration report



Name SERCtitrate **Serial number** B433917554

Responsible User

Time	Login name	Full name
2/21/2024 11:22:53 AM	SystemInternal	System internal account

- (1) Modified
- (9) Standard evaluation used
- (22) New (Manual)